



GC SMART GUIDE

CLICK TO SELECT A PRODUCT CATEGORY




Preventive &
Remineralisation



Glass Hybrid
Restorative System




Conventional & Resin-
Modified Glass Ionomers

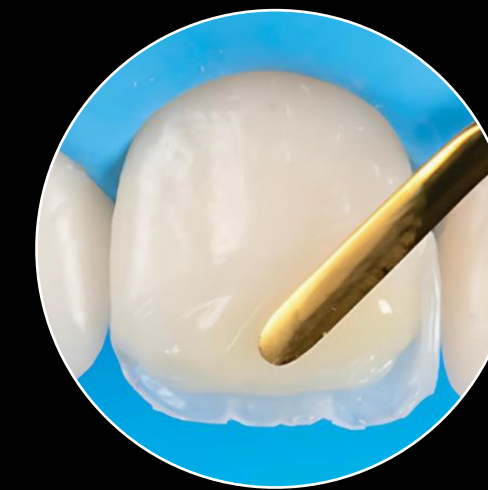

Composites, Bonds
Finishing & Polishing



Luting
Cements

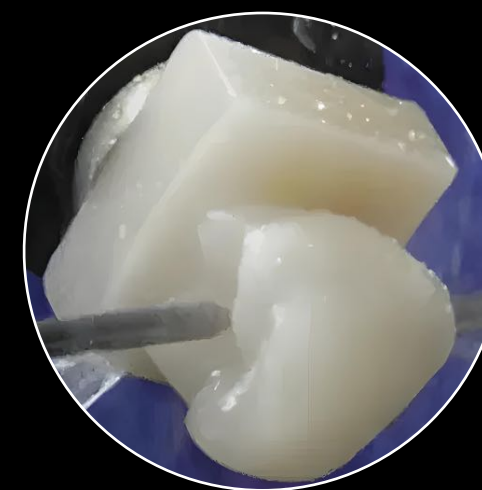

CAD/CAM &
Temporisation



Impression &
Reline



GC Smart
Guide Instructions



Contact a GC
Product Specialist



CPD Newsletter



GC Australasia
Website



About GC

Follow GC On Social Media



GC PREVENTIVE PRODUCTS

Click a product image to view product information



Tri Plaque ID Gel



Dry Mouth Gel



GC Tooth Mousse & Tooth Mousse Plus



RECALDENT® Gum



GC Ruscello Floss

ALWAYS READ THE LABEL AND FOLLOW THE DIRECTIONS FOR USE




GC Tri Plaque ID Gel



Clinical Indications

- Differentiate between new, old and mature plaque within the biofilm
- Identify active **high risk** acid producing bacteria where pH is below 4.5
- Improve patient motivation and education
- Compatible with guided biofilm removal protocols
- Identify high risk acid producing bacteria where pH is below 4.5
- Enables more targeted oral hygiene instruction
- Easy to apply and wash off

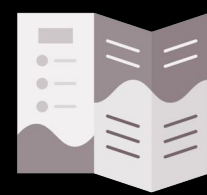


[View product brochure](#)




[View instructions for use](#)

How does Tri Plaque ID Gel work?

GC Tri Plaque ID Gel contains sucrose and pigments that are able to penetrate and stain the plaque biofilm

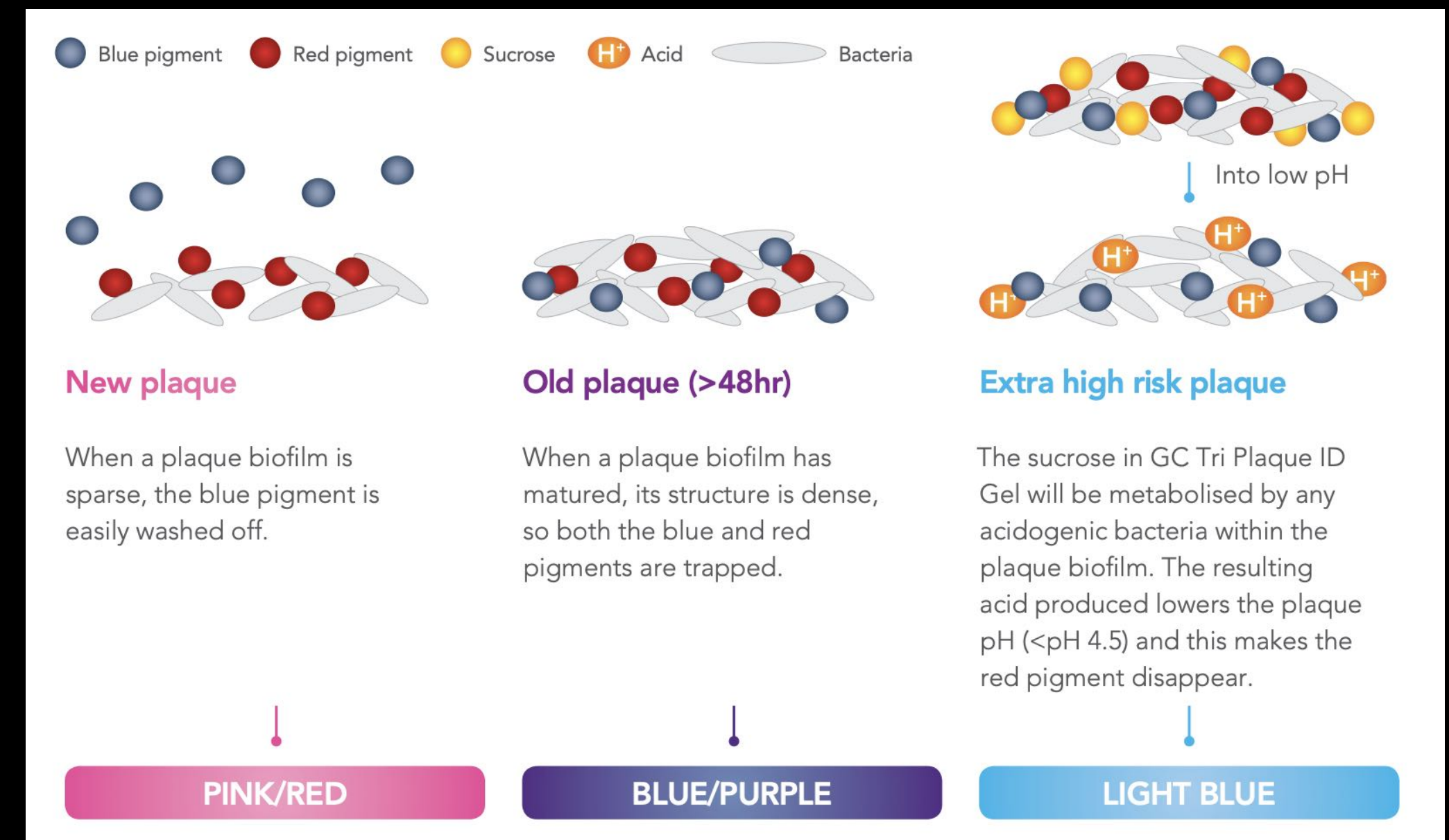

Prof. Ian Meyers

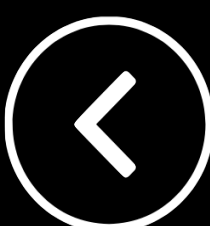


GC Dry Mouth Gel

Dry Mouth Gel is a unique, sugar free product designed to help relieve dry mouths and provide long lasting comfort and a soothing effect for patients.

Uniquely, unlike most saliva substitutes, Dry Mouth Gel has a neutral pH, meaning that effective symptomatic relief can be provided, whilst preserving oral pH within the safe range to prevent demineralisation.

Clinical Indications

- Provision of comfort and relief for dry mouth conditions
- Provides a protective coating for teeth and oral tissues
- Control of halitosis
- Up to 4 hours relief from symptoms of dry mouth
- Sugar free & near neutral pH
- Five flavours - Raspberry, Mint, Orange, Lemon and Fruit Salad
- Doesn't interfere with fluoride products or GC Tooth Mousse

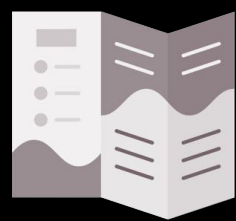

[View product brochure](#)



[View instructions for use](#)



Dry Mouth Risk Factors Examples

Dehydration, ageing, anxiety / stress, mouth breathing, snoring, tobacco

Systemic Conditions Examples

Diabetes, stroke, Alzheimer's disease, Sjogren's syndrome, Salivary gland hypofunction, rheumatoid arthritis, nerve damage to your head or neck

Medication Related Examples


Chemotherapy or radiotherapy, antidepressants, antihistamines, decongestants, diuretics



[View On-Demand Webinar \(1 CPD\)](#)



[View Dry Mouth Gel Patient Video](#)

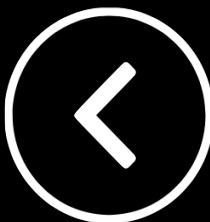

GC Tooth Mousse™ & GC Tooth™ Mousse Plus



RECALDENT and RECALDENT logo are trademarks of Mondelēz International group used under license.

Clinical Indications, Features and Benefits

- Rapidly acid neutralisation & buffering from plaque bacteria, reflux, food and drinks
- Increases salivary concentration of bioavailable calcium, phosphate and fluoride
- Promotes remineralisation to protect, strengthen and brighten enamel
- Prevents tooth decay & erosion
- Prevents dentine sensitivity from erosion, gingival recession and tooth whitening
- Prevents and reverses white spot lesions during *and* after orthodontic treatment
- Contains 10% RECALDENT™ (CPP-ACP) - derived from milk
- Tooth Mousse is suitable for children and adults of all ages
- Tooth Mousse Plus contains 900 ppm fluoride (suitable for adults and children over the age of 6)

Seven delicious flavours and up to 160 applications per tube provides patients with an economical and tasty way to protect their oral health



View Tooth Mousse patient brochure



View Tooth Mousse Plus patient brochure



View Tooth Mousse instructions for use



View Tooth Mousse Plus instructions for use



View Tooth Mousse FAQ's



View Tooth Mousse Patient Video



View Tooth Mousse Plus Patient Video



View Tooth Mousse Plus For Ortho Patient Video



View On-Demand Clinical Presentation (1 CPD)

The higher the risk, the greater the protection that is required



Min. once a day after brushing
e.g. prevention, children, adolescents, adults, elderly,



Twice a day after brushing
e.g. orthodontics, poor diet, poor oral hygiene, pregnancy, tooth whitening, caries risk, WSL



As often as required
e.g. altered saliva flow, xerostomia, medications, erosion, poly-pharmacy, radiotherapy, chemotherapy, periodontal disease, high caries risk, hypomineralisation

Contraindications - not to be used by patients with known or suspected allergy to

- Milk Protein (casein only)
- Benzoates (preservative)
- Soy bean (doesn't apply to Milkshake Flavour)

www.toothmousse.com.au



Recaldent® Gum

Chewing Recaldent® Advanced sugar-free gum with long lasting flavours helps stimulate saliva flow. Saliva keeps our mouths healthy and is nature's way of protecting our teeth.

Benefits

Stimulates saliva flow

Sugar-free


With RECALDENT (CPP-ACP) casein phosphopeptide-amorphous calcium phosphate



Do not use if you have allergy to cow's milk proteins or soybean derived products.

Recaldent® Advanced chewing gum contains PHENYLALANINE.


Excess consumption may have a laxative effect.

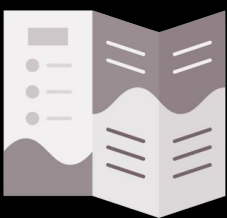

[View product brochure](#)



RECALDENT and RECALDENT logo are trademarks of Mondelēz International group used under license.

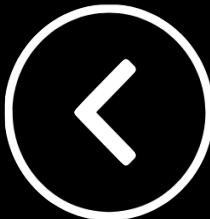

GC Ruscello Floss

Features and Benefits

GC Ruscello - brighten your day with its colours and freshness and look after your teeth with its wedge shaped fibre innovation

GC Ruscello floss uses micro-fibre technology - each fibre is approximately 1/100th the thickness of human hair

Ruscello has 5 times as many fibres as existing floss, making it easy to reach every corner in the mouth

Ruscello fibres are wedge-shaped: plaque is easily taken up into its micro-porosities and wiped away

Ruscello floss has exceptionally soft texture - it does not hurt even when attached to the gingiva

A waxed coating enables the floss to slide effortlessly between teeth

GC Ruscello floss has a great mint flavour for a refreshing finish

Two options - waxed and unwaxed

Also available in 30m and 200m




30m Mint Waxed



200m Mint Waxed



200m Mint Unwaxed

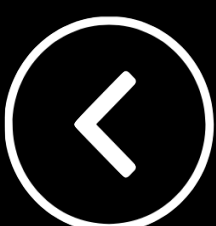

GC GLASS HYBRID RESTORATIVE SYSTEM

Click a product image to view product information



EQUIA Forte HT Fil



EQUIA Forte Coat

ALWAYS READ THE LABEL AND FOLLOW THE DIRECTIONS FOR USE



EQUIA Forte HT Fil



Clinical Indications

- Class I restorations
- Stress bearing Class II restorations
- Non-stress bearing Class II restorations
- Intermediate restorative
- Class V and root surface restorations
- Core build up
- Amalgam replacement
- Molar Hypomineralisation
- High caries risk patients

EQUIA Forte® HT Fil is a long-term bulk fill glass-hybrid restorative system with enhanced mechanical properties, superior fluoride release, sound marginal seal and excellent handling

Glass hybrid systems offer the strength of the latest innovative glasses combined with the wear resistance of a nano-filled coating (Equia Forte Coat). The unique synergy between the coating and the restorative material improves the durability and wear resistance of the restorations.



Dr. Zeynep Bilge Kütük, Turkey

Clinical Features & Benefits

- Packable, non-sticky and moisture-tolerant - fast and easy bulk placement
- Enables posterior restoration placement with more natural aesthetics
- Excellent biocompatibility & biomimetic properties
- Hydrophilic - withstand the humidity of the oral environment and the pulpal fluid flow
- Chemically bonds to dentin, enamel and cementum - strong, stable and chemically-fused seal
- Pulp-friendliness - suitable for use in deep cavities
- Ion exchange - uptake from saliva & aqueous solutions applied during brushing and tooth paste
- Co-efficient of thermal expansion close to tooth structure and the low shrinkage behaviour
- Good marginal adaptation of glass hybrids - virtually no post-operative tooth sensitivity
- Use in conjunction with Equia Forte Coat to increases surface hardness and wear resistance



[View product brochure](#)



[View comprehensive guide](#)



[View instructions for use](#)



[View Technique Guide](#)



[View Schwendicke Paper](#)




[View Schwendicke Article](#)





EQUIA Forte Coat®

EQUIA Forte® Coat incorporates a new multi-functional monomer that produces a tougher resin matrix for better protection. Used together with EQUIA Forte® HT Fil, the result is a strong bulk fill glass hybrid restorative system with exceptional physical properties and aesthetics.

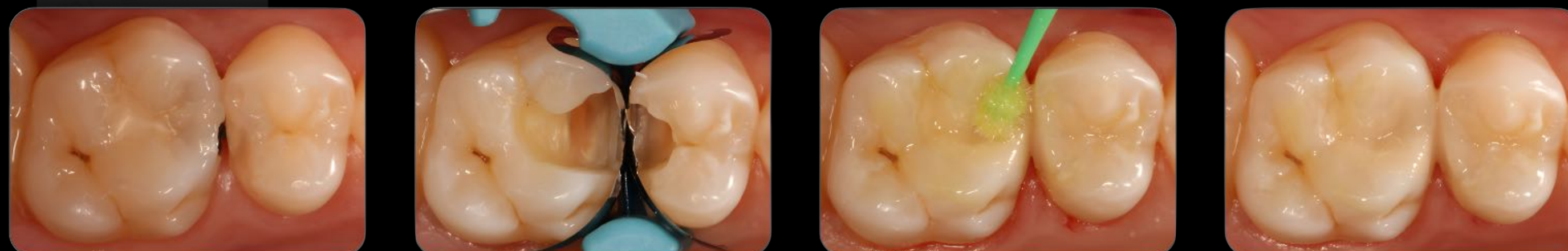

High caries risk patient - restored with Equia Forte HT and Equia Forte Coat®

Courtesy of Dr. Zeynep Bilge Kütük, Turkey

Features and Benefits


Nano-filled, self-adhesive resin coat

Single dispersion nano-filler and reinforced resin technologies

Improves strength and hardness of glass hybrid restoratives

Provides smooth glossy surface


Helps to maintain high wear resistance over time

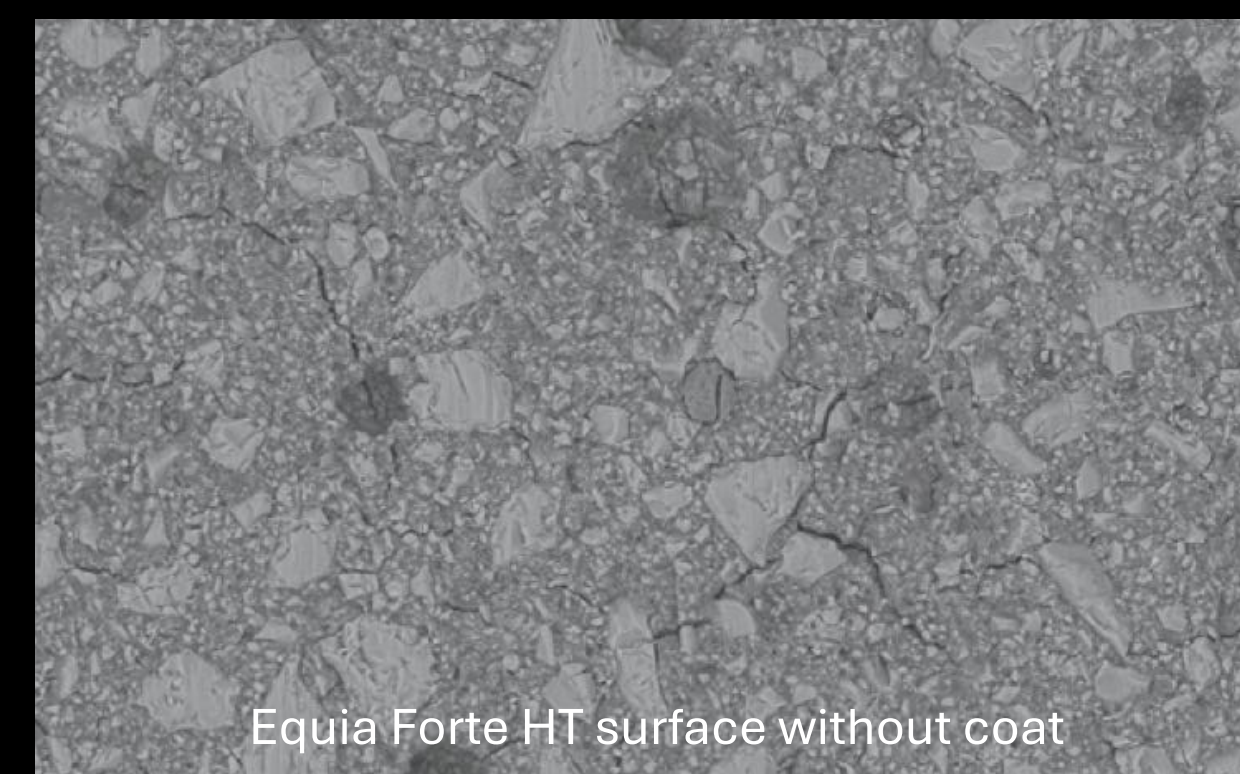

Equia Forte HT surface without coat




With Equia Forte Coat

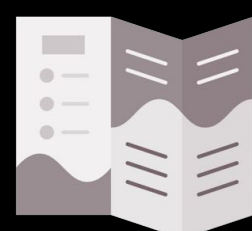


[View product brochure](#)



[View instructions for use](#)

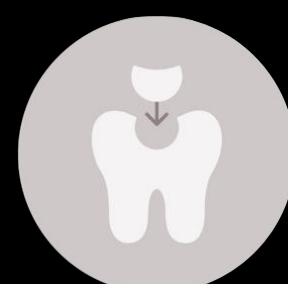

[View Technique Guide](#)

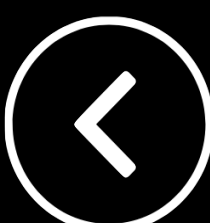

GC CONVENTIONAL & RESIN MODIFIED GLASS IONOMERS

Click a product image to view product information



DENTIN CONDITIONER



CAVITY CONDITIONER



EQUIA Forte Coat



COCOA BUTTER



GC Fuji II LC



GC Fuji VII




GC Fuji VII EP



Fuji IX - GP, FAST.
& EXTRA



GC Fuji BULK

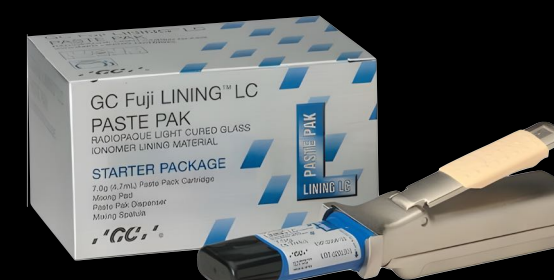

GC Fuji LINING LC



GC Fuji ORTHO BAND
PASTE PAK



View cement
usage guide

ALWAYS READ THE LABEL AND FOLLOW THE DIRECTIONS FOR USE



DENTIN CONDITIONER

Dentin Conditioner is a mild polyacrylic acid solution designed to remove the dentinal smear layer and to condition dentine and enamel before the application of your glass ionomer restorative.

It will increase the bond between your glass ionomer cement and tooth structure for added longevity.

The 10% polyacrylic acid solution cleans tooth surfaces to improve adhesion of glass ionomer cements

Applied on the tooth for 20 seconds

Step-by-step

1. Using a cotton pellet or sponge, apply dentin conditioner to the bonding surfaces for 20 seconds
2. Rinse thoroughly with water
3. Dry by blotting with a cotton pellet or gently blowing with an air syringe. DO NOT DESICCATE
4. Best results are obtained when prepared surfaces appear moist (glistening)
5. Apply glass ionomer as a bonding agent, base/liner or restorative material
6. Close the bottle immediately after use



Clinical Indications

Improves the marginal seal

Longevity of the restoration


Leaves smear plugs in the tubules

Less risk of post-operative sensitivity

Blue tinting




Dentine smear layer

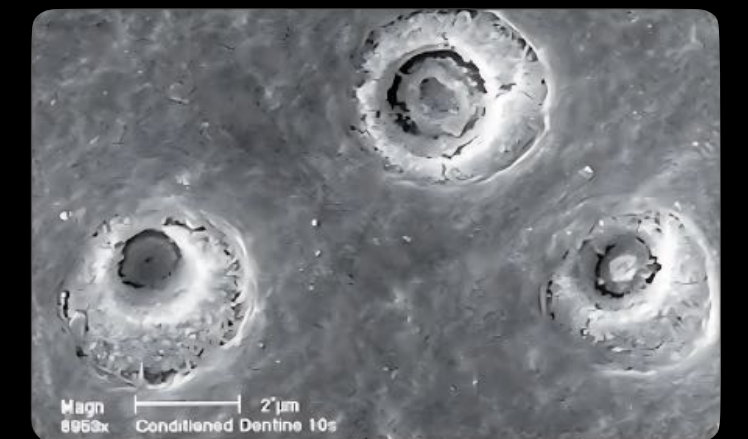

Conditioned dentine



[View instructions for use](#)

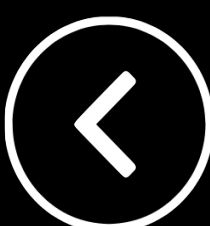

CAVITY CONDITIONER

Cavity Conditioner is a 20% polyacrylic acid solution designed to remove the dentinal smear layer and to condition dentine and enamel before the application of your glass ionomer restorative.

It will increase the bond between your glass ionomer cement and tooth structure for added longevity.

Applied on the tooth for 10 seconds

Clinical Indications

Cleans the surface effectively | Achieve optimal chemical adhesion

Creates mechanical retention

Less risk of post-operative sensitivity

Enhances the marginal seal

Blue tint




[View instructions for use](#)



Step-by-step


1. After tooth preparation, apply CAVITY CONDITIONER to the bonding surfaces for 10 seconds using a cotton pellet or sponge
2. Rinse thoroughly with water and dry. DO NOT DESICCATE
3. Apply GC Fuji BOND LC, GC Fuji II LC (improved) or other glass ionomer/glass hybrid as a bonding agent, base/liner or restorative material
4. Close the bottle immediately after use

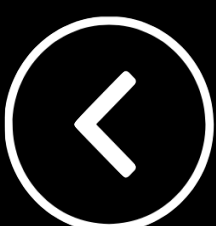



EQUIA Forte Coat®

EQUIA Forte® Coat incorporates a new multi-functional monomer that produces a tougher resin matrix for better protection. Used together with EQUIA Forte® HT Fil, the result is a strong bulk fill glass hybrid restorative system with exceptional physical properties and aesthetics.

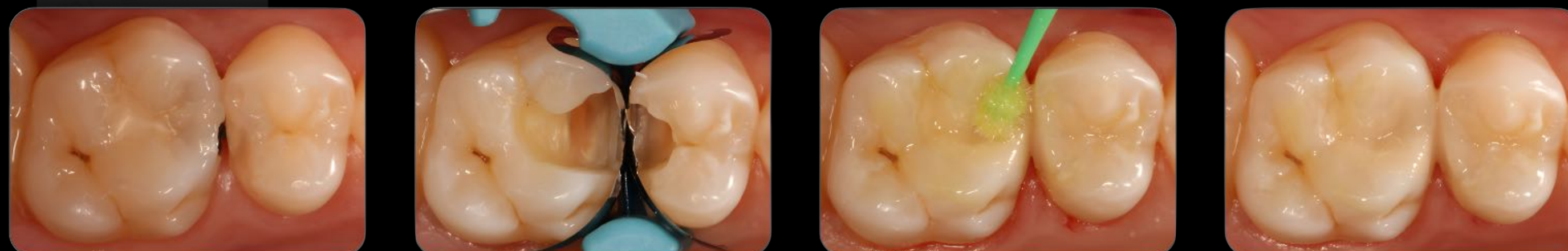

High caries risk patient - restored with Equia Forte HT and Equia Forte Coat®

Courtesy of Dr. Zeynep Bilge Kütük, Turkey

Features and Benefits


Nano-filled, self-adhesive resin coat

Single dispersion nano-filler and reinforced resin technologies

Improves strength and hardness of glass hybrid restoratives

Provides smooth glossy surface


Helps to maintain high wear resistance over time

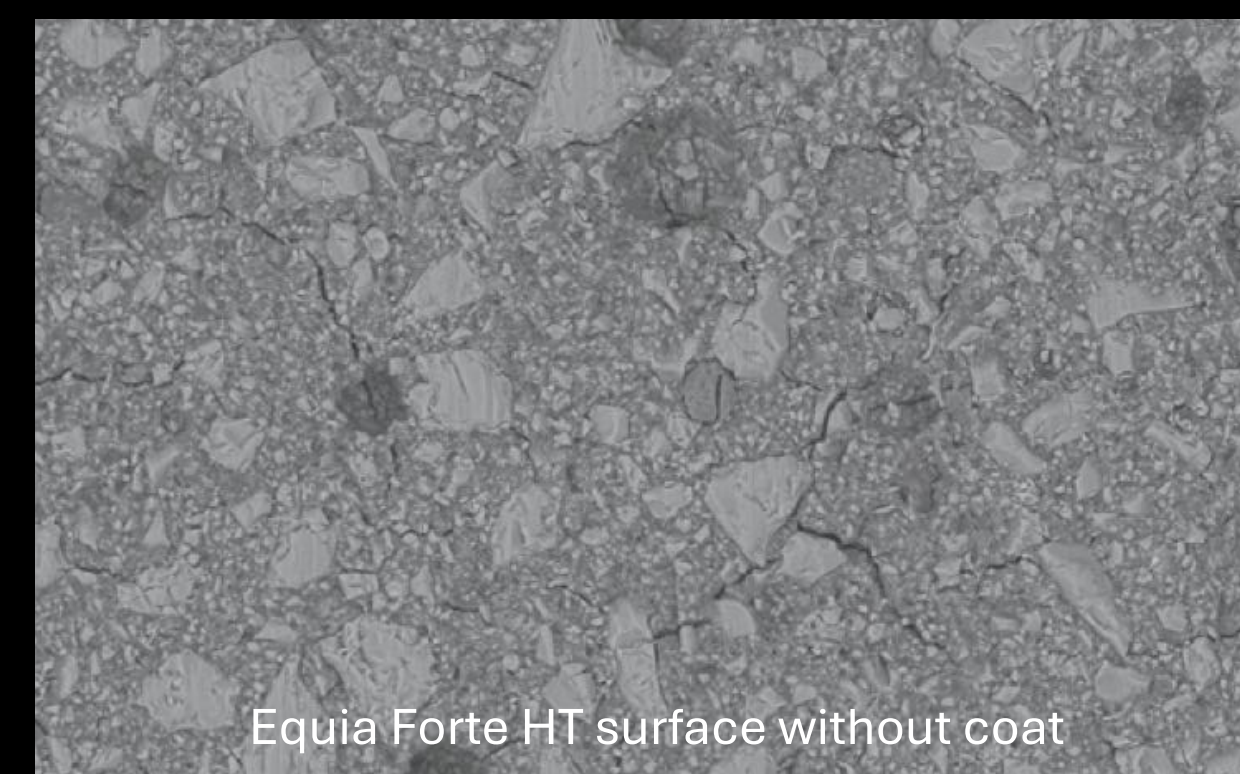

Equia Forte HT surface without coat



With Equia Forte Coat

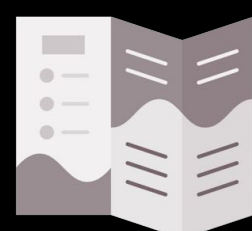


[View product brochure](#)



[View instructions for use](#)



[View Technique Guide](#)

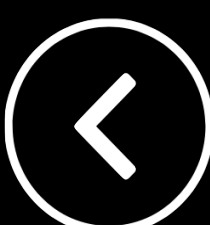

COCOA BUTTER

Indications

For the prevention of water and saliva contamination to exposed surfaces of glass ionomer cements during the first 24 hours after placement

For moisture proofing of luted lines

Clinical Features & Benefits

Excellent separation of GIC and metal instruments


Smooth filling contours

Prevents cement from sticking to the adjacent teeth

Easy to use



[View instructions for use](#)

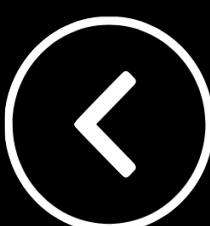

GC Fuji II LC



GC Fuji® II LC Capsule is the world's first resin-reinforced glass ionomer and has remained the benchmark for light cured glass ionomer cements. GC Fuji® II LC Capsule is a self-adhesive restorative material that is moisture tolerant, simple to use, quick to apply and is ready for final finishing immediately following 20 seconds light curing

Clinical Indications


- Class III and V restorations
- Root surface caries and cervical lesions
- Cervical lesions - Carious, Erosion, Abrasion & Abfraction
- Proximal lesions - primary and secondary teeth
- Reinforced base in sandwich technique
- Lesions where cavity margins are in dentine or
- When enamel cannot be adequately etched
- Base/liner under composite in the sandwich technique
- Core build up material
- Fissure sealant
- Adhesive for bonding orthodontic brackets

Clinical Features & Benefits

- Excellent aesthetics and polish-ability
- Abrasion resistance, ideal for cervical areas
- Good radiopacity
- Triple cure mechanism gives guaranteed cure for immediate finishing after light curing
- Available in 11 translucent Vita® shades
- Fluoride releasing and rechargeable uptake from external sources



(Dr. G. Knight)

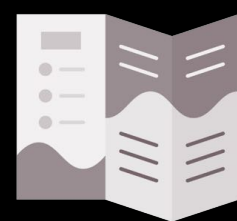

[View product brochure](#)



[View instructions for use](#)



[View capsule activation guide](#)




GC Fuji VII



Clinical Indications

- Surface and fissure protection
- Root surface protection
- Protection of erupting molars
- Hypersensitivity prevention and control
- Intermediate endodontic sealing
- Intermediate restorations
- Caries stabilisation and indirect pulp capping
- Cementation of stainless steel crowns
- Fissure protection

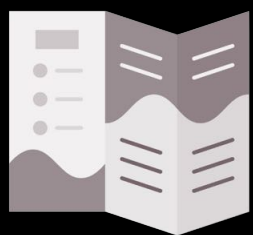

[View product brochure](#)



[View instructions for use](#)




[View capsule activation guide](#)


GC Fuji VII is the revolutionary glass ionomer sealant and surface protectant. It is fantastic at keeping young patients' teeth cavity-free. There is no isolation or bonding agent required. It works in a moist field. With GC Fuji VII, you will never have to worry about sealing over immature enamel or non-cavitated lesions. The self-bonding, high fluoride releasing GC Fuji VII creates a strong, acid-resistant fused layer.

Clinical Features & Benefits

- Flowable consistency to ensure effective wetting & intimate adhesion to tooth surfaces
- Moisture tolerant- can be applied in situations where saliva control is difficult
- Available in white and special pink shades for clear identification

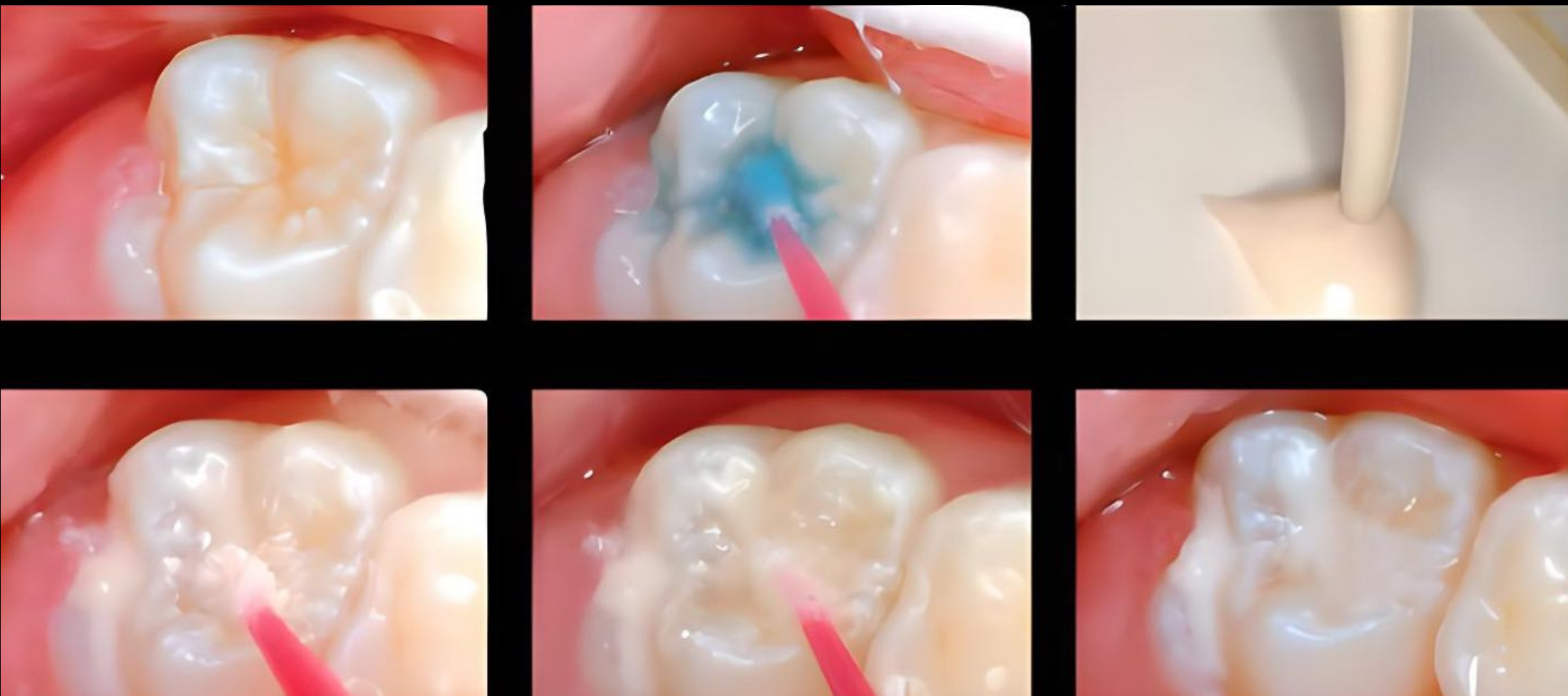

(Dr J. Lucas)

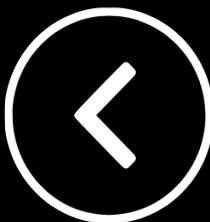


GC Fuji VII EP



Clinical Indications

- Higher risk patients
- Root surface protection
- Fissure protection
- Hypersensitivity prevention and control
- Protection for immature enamel (erupting molars)
- Temporary filling of endodontic access
- Intermediate restorations
- Protection of teeth with Molar Hypomineralisation

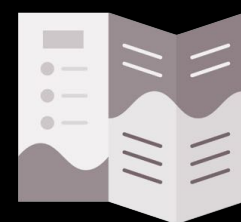

[View product brochure](#)



[View instructions for use](#)





[View capsule activation guide](#)

GC Fuji VII EP is a clinically proven glass ionomer cement formulation enhanced with RECALDENT® (CPP-ACP) for superior protection and remineralisation.

Once applied, GC Fuji VII EP forms a mechanical barrier to protect the covered and surrounding tooth surfaces and to minimise plaque accumulation. GC Fuji VII EP is essential when treating higher risk patients, including those with high acid exposure, reduced saliva flow, compromised cleaning and plaque removal capability.

Clinical Features & Benefits

- 3% RECALDENT® (CPP-ACP) enhances protection to adjacent dentine and enamel
- Remineralisation of tooth structure - prevents caries lesions from progressing
- Enhanced biofilm resistant surface
- Enhanced mineral release - calcium, phosphate, strontium and fluoride
- Superior access and delivery through GC Elongation Tips
- Available in pink and white shades

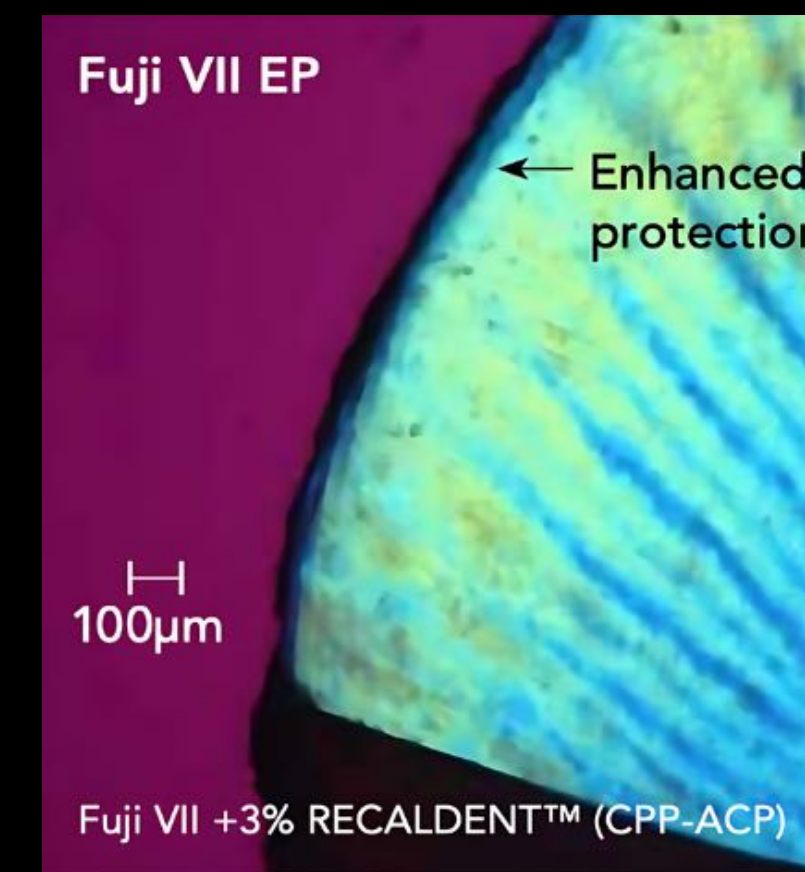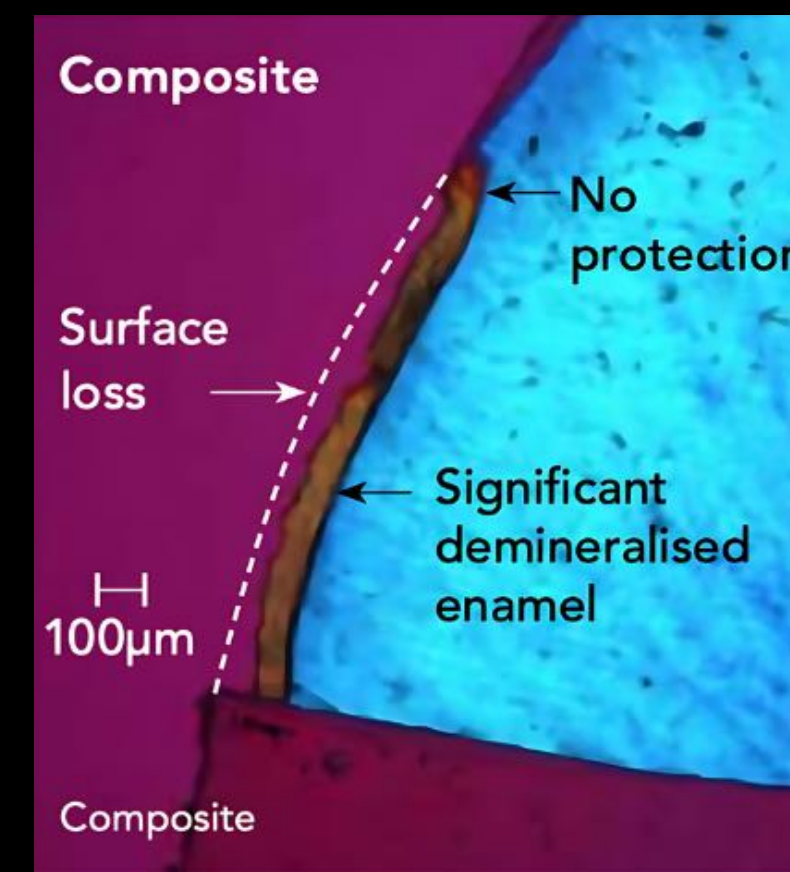

Polarised light micrographs of lactic acid (pH 4.8) demineralised lesions in enamel adjacent to Fuji VII EP and resin composite



Fuji IX

GP | GP FAST | GP EXTRA

Clinical Features & Benefits

- Chemical adhesion to tooth structure allows minimal cavity preparation
- No etching, no bonding, moisture tolerant
- High compressive and flexural strength for posterior indications
- Stronger for improved longevity and wear resistance
- Good handling characteristics, fluoride release
- Single step placement - no need for layering technique
- Good radiopacity

Clinical Indications

- Class I and II restorations in deciduous teeth
- Non-load bearing Class I and Class II restorations in permanent teeth
- Intermediate restorative and base material for heavy stress situation in Class I and Class II cavities using sandwich laminate technique
- Class V and root surface restorations
- Core build-up




Fuji® IX GP is one of the world's leading glass ionomer restoratives. A conventional, self-curing glass ionomer, it is easy to pack and contour in the cavity. If speed is of the essence, you will want to use

GC Fuji® IX GP EXTRA is the material of choice for better aesthetics. The special glass filler used give a very high translucency that changes little over time. The result is natural aesthetics never before achievable with conventional glass ionomer

GC Fuji® IX GP FAST enables placement of long lasting restorations in only half the time - only 3 minutes from start of mixing to finish

Fuji® IX GP EXTRA (Auto-cure GIC)	Manipulation time 2 min 30 sec	- Shades: A2, A3, A3.5, B1, B3, C4 - Excellent translucency
Fuji® IX GP FAST (Auto-cure GIC)	Manipulation time 3 min	Shades: A2, A3, A3.5, B2, B3, C4

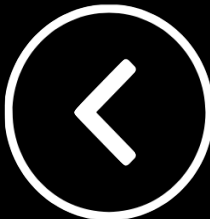

[View product brochure](#)

[View GIC usage guide](#)

[View capsule activation guide](#)



[View Fuji IX GP instructions for use](#)



[View Fuji IX GP Fast instructions for use](#)



[View Fuji IX GP Extra instructions for use](#)



GC Fuji BULK



Clinical Indications

- Class I and II restorations in deciduous teeth
- Non-load bearing Class I and Class II restorations in permanent teeth
- Base material for load bearing situations in Class I and Class II cavities using a sandwich laminate technique
- Transitional and intermediate restorations
- Class V and root surface restorations
- Core build-up


Perfect choice for

- Paediatrics & Geriatrics
- Dry mouth patients and patients with erosion
- Transitional restorations
- Caries-active patients and exposed root surfaces


Fuji® BULK is a robust, rapid-setting, multipurpose auto-cure glass ionomer cement (GIC) designed to meet the challenge of bulk placed restorations in more acidic oral environments. Truly remarkable in its formulation, Fuji BULK provides a unique balance between restoring function and protecting surrounding tooth surfaces from acid challenges.

Clinical Features & Benefits


- Designed for bulk placed restorations in more acidic oral environments
- Balance between restoring function and providing maximum protection to surrounding tooth surfaces from acidic challenges
- Fast setting auto-cure
- Seals and restores and is ready for final finishing just 2 minutes after start of mix

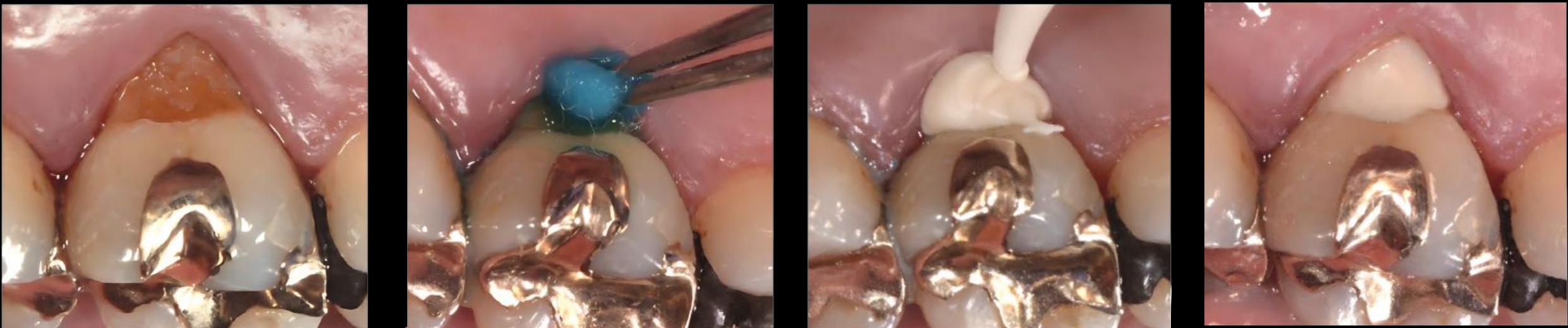

Dr E. Stephen Vouliotis

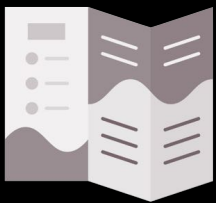

[View product brochure](#)

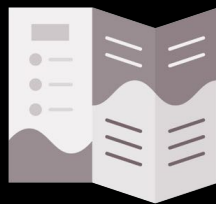


[View GIC usage guide](#)



[View instructions for use](#)



[View capsule activation guide](#)

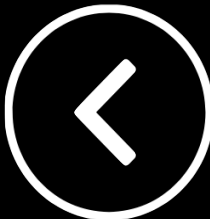

GC Fuji LINING LC




Clinical Indications

Light-cured glass ionomer cement
Liner or a base under composite or amalgam restorations

Clinical Features & Benefits

Easy adaptation to cavity walls without slumping
Excellent bonding to tooth structure offering a durable marginal seal
Controlled setting with light curing
Low water solubility
Insensitive to humidity
Radiopaque for easy post-operative diagnosis
Simple, accurate dispensing
Very easy bubble free mixing

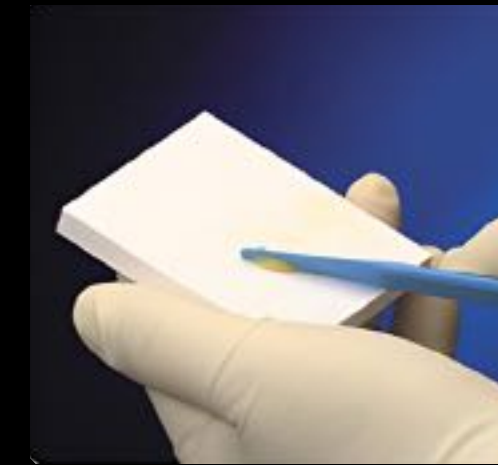

Easier Mixing

The paste-paste formulation always mixes to the ideal consistency because the exact proportions are dispensed each time. The result is a perfect mix in just 10 seconds, without any air bubble incorporation.



No dripping from the placement instrument

Once mixed, the lining cement is easily transferred from mixing pad to instrument where it holds without dripping, ready for placement exactly where you want it.



No slumping in the mouth

GC Fuji LINING LC has excellent flow, without running or slumping. You can easily control the even spread of cement over cavity walls and floor to form a strong, protective lining in even the shallowest of cavities.



Easy to identify after placement


The Universal shade balances the need for aesthetics, especially when lining under a shallow composite restoration, with the need for a visual contrast to help guide placement.



[View product brochure](#)



[View instructions for use](#)

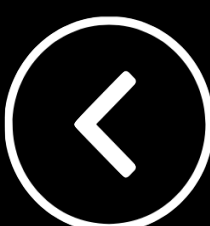

GC Fuji ORTHO BAND PASTE PAK



Indications

Cementation of orthodontic bands

Cementation of appliances retained with stainless steel crowns
(e.g., Herbst or expansion appliances)

Clinical Features & Benefits

Self-curing resin-reinforced glass ionomer for bonding of bands and appliances

Paste Pak cartridge delivery - mixed for 15sec for optimal consistency

Film thickness of just 3 μ m

Optimal handling and working time of 3mins 30sec

High fluoride release

Blue colour

Blue colour for easy identification during application, removal and after debonding

Fuji ORTHO BAND Paste Pak is a blue, resin modified glass ionomer for bonding orthodontic bands. It can also be used for the cementation of appliances retained with stainless steel crowns (ex: Herbst or expansion appliances).

The Paste Pak delivery saves time and reduces waste. It delivers consistent mixes every time, no matter the amount that the mix requires. Its blue colour distinguishes the cement for easy clean-up and removal.



[View product brochure](#)



[View instructions for use](#)



GC COMPOSITES & BONDS

Click a product image to view product information



G-Premio BOND



GC Modeling Liquid



Composite Manipulation Brush



Essentia Masking Liner



G-aenial Universal Injectable



G-aenial Anterior / Posterior



G-aenial Universal Flo



G-aenial Flo X



G-aenial A'Chord



Gradia Direct X



Essentia



Essentia Universal



everX Flow



EverStick C&B, Perio, Ortho



Epitex & GC New Metal Strips

ALWAYS READ THE LABEL AND FOLLOW THE DIRECTIONS FOR USE



G-Premio BOND

G-Premio BOND™ is a one-bottle universal bonding agent compatible with all etching modes. It can be used for direct bonding and repair.

G-Premio BOND™ offers the advantages of a **universal adhesive**, with great performance in all situations, thanks to a unique combination of three functional monomers (4-MET, MDP, MDTP). **G-Premio bond is HEMA-free**

With an easy application of G-Premio BOND™ achieves immediate chemical adhesion and high bond strengths that are ideal clinically.


Clinical Indications

- Direct bonding of light-cured composites and compomers to tooth structure
- Bonding of dual-cured core build up composites to tooth structure, as long as these materials are light-cured
- Intraoral repair of composite, metal-based and zirconia/alumina-based restorations
- Intraoral repair of ceramic restorations in combination with a silane coupling agent (like G-Multi PRIMER)
- Sealing of tooth preparation (cavity or abutment) for indirect restorations
- Cavities where onlays and inlays would also be indicated



Features & Benefits

- Fast and simple, no rubbing or complex application procedures
- Immediate chemical bonding, high bond strengths without needing to rub the surface.
- Chemically bonds within the first few seconds of application
- The drying step is highly reproducible: 5 seconds, MAXIMUM air pressure
- Fast 25 sec self-etch application and chemical adhesion
- Perfect for paediatric patients or where moisture control is challenging
- Yellow colour prior to light-curing gives visual control of application
- Visible during application, invisible after light curing
- Wets the tooth just like water
- Perfect margins, 3um film thickness - avoids bond pooling
- Eliminates patient discomfort

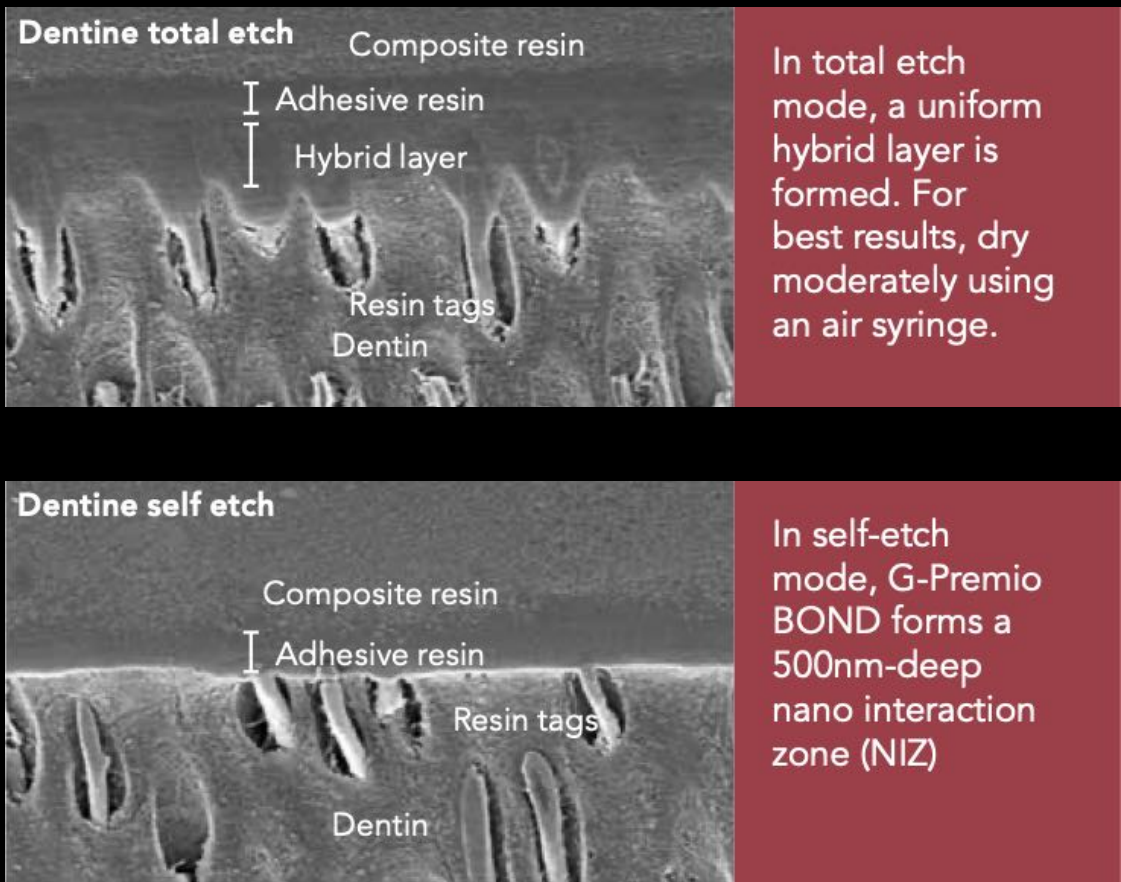

- 
View product brochure
- 
View instructions for use
- 
View Technique Guide
- 
View Repair Technique Guide
- 
View clinical article



Self-etch

Selective etch

Total etch



GC Modeling Liquid

Modeling Liquid is GC's solution for an easy, fast and efficient application of all direct composite materials. Applied on a brush and/or instrument, it eases the placement and shaping of composites and helps to achieve good morphology with a smooth finish.



Features & Benefits

- Can be used with an instrument or brush for quick and easy application and shaping of composites
- A smooth finish for a shorter finishing procedure
- Adapting the final layer of composite using a wetted brush aides in achieving a smoother surface
- Invisible for great aesthetic results
- GC Modeling Liquid becomes fully transparent after light-curing

Clinical Indications

- Modeling of composite materials for direct restorations
- Use with an instrument or the brush technique
- Pre-treatment of fibres in direct and indirect procedures



Dispense one drop of GC Modeling Liquid



Protect it from the light with a cover



Touch the drop with the brush or instrument in order to moisten it



Remove the excess thoroughly using a gauze swab



Use your instrument or brush to apply the composite



Use the flat brush to adapt and shape the composite



Dispense one drop of GC Model the dentin mamelons using a round brush



Shape and smoothen the final layer of composite using a wetted flat brush or instrument.



Apply the composite using your wetted instrument



The brush can be used to adapt the material cusp by cusp



GC Modeling Liquid and the brush aide to recreate the fissures and refine anatomy




View product brochure



View instructions for use



Watch How to use GC Modeling Liquid video

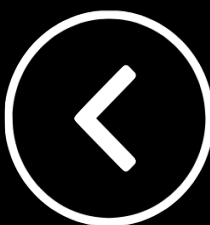

GC Composite Manipulation Brush

To aid and assist with the manipulation of composites. Ideal for composite stains and tints and use with GC Modeling Liquid



Dispense one drop of GC Modeling Liquid



Protect it from the light with a cover




Touch the drop with the brush or instrument in order to moisten it



Remove the excess thoroughly using a gauze swab



Use your instrument or brush to apply the composite

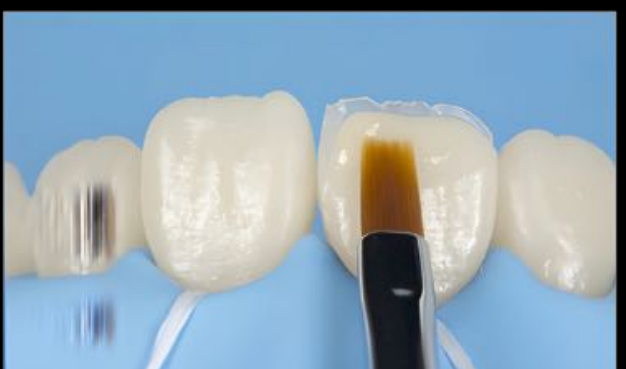

Use the flat brush to adapt and shape the composite



Dispense one drop of GC Model the dentin mamelons using a round brush



Shape and smoothen the final layer of composite using a wetted flat brush or instrument.



Apply the composite using your wetted instrument




The brush can be used to adapt the material cusp by cusp



GC Modeling Liquid and the brush aide to recreate the fissures and refine anatomy



Watch How to use GC Modeling Liquid video

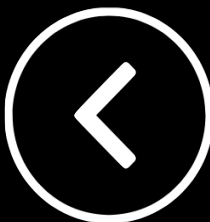

Essentia™ Masking Liner



Essentia Masking Liner is an opaque, heavily-filled injectable composite that will block discolouration without increasing value, creating the base for great aesthetics.

Clinical Indications

Suitable for anterior and posterior teeth

Covers amalgam staining

Blocks discolouration to improve aesthetic outcomes

Masking dark cores under veneers, crowns and direct restoratives

Wet and adapt well to the base of cavities and is well suited to placement in very thin layers




[View product brochure](#)



[View Essentia Clinical Guide](#)



Dr Chee Chang

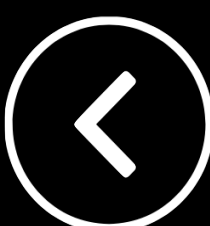


G-ænial Universal Injectable



Clinical Indications

- Direct restorative for Class I, II, III, IV and V cavities
- IMT (injection moulding technique)
- Fissure sealant
- Sealing hypersensitive areas
- Repair of direct and indirect aesthetic restorations
- Blocking out undercuts
- Liner or base
- Splinting of teeth in combination with fibres such as GC everStick fibres

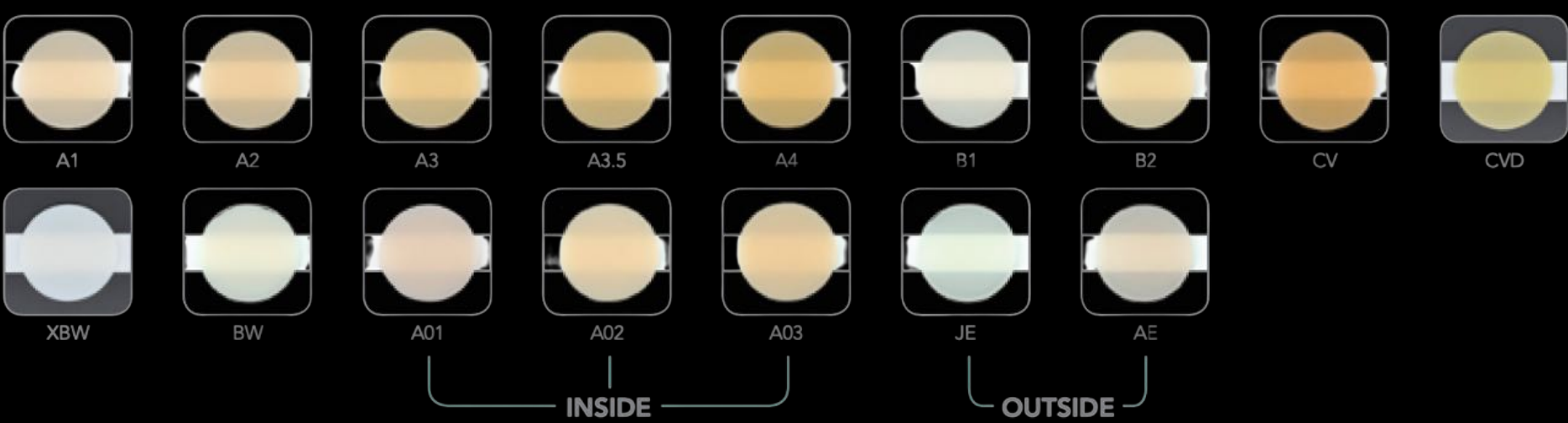


G-ænial™ Universal Injectable is a universal restorative composite displaying exceptional strength & resistance for long-lasting aesthetic restorations. It is incredibly strong - even stronger than leading paste composites on the market.

Features & Benefits

- Exceptional strength & wear resistance thanks to the combination of ultra-fine barium fillers
- GC's Full-coverage Silane Coating (FSC) technology
- Suitable for any cavity class without size limitation
- Strong enough without the need to cover with a conventional composite
- Unique thixotropic viscosity, optimal for free-hand build-up of cusps
- Shorter finishing steps as the anatomy is already created
- New syringe design featuring easy extrusion & no uncontrolled flow-out
- Optimal viscosity to easily separate the material from the tip and does not slump
- New tips with a long bendable needle for unparalleled access to difficult posterior cavities
- 16 shades in three levels of translucency, Excellent polish-ability and gloss retention
- High radiopacity (252%) - follow up your restorations and to easily detect any secondary caries



Injection Moulding Technique (Dr Anthony Tay)

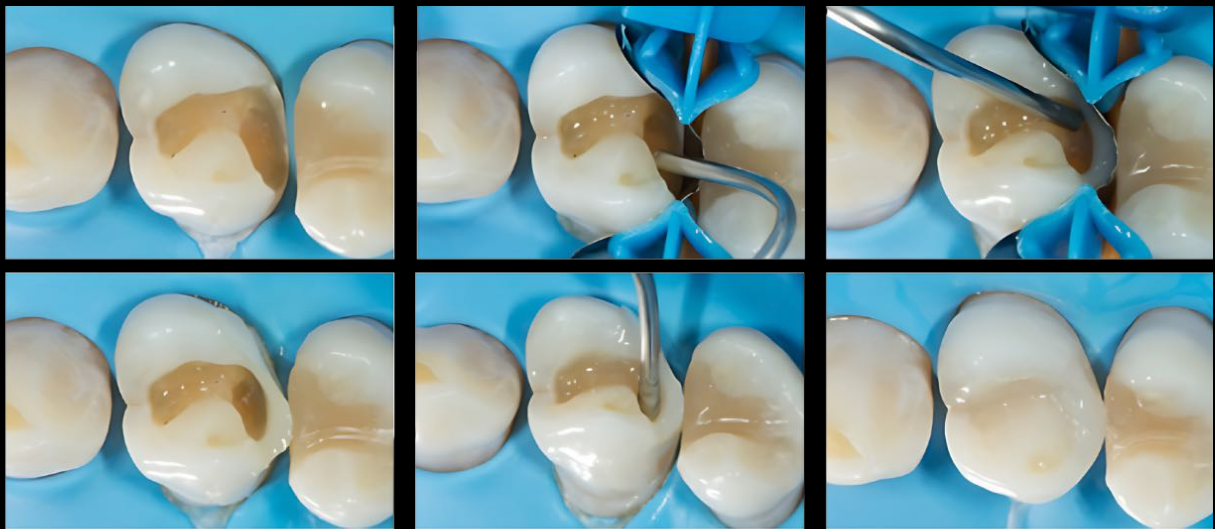

Posterior Restorations (Dr Javier Tapia Gaudix)



View Cusp Build-up Video



View Anterior Injection Moulding Video



View Posterior Injection Moulding Video




View product brochure



View instructions for use



View Injection Moulding Guide

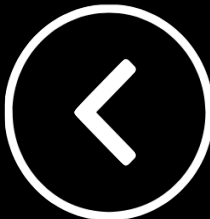

G-ænial Anterior & Posterior



Clinical Indications

G-ænial ANTERIOR

Direct restorative for Class II, IV, V cavities
For wedge-shaped defects and root surface cavities
For veneers and diastema closure

G-ænial POSTERIOR

Direct restorative for Class I and II cavities

Features & Benefits

- Natural reflectivity creates invisibility - excellent light scattering properties
- Forward-thinking shading concept
- High strength, low shrinkage, easy handling
- Anterior enables sculpting , Posterior firmer and more packable
- Single shade simplicity
- Good polishability and high gloss retention




[View product brochure](#)



[View instructions for use](#)

G-ænial Anterior and Posterior is the material of choice when you are looking for a restorative that will meet the aesthetic expectations of all your patients - in most cases with only one shade. Being a forgiving composite that offers beautiful and natural looking high gloss restorations whilst being extremely easy to use, makes it a composite suited to everyday use. G-ænial also gives you more working time so you are free to shape and sculpt it to obtain the ideal anatomical form and aesthetics with ease.

Anterior Shade Range

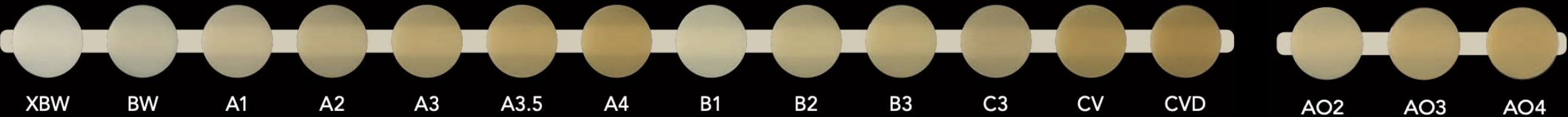

JE – Junior Enamel
for young patients

AE – Adult Enamel
for adult patients


SE – Senior Enamel
for senior patients

TE – Translucent Enamel
simulates a depth-effect within restorations

IE – Incisal Enamel
can be applied on the enamel occlusal third and proximal surfaces, particularly with young patients

CVE – Cervical Enamel
offers the right translucency to let the darker cervical dentine shine through

Posterior Shade Range

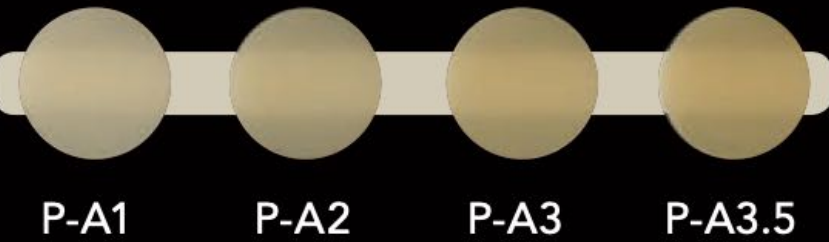


P-JE – Junior Enamel
for young patients

P-IE – Incisal Enamel
can be applied on the enamel occlusal third and proximal surfaces, particularly with young patients



Dr P Sheridan

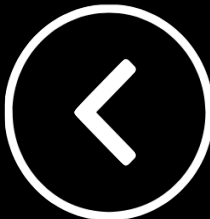

G-ænial Universal Flo




With GC G-ænial™ Universal Flo you get an exclusive blend of injectable viscosity for easy and convenient placement even in deep cavities, combined with outstanding physical properties for safe and durable restorations

G-ænial™ Universal Flo provides the answer to many difficult cases where you do not want to compromise between perfect placement, aesthetics and strength

Its unique combination of strength and flowability makes it also perfectly adapted to be used as a light-cured composite luting cement

Clinical Indications

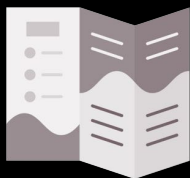

View product brochure



View instructions for use



View Technique Guide

True innovation in filler size and composition, in filler silanation and in manufacturing techniques has given rise to a remarkable composite resin. G-ænial Universal Flo is the first truly universal Injectable Composite, a material that clinicians can place with ease for beautifully adapted, beautifully finished, fast, aesthetic and tough composite restorations

Features & Benefits

Injectable composite, feels like a flowable, behaves like a conventional restorative

Superior polish retention like a nano-fill

15 shades based on Vita guide - 3 levels of translucency

Single shade or multi-shade layering application options




High strength, strong, polishable and high wear resistance

Fine 200nm filler size and composition

Full coverage silane technology



Dr S Koide

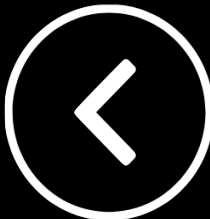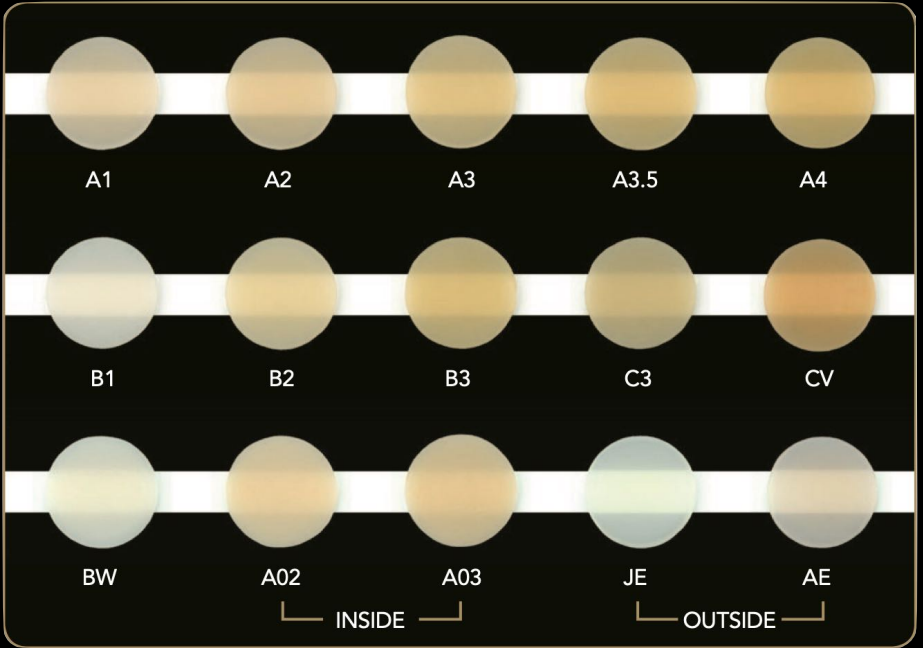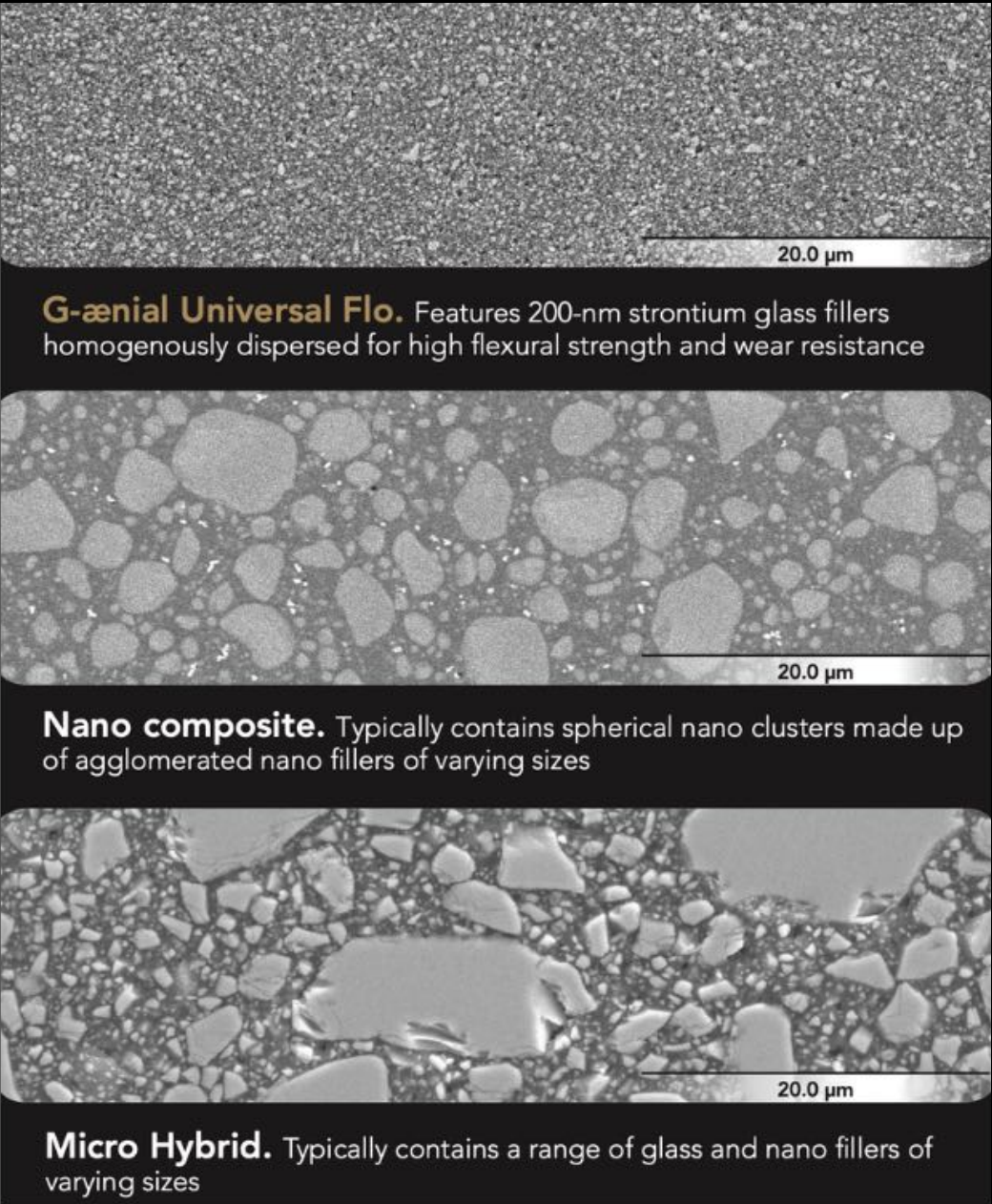

G-ænial Flo X



G-ænial™ Flo X is a high radiopacity and high strength flowable composite, designed specifically for optimum handling in lining, micro cavity and fissure sealing applications.

The best handling "conventional" flowable - yet with greater strength than most conventional composites.

Clinical Indications

Ideally suited for liner applications, sealant or where more flow is required

High radiopacity

Excellent wetting and adaptation following bond application

A full strength flowable composite for traditional flowable applications



Features & Benefits

A flowable with impressive strength



Innovation in filler technology gives you a high strength composite in a flowable consistency. Designed to be strong enough for application in any classification of cavity. G-ænial Flo X and G-ænial Universal Flo deliver similar high levels of strength while offering different viscosities tuned to the variation in composite handling preferences of individual clinicians.

High radiopacity

Perfect for defining the cavity outline in radiographic images for direct or indirect procedures. G-ænial Flo X shows one of the highest levels of radiopacity among flowable composites

Enhanced fluidity


Ideal handling for lining application G-ænial Flo X gives you total control over the flow. It's designed to wet and adapt yet maintain control for complete ease of application.

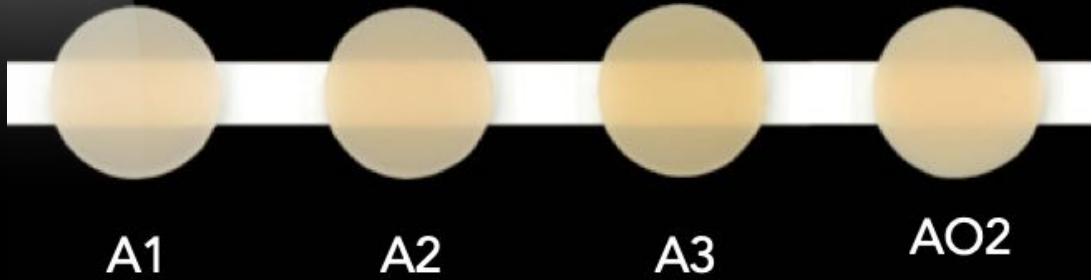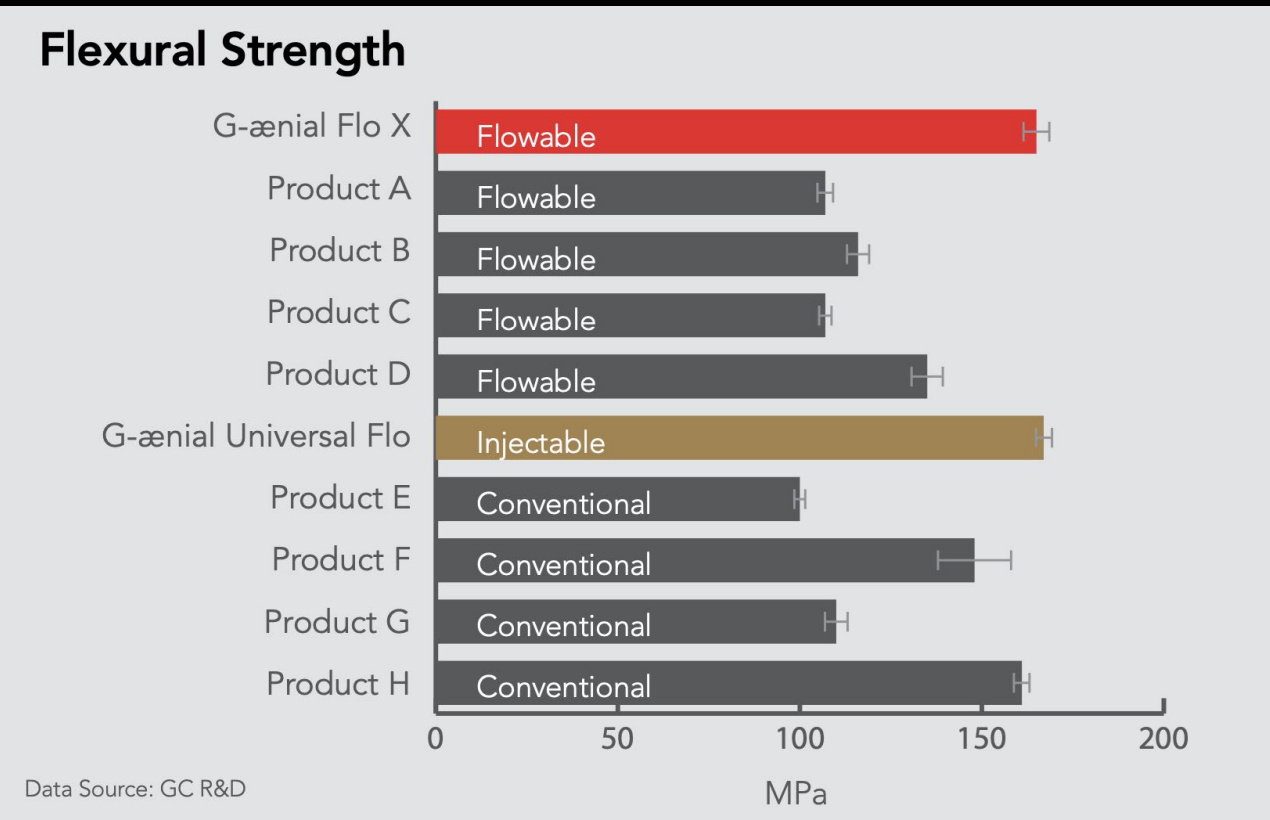

[View product brochure](#)



[View instructions for use](#)

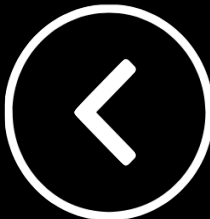

G-ænial A'Chord

G-ænial® A'CHORD is a universal composite restorative with a simplified unishade system that achieves aesthetics of 16 Classic VITA® shades in just 5 shades. Building on 10 years of clinical success, it balances simplicity, aesthetics, and advanced technology for aesthetics and strength for anterior and posterior restorations.

The composite's non-sticky consistency, achieved through a Bis-MEPP monomer and optimised filler-monomer, allows easy sculpting with hand instruments or brushes. GC's proprietary Full-Coverage Silane Coating (FSC) and High Performance Pulverised CERASMART (HPC) technologies enhance the material's strength, wear resistance, gloss retention, low discolouration, and excellent radiopacity.

Features and Benefits

- Natural fluorescence: beautiful invisible and naturally fluorescent end results under any lighting conditions (UV and even near-UV)
- Long-lasting gloss: exceptional polishing and high gloss retention thanks to specialised proprietary technology
- Extreme challenges: extended portfolio for extreme challenges: 2 cervical shades for elderly patients, 3 opaque for masking discolouration & preparation lines, 2 enamel for life-like translucency and 2 bleach shades
- High stain and war resistance: Specific mix of monomers ensure durable shade stability for long lasting beautiful restorations
- Excellent radiopacity: radiopacity (318% Al) optimises radiographic follow-up over time





Clinical Indications

- Direct restorative for Class I, II, III, IV and V cavities
- Direct restorative for wedge-shaped defects and root surface cavities
- Direct restorative for veneers and diastema closure



Dr A. Mak

Dr Y. Choi

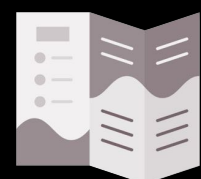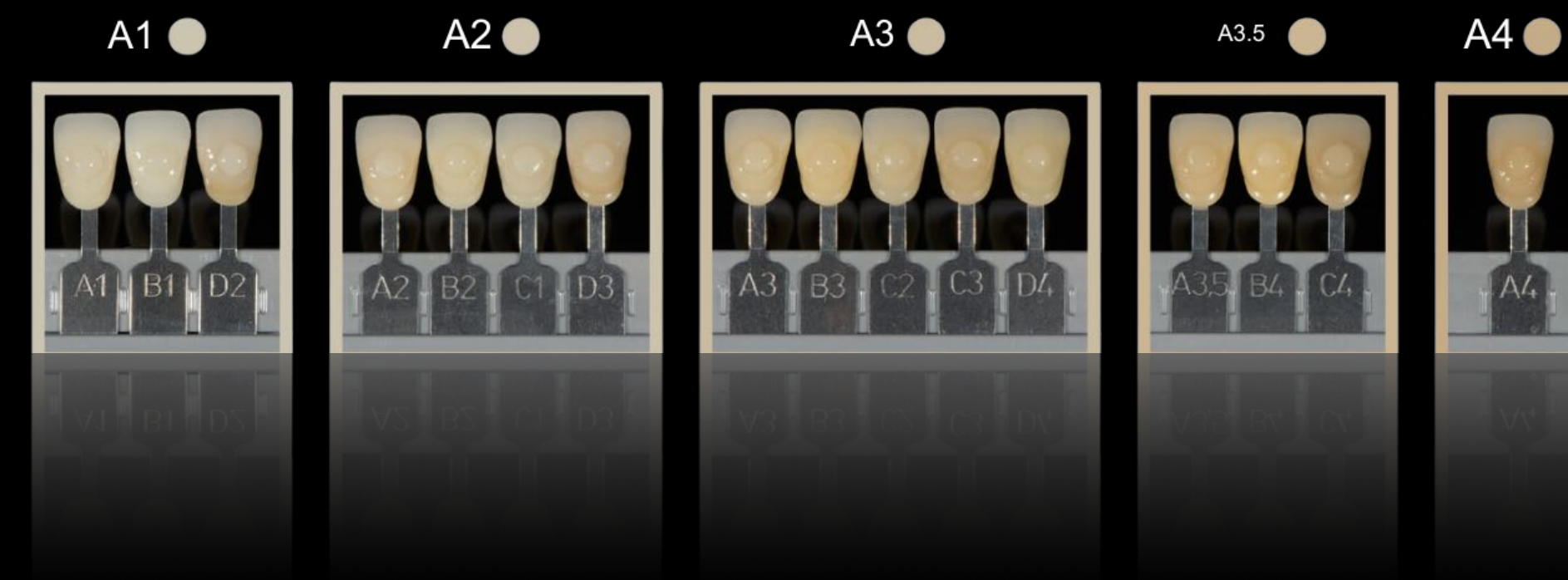

View product brochure



View Gaenial A'Chord Shade Guide



View Unitip instructions for use



View Syringe instructions for use



View Sticky or Not Video



View Gaenial A'Chord Polishing Video

Gradia® Direct X

GRADIA® DIRECT X has the same great handling you've come to love from GRADIA DIRECT with increased radiopacity.

GRADIA® DIRECT X simplifies your aesthetic restorations by combining smooth handling, durability and high polishability into a single, universal composite.

This product has been designed to provide excellent wear resistance and 'natural feel' in the mouth without wearing away opposing enamel.

Clinical Indications



- Direct restorative for Class I, II, I, IV, V cavities.
- Direct restorative for wedge-shaped defects and root surface cavities.
- Direct restorative for veneers and diastema closure.




Features and Benefits

- GRADIA DIRECT X achieves unmatched aesthetic invisibility through a remarkable chameleon effect
- 90% of restorations can be completed using a single shade
- Contains nano fillers, glass fillers and pre-polymerised fillers that imitate the complex light transmission, diffusion and reflection characteristics found in natural teeth
- Strong, radiopaque anterior and posterior composite
- Simple shade selection
- Optimal combination of handling and finishing characteristics
- A simplified shade range - majority of restorations can be completed using a single layer technique

A simplified shade range allowing the majority of restorations to be completed using a single layer technique.

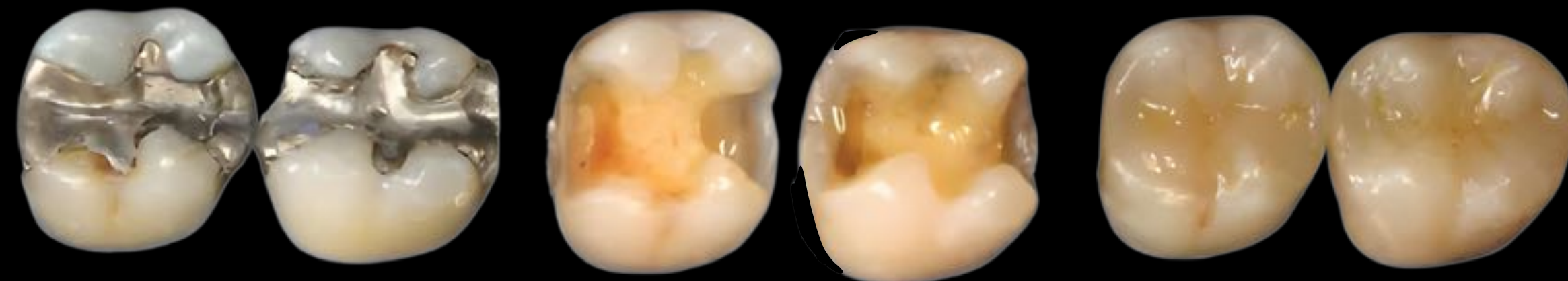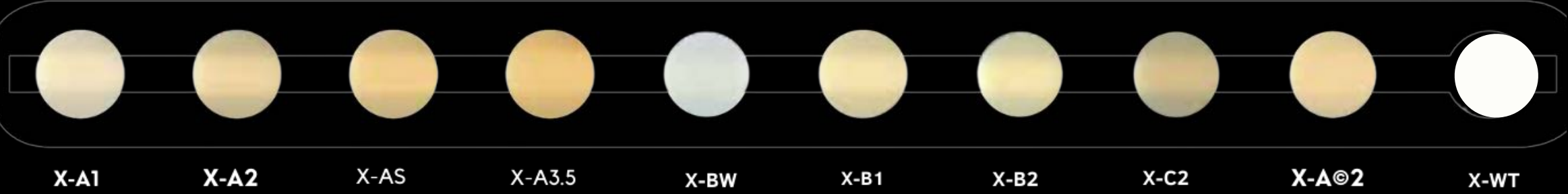

Dr L Harris

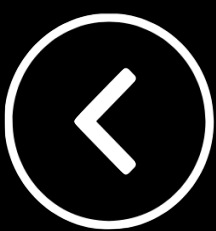

View product brochure



View composite brochure



View instructions for use



Essentia®

Essentia represents the pure essence of composite layering. With just seven syringes, you now have a perfect aesthetic solution for all your restorations. Traditional shade systems rely on different hues (colours), known as A, B, C and D. However, it is known that enamel and dentin have different characteristics which are also changing with time, becoming not only stronger in colour but also more translucent.

By breaking away from traditional monochromatic shades and focusing on reproducing the aging characteristics of enamel and dentin, you can obtain a much more natural result in just two simple layers.


Features and Benefits

- A simplified range of 7 shades - 3 dentine, 2 enamel and 1 universal shade
- High aesthetic results made simple with a bilaminar concept
- Mimic the natural ageing of teeth
- High gloss retention and wear resistance
- Ultra-fine particles used in the enamel shades
- Handling of each shade optimised to its use
- Characterisation options with Essentia™ modifiers for high-end aesthetics
- Effective masking liner to block discolouration without increasing value




Clinical Indications

- Direct restorative for class I, II, III, IV, V cavities
- Direct restorative for wedge-shaped defects and root surface cavities
- Direct restorative for veneers and diastema closure

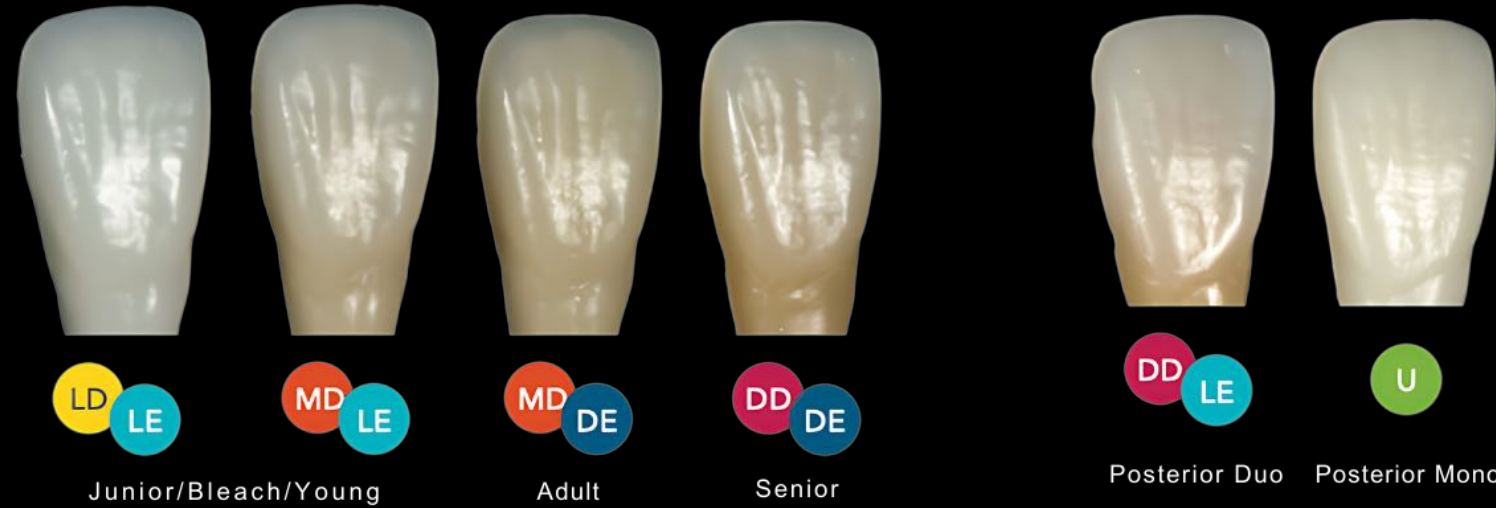

Dr Javier Tapia Guadix

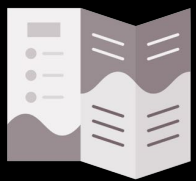

View product brochure



View instructions for use



View Clinical Guide



Essentia® Universal

Features and Benefits

The solution for all posterior restorations Essentia® Universal shade exhibits an excellent chameleon effect and seamlessly blends with the surrounding tooth structure, regardless of the shade of the cavity. Just one shade with amazing blendability

Available in three viscosities - all strong enough to be used in all cavity classes- Essentia® Universal shade is the ultimate solution for posterior restorations.

Essentia® Universal is a full-strength paste composite with a simplified one shade concept. Its Universal shade displays great blending effect together with a packable consistency for easy placement and shaping.

Essentia® HiFlo displays high wettability and adaptation, thanks to its fluid viscosity. Based on 700m barium glass fillers, this composite features a high radiopacity of over 300% and excellent physical properties, making it possible to be used in any cavity class.

Essentia® LoFlo is a full-strength injectable composite with low flowability, optimal when a thixotropic material is needed. Based on 200nm ultra-fine strontium fillers, this composite features excellent physical properties, making it possible to be used in any cavity class.



Clinical Indications

Direct restorative for Class 1, 11, III, IV and V cavities

Liner or base | Blocking out undercuts | Fissure sealant

Sealing hypersensitive areas

Repair of direct and indirect aesthetic restorations, temp crown & bridges and defective margins when margins are in enamel

Dr Bojidar Kafelov



Essentia® Universal



Prof Joseph Sabbagh




Essentia® HiFlo



Dr Bojidar Kafelov



Essentia® LoFlo

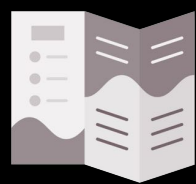


[View product brochure](#)



[View instructions for use](#)



[View HIFLO LOWFLO instructions for use](#)

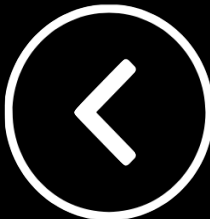

everX Flow




Clinical Indications

- Dentin replacement material for all direct composite restorations
- Large posterior cavities, deep cavities and endo-treated teeth, cavities with missing cusps or after amalgam removal
- Alternative to cavities where inlays and onlays would also be recommended
- Core build-up - crowns, inlays and onlays

Features & Benefits

- Optimal to strengthen large restorations thanks to the reinforcing effect of fibres
- Excellent fracture toughness reducing the risk of catastrophic failures
- Perfect thixotropy for easy placement and adaptation
- No slumping even in upper molars
- Excellent aesthetic results with the Dentin shade (in layers of 2mm)
- Indicated both for dentin replacement and for core build-up preparations

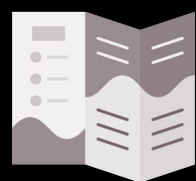

View Product Brochure



View instructions for use



View Technique Guide



View Garoushi Paper



View Magne Paper




View Ozsevik Paper

everX Flow™ is a short-fibre reinforced flowable composite indicated for dentine replacement in direct restorations (together with a conventional composite as enamel layer), and for core build-up.


Thanks to the short fibres it contains, ever Flow efficiently reinforces restorations and displays an exceptionally high fracture toughness. Fibres also help to redirect cracks and avoid catastrophic failures, which makes ever Flow an optimal material to use in weakened or cracked teeth, for instance after amalgam removal.

everX Flow has a very thixotropic viscosity which allows it to flow and adapt perfectly to the cavity floor, but without any slumping - even when used in upper molars. This optimal consistency makes your restorative treatment both faster and easier. everX Flow is available in two unique shades to answer all your clinical needs.


The Bulk shade displays a depth of cure of 5.5mm and is perfect for deep cavities, or whenever you want to speed up the treatment. The Dentine shade has a higher opacity and requires to be layered, and will be the best choice when you are looking for the best possible aesthetic results.

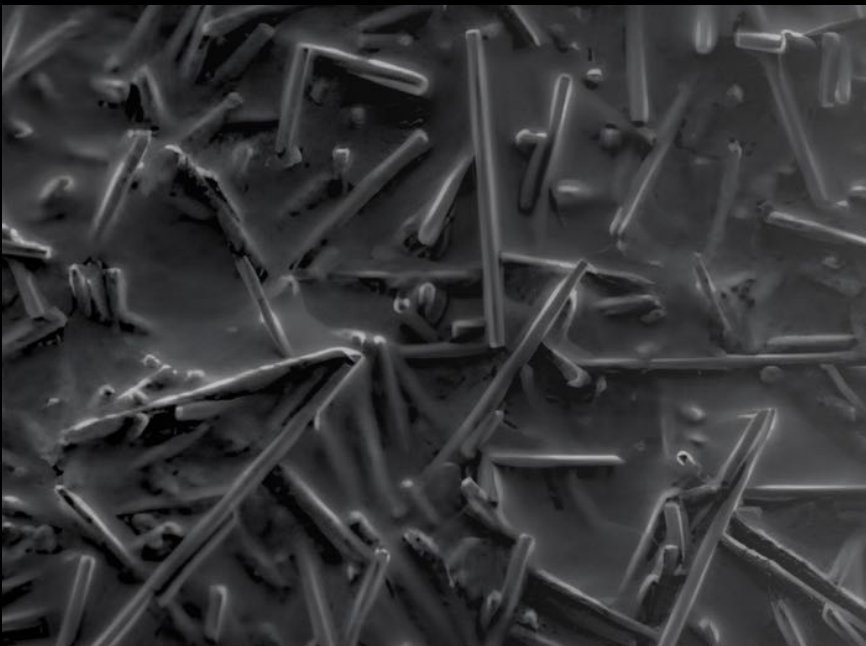

E-glass fibres in everX Flo

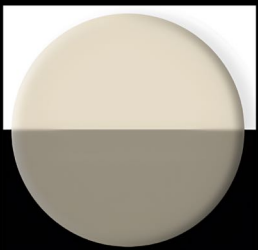

Bulk Shade for Deep Cavities and Cores - 5.5mm depth of cure

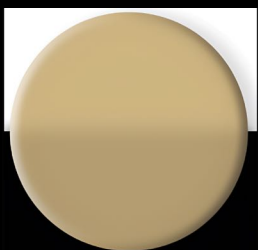

Dentine Shade for cores and aesthetics - 2.0mm depth of cure




Posterior restoration with everX fibre

Courtesy: Dr. Yo Han Choi



Watch On-Demand Webinar

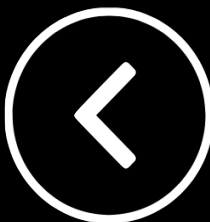

everStick



Clinical Indications

everStick®C&B

- Surface retained bridges
- Inlay and onlay bridges
- Hybrid bridges
- Full cover crown bridges
- Implant-supported bridges
- Temporary bridges
- Immediate bridges

everStick®NET

- Labial splinting of traumatised teeth
- Labial periodontal splinting
- Repair and reinforcing of veneers

everStick®ORTHO

- Aesthetic retention phase after active orthodontic treatment

everStick®PERIO

- Surface-retained or intra-coronal splinting
- Combined periodontal splints and surface-retained bridges

everStick® fibre reinforcements are made of silanated glass fibres in thermoplastic polymer and light curing resin matrix. everStick products address the advantages of minimally invasive dentistry where the patient's own healthy tooth tissue is saved for as long as clinically possible. This also means that other treatment options remain available should the patient ever need them in the future.

Features & Benefits

- Minimally invasive and reversible; leave other treatment options available
- Superior mechanical properties
- Unique patented bonding
- Elasticity close to that of dentin
- Extensive research data



View product brochure



View Ortho brochure



Instructions for use & Technique Guides



everStick C&B



everStick Perio




everStick Post



everStick Ortho




everStick NET

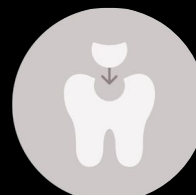


everStick C&B
Technique Guide



everStick Perio
Technique Guide

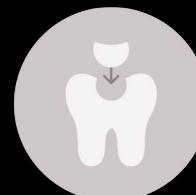

everStick Post
Technique Guide

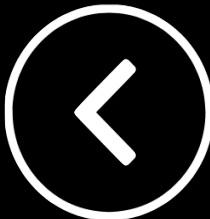

Epitex



Clinical Indications

Epitex finishing and polishing strips are for contouring, finishing and polishing of composite, glass ionomers and metal restorations. Epitex finishing and polishing strips are ultra-thin because the abrasive particles are not bonded to the strip by an adhesive. Thanks to this feature, Epitex allows easier access in tight contact points, helping to minimise gingival damage.

Features and Benefits

- Extremely thin and flexible material
- Excellent adaptability even between tight interproximal contacts
- Very tear resistant
- Abrasive particles are firmly embedded into the strip surface
- Handy strip dispenser
- Four grades:
 - COARSE - preparation
 - MEDIUM - final preparation
 - FINE - finishing
 - EXTRA FINE - polishing
- Available in 10m length individual pack or as a starter kit

GC New Metal Strips



Clinical Indications



Colour-coded polishing strips for preparation, contouring, finishing and polishing of composite, glass ionomer or metal restorations.

New Metal Strips are electro-plated with aluminium-oxide particles.

The strips will not stretch or break, and the excellent wear resistance and superior abrasive bonding allow wet finishing, which assures long-lasting and repeated use.

Features and Benefits

- Electroplated aluminium oxide particles
- Very flexible and tear resistant
- Autoclavable and easy to clean
- Able to reach between tight interproximal contacts
- Wide (4mm) and Narrow (2.6mm & 3.3mm)
- Individual packs of 12 or assorted pack of 12

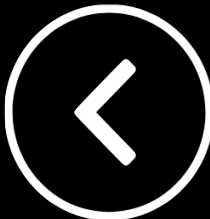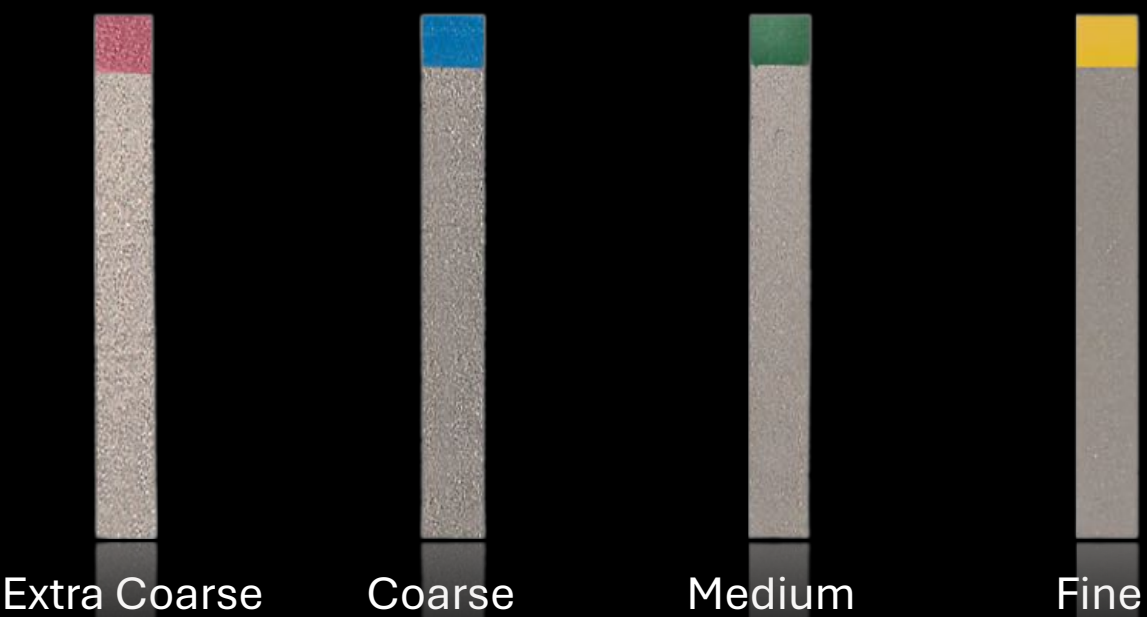

GC LUTING CEMENTS

Click a product image to view product information



G-CEM ONE



G-CEM Capsules



G-CEM Linkforce



G-Multi-PRIMER




G-CEM Veneer



GC Fuji PLUS CAPSULE



Click to download GC Luting Guide App

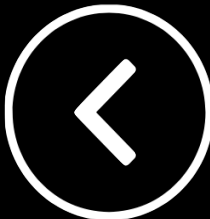

ALWAYS READ THE LABEL AND FOLLOW THE DIRECTIONS FOR USE



G-CEM ONE™

Clinical Indications

- Cementation of all types of all ceramic inlays, onlays, crowns & bridges
- Cementation of resin and metal-based inlays, onlays, crowns & bridges
- Cementation of metal, ceramic, fibre posts, and cast post and cores
- Cementation of all ceramic and composite veneers
- Final cementation of crowns and bridges on implant abutments

Features & Benefits



- Excellent bond strength to enamel, dentin and to all indirect materials
- Optimal bond strength for retentive and non-retentive preparations
- Optional Adhesion Enhancing Primer
- Exceptional self-curing for your peace of mind when thick or opaque restorations are used
- Tack cure feature for easy excess removal
- Increased moisture and saliva tolerance when Adhesive Enhancing Primer is used
- Optimised flow of the paste, virtually no post-operative sensitivity
- Invisible, wear resistant margin for aesthetic outcomes
- 4 stable shades: A2, Translucent, AO3, BO1 (White Opaque)



G-CEM ONE™ is a truly universal, non-technique sensitive, versatile and reliable product that gives the flexibility of being effective in all cementation procedures for any type of restorations; from metal-based to resin and all-ceramic inlays, onlays, crowns, bridges, and posts. It demonstrates excellent bond strength to enamel, dentine and all indirect restorations.

In addition, the optional Adhesive Enhancing Primer (AEP) and its innovative GC Touch Cure feature provides optimal bond strength in challenging clinical situations, such as low retentive surfaces. The chemical initiator in the AEP accelerates the chemical cure of the cement to ensure optimal bonding. The AEP is easy to apply; no additional light-curing is needed - apply and air-dry.


G-CEM ONE can be combined with the dedicated Adhesive Enhancing Primer (AEP). This innovative Primer induces what we call "touch cure mechanism" for an accelerated and very effective dark-cure (self-cure) reaction of G-CEM ONE.

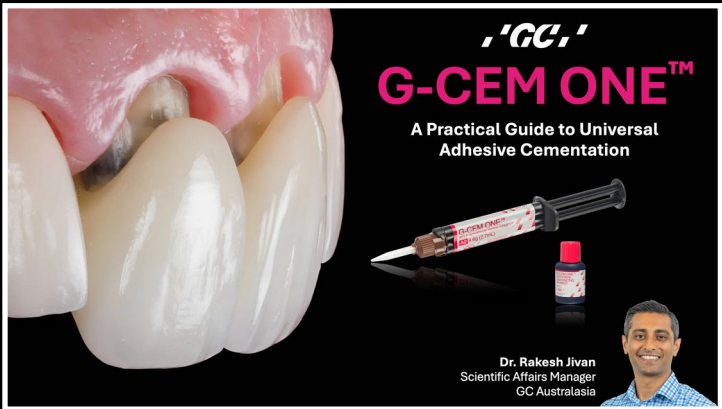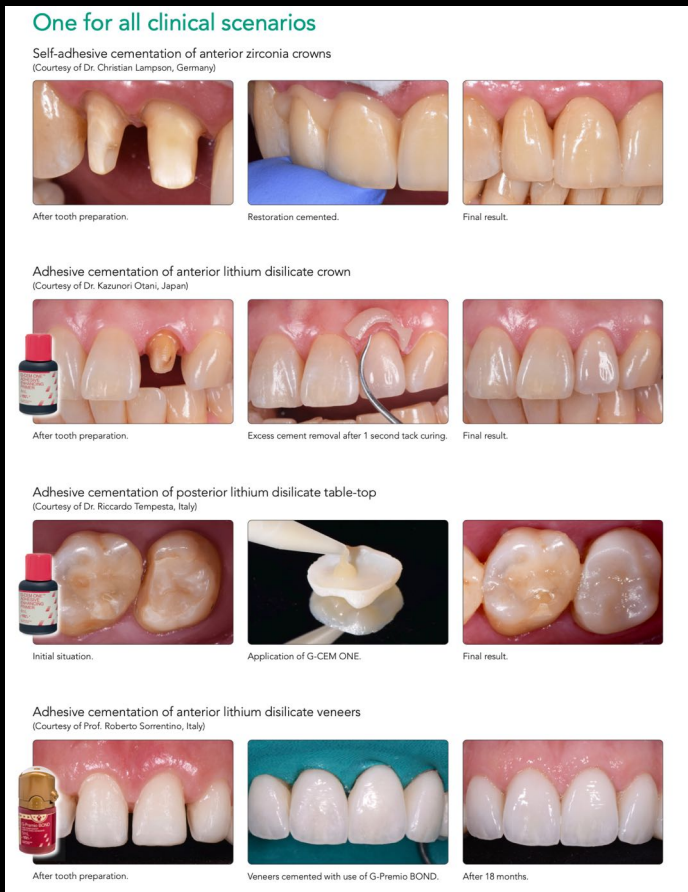

View On-Demand G-CEM One Clinical Presentation (1 CPD)



View On-Demand Zirconia Bonding Webinar (1 CPD)

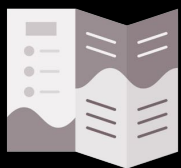

View product brochure



View Luting Guide brochure



View instructions for use



View Technique Guide



View Dr Miles Cone Video Part 1



View Dr Miles Cone Video Part 2



View Dr Miles Cone Video Part 3



View Dr Miles Cone Video Part 4




View Dental Advisor Video



View Simply The One Video



View Luting Non-Retentive Lithium Disilicate Cementation

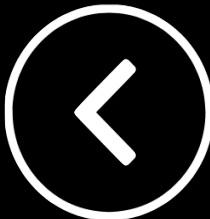

G-CEM Capsule




Clinical Indications


Cementation of all types of all ceramic, resin and metal-based inlays, onlays, crowns and bridges (max. 2 abutments)

Features & Benefits

- Excellent bonding property to restorative materials, including Zirconium
- High bond-strength to enamel in both self-cure and light-cure modes
- No need for pre-conditioning on the tooth structure
- Direct delivery into the root canal with special elongation tip
- Long-term stability and colour stability since there is no dissociation after curing
- Lowest linear expansion

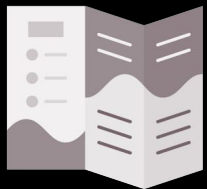

View product brochure

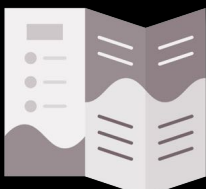

View Luting Guide brochure



View instructions for use



View capsule activation guide


G-CEM is a dual-cure self-adhesive universal resin cement delivered in capsules, designed for the adhesive luting of all-ceramic, zirconia, alumina, metal or composite indirect restorations.

G-CEM combines the improved handling and self-adhesion of conventional cements with the superior mechanical properties, adhesion and aesthetics of resin cements.



Light cure for 2-4 seconds until the cement has reached a gel phase. The excess cement is then easily removed using a scaler tip.

Cementation indications	Surface pre-treatment†		G-CEM
- Inlays, Onlays	Sandblast or etch	Ceramic Primer	Self-adhesive luting cement
Metal, PFM			
- Inlays, Onlays	✓		Y
- Crowns	✓		Y
- Bridges	✓		Y
- Endodontic Posts	✓		Y
Ceramics (less than 600 MPa flexural strength)			
- Inlays, Onlays	✓	✓	Y
- Crowns	✓	✓	Y
- Bridges	✓	✓	Y
High Strength Ceramics (greater than 600 MPa flexural strength)			
- Crowns	✓		Y
- Bridges	✓		Y
- Endodontic Posts	✓		Y

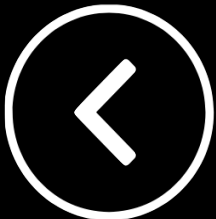

G-CEM LinkForce




Clinical Indications


Inlays, Onlays, Veneers, CAD/CAM, Crowns, Bridges

From inlays/onlays, overlays, veneers and tabletops to CAD/CAM prosthetics, the choice between different types of restorations has become wider. Add the introduction of new substrates such as zirconia, lithium disilicate and hybrid ceramics and it is no surprise that mastering all the cementation procedures for all indications is now a difficult challenge. Is there one solution that can adapt to all situations?


GC has developed a universal adhesive resin cement for you that offers predictable results and allows you to work in a standardised way. G-CEM LinkForce® from GC - the universal and powerful solution to all your adhesive cementation challenges that ensures a strong bond to all substrates in all indications, without compromises.

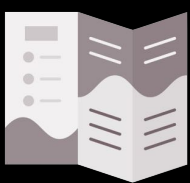

View product
brochure

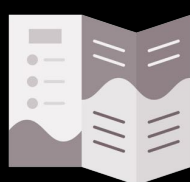

View Adhesive Resin
Cement brochure

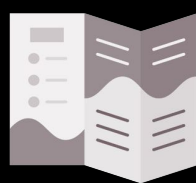

View Luting Guide
brochure



View G-Cem Linkforce
technique guide



View G-Cem Try-in Paste
instructions for use



View G-Cem LinkForce
instructions for use

Features & Benefits

Secured adhesion in all situations with only one system with three base elements

G-Premio BOND bonds to ALL preparations with no compromises: bonds to teeth, metal abutments and composite core build-ups

G-Multi Primer ensures a stable adhesion to ALL restorations: stable chemical adhesion on all surfaces, even on precious metals

G-CEM LinkForce provides a strong link in ALL indications: universal without limitations in indications or substrate

Light-cure, dual-cure: the choice is yours

High bond strength thanks to an efficient light-curing of the bond: very thin 3um film thickness does not interfere with crown placement


Efficient self-cure mode: particularly useful when luting opaque or thick restorations

Optimal light-curing of the resin cement: ideal for luting veneers

Aesthetic & stable over time

Shades to match all needs, accompanied by their corresponding try-in pastes

Tooth-like fluorescence, colour stability over time

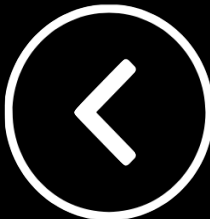

G-Multi PRIMER



Clinical Indications

Priming of the adherent surface of indirect restorations made of glass ceramics, hybrid ceramics (e.g. CERASMART), zirconia, alumina, composites and metal before luting procedures


Intra-oral and extra-oral repair of fractured indirect restorations made of ceramics, hybrid ceramics, zirconia, alumina, and composites

Features & Benefits

- Silane will provide adhesion to glass ceramics, Lithium Disilicate, Hybrid, ceramics and composites
- MDP ensures adhesion to Zirconia and non-precious metal alloys
- MDTP will ensure adhesion to precious metal alloys
- No refrigeration: 2-year shelf life

G-Multi PRIMER™ is a universal primer designed to pre-treat virtually any substrate for indirect restorations. G-Multi PRIMER™ uses three different chemical agents to ensure perfect adhesion in all situations. By adding silane to a universal primer rather than a bonding agent, the stability of the adhesion is assured.

Substrate	Pre-treatment at lab or chairside	Step 1	Step 2	Step 3	Step 4	Step 5 In case contamination with saliva/ blood before primer application	Step 6 Primer Application	Remark
Fedspathic ceramics, Leucite-reinforced ceramics	Etching at chairside	Try-in	Rinse thoroughly & dry	Etch inner surface for 60 secs with 5% HF acid	Rinse & dry	Clean with alcohol & dry	G-Multi PRIMER Apply & dry	Adhesion of G-Multi PRIMER through silane
Lithium disilicate	Etching at chairside	Try-in	Rinse thoroughly & dry	Etch inner surface for 20 secs with 5% HF acid	Rinse & dry	Clean with alcohol & dry	G-Multi PRIMER Apply & dry	
Zirconia Alumina	Sandblasting at chairside	Try-in	Rinse thoroughly & dry	Sandblast	Rinse & dry	New sandblasting or clean with Ivoclean*	G-Multi PRIMER Apply & dry	Do not clean the Zr oxide surfaces with phosphoric acid.
	Already sandblasted by lab	Try-in	Rinse thoroughly & dry			New sandblasting or clean with Ivoclean*	G-Multi PRIMER Apply & dry	Adhesion of G-Multi PRIMER through MDP
Metal Composite Hybrid Ceramics	Sandblasting at chairside	Try-in	Rinse thoroughly & dry	Sandblast**	Rinse & dry	Clean with alcohol & dry	G-Multi PRIMER Apply & dry	Adhesion of G-Multi PRIMER through silane (to glass fillers), MDP (for non-precious metal, resins) and MDTP (for precious metal)
	Already sandblasted by lab	Try-in	Rinse thoroughly & dry			Clean with alcohol & dry	G-Multi PRIMER Apply & dry	
Fibre Post	At chairside	Try-in	Rinse thoroughly & dry			Clean with alcohol & dry	G-Multi PRIMER Apply & dry	

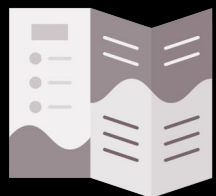

View product brochure



View instructions for use



View Technique Guide



Silane


Adhesion to glass ceramics, hybrid ceramics and composites

MDP

Adhesion to zirconia, alumina and non-precious metals

MDTP

Adhesion to precious metals

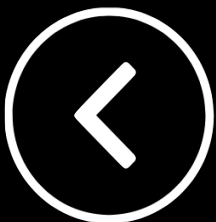

G-CEM Veneer




Introducing G-CEM Veneer: predictable aesthetics, high mechanical properties and unique flow.* G-CEM Veneer has all the ingredients to deliver a stunning result for a confident smile. The secret is in its optimally balanced formula with Full-coverage Silane Coating (FSC) technology. This innovative filler treatment allows a high filler rate to be achieved with homogeneous repartition of fillers in the matrix. The result: unique handling and excellent physical properties.

Clinical Indications


Adhesive luting of ceramic and composite veneers, inlays and onlays with a thickness (< 2.0mm) and translucency that enables the complete light-curing of the cement



View product brochure

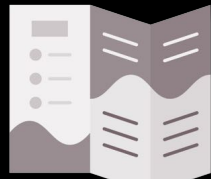

View Luting Guide brochure

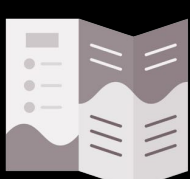

View Adhesive Resin Cement brochure




View instructions for use




View Technique Guide

Features & Benefits of G-CEM Veneer

- Stays in place when untouched and it is fluid while put under pressure
- Covers the entire surface uniformly without voids, while the film thickness is kept to a minimum
- Not sticky; it is easy to apply with high precision and it allows easy excess removal
- Unique consistency allows for an optimally balanced flow, without the need for preheating
- Well suited for the cementation of aesthetic restorations and offers high colour stability
- G-CEM Try-in Pastes, for trial-fit of your indirect restorations
- The shades of G-CEM Try-in Paste perfectly match the four shades of G- CEM Veneer and G-Cem LinkForce
- Use of G-CEM Try-in pastes will aid in the verification of the selected shade for the desired aesthetic result
- Water-based for easy removal after application

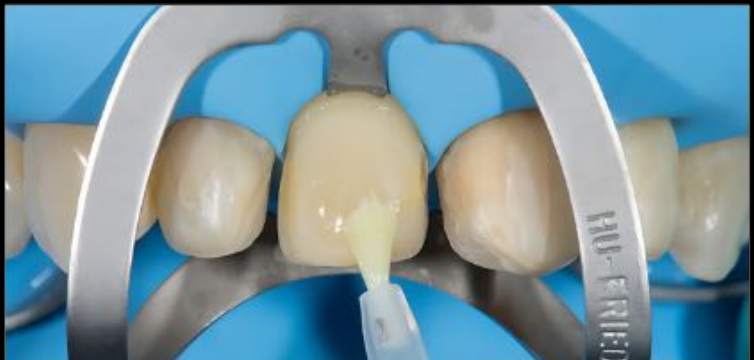

Dr. J Tapia Guadix,

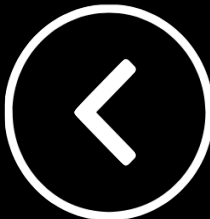

GC Fuji PLUS CAPSULE



Fuji PLUS has been designed to offer you simple and error-free cementation procedures and a strong and durable bond you can rely on. The outstanding performance and the availability in various shades make GC Fuji PLUS a multi- purpose cement.

Fuji PLUS is a resin modified glass ionomer luting cement that has consistently delivered stress-free cementation for routine crown and bridge restorations. It provides the highest level of confidence through exceptional long-term clinical results. Fuji PLUS has an extended working time for luting long-span bridge work, combination work and luting of several restorations in one step.


Clinical Indications - Cementation of

All types of metal, porcelain fused to metal, and resin crowns, inlays, onlays and bridges

Porcelain inlays

All types of core material including prefabricated pins, posts and cast posts

Cementable reinforced all ceramic crowns (e.g. Procera)

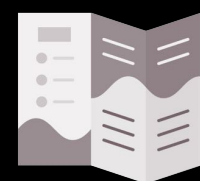

[View product brochure](#)



[View instructions for use](#)



[View capsule activation guide](#)



[View Luting Guide brochure](#)

Features & Benefits

Easy mixing and handling like a conventional cement

Similar mechanical properties to resin cements

Intrinsic adhesion to tooth structure

Minimal film thickness

No post-operative sensitivity

Radiopaque

In combination with the GC Fuji PLUS Conditioner

Optimal marginal seal

Continued high fluoride release

Step-by-step technique

1. Clean the prepared teeth with pumice and water. As an option for increased retention - apply Fuji PLUS conditioner for 20 seconds. Rinse and gently dry.




2. Mix GC Fuji PLUS. Coat the internal surface of the restoration with 1mm of cement and seat within 30 seconds.



3. Remove excess cement at the first formation of gel stage (approx. 1 minute after seating).



4. Maintain isolation until the final set is verified (approx. 4 minutes).

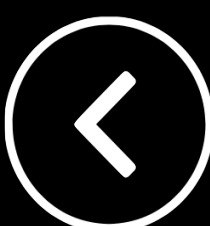

GC CAD/CAM & TEMPORISATION

Click a product image to view product information



GC Initial LiSi Block



Cerasmart 270



OPTIGLAZE™ color



Tempsmart



Cavition

ALWAYS READ THE LABEL AND FOLLOW THE DIRECTIONS FOR USE



GC Initial[®] LiSi Block



Lithium Disilicate CAD/CAM Block for Chairside Solutions

A fully crystallised lithium disilicate block that delivers optimal physical and aesthetic properties without firing.

GC Initial LiSi Block is a fully crystallised lithium disilicate block that delivers optimal physical and aesthetic properties without firing.

This unique block features GC's proprietary High Density Micronisation (HDM) technology for CAD/CAM dentistry to deliver high wear resistance, smooth margins and aesthetic final results. This makes it an ideal time saving solution for single visit chairside treatments.

The unique ultra-fine crystal makes it easy to grind and can be quickly milled in its fully crystallised stage. Thanks to GC Initial LiSi Block's exceptional properties, high gloss and natural opalescence can be obtained in few minutes by polishing only and restoration is then ready for luting.

Optimised acid and wear resistance help to preserve the aesthetics of your restorations over time. GC Initial LiSi block is available in four aesthetic shades and two translucencies. When more characterisation is required, the GC Initial[®] ceramic line comprises fully synergistic aesthetic solutions.

Why choose GC Initial[®] LiSi Block?

Reduced chair time

Just mill, polish and place!

Initial LiSi Block offers reduced process time

No need to fire because it is fully crystallised - saving 40%+ of chairside time

Durable aesthetic & smooth margins

Optimised acid and wear resistance | helps preserve the aesthetics of your restorations over time

Excellent edge stability for smooth margins

Natural opalescence

Available in both high and low translucencies


GC Initial LiSi Block offers a natural opalescence in any light

Choose your preferred finishing procedure:

Superior gloss can be obtained in a few minutes by polishing only

For sophisticated aesthetic cases, remarkable results can be achieved with

GC Initial Lustre Pastes NF and Initial Spectrum Stains

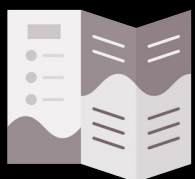

View product brochure



View instructions for use



View Milling Parameters




CEREC SPEEDFIRE Guide




View clinical paper



Scientific References

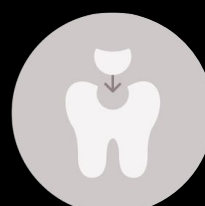

Watch LiSi Polishing Video

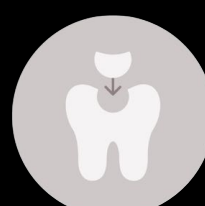

Watch Dr Javier Tapia Gaudix Game Changer Presentation



Watch One & Done Webinar



CERASMART 270

CERASMART 270 is an aesthetic hybrid block made to provide fast and long lasting clinical success.

Clinical Features & Benefits

Strength and flexibility blend together

Full Coverage Silane Coating (FSC) nano-filler technology

Highest flexural strength in its category

High flexibility to buffer masticatory pressure

Precision and durability unite

Fast and precise milling, sharp marginal adaptation

Long-lasting polish and less wear to the opposing dentition

High radiopacity for an easy follow-up

Aesthetics and CAD/CAM connect

Balanced fluorescence and opalescence

Characterise with OPTIGLAZE (Pre-treat, apply, light cure)

Easy to obtain a high gloss surface

CERASMART™ 270 for CEREC - refill of 5 blocks in sizes 12, 14 & 14L

HT shades (A1, A2, A3, A3.5, B1)

LT shades (A1, A2, A3, A3.5, B1)

Bleach shade (BL)




Clinical Indications

Quadrant restoration - Partial and full-coverage crowns



Dr. Dimitrov

Composite resin replacement - Inlay and Onlays

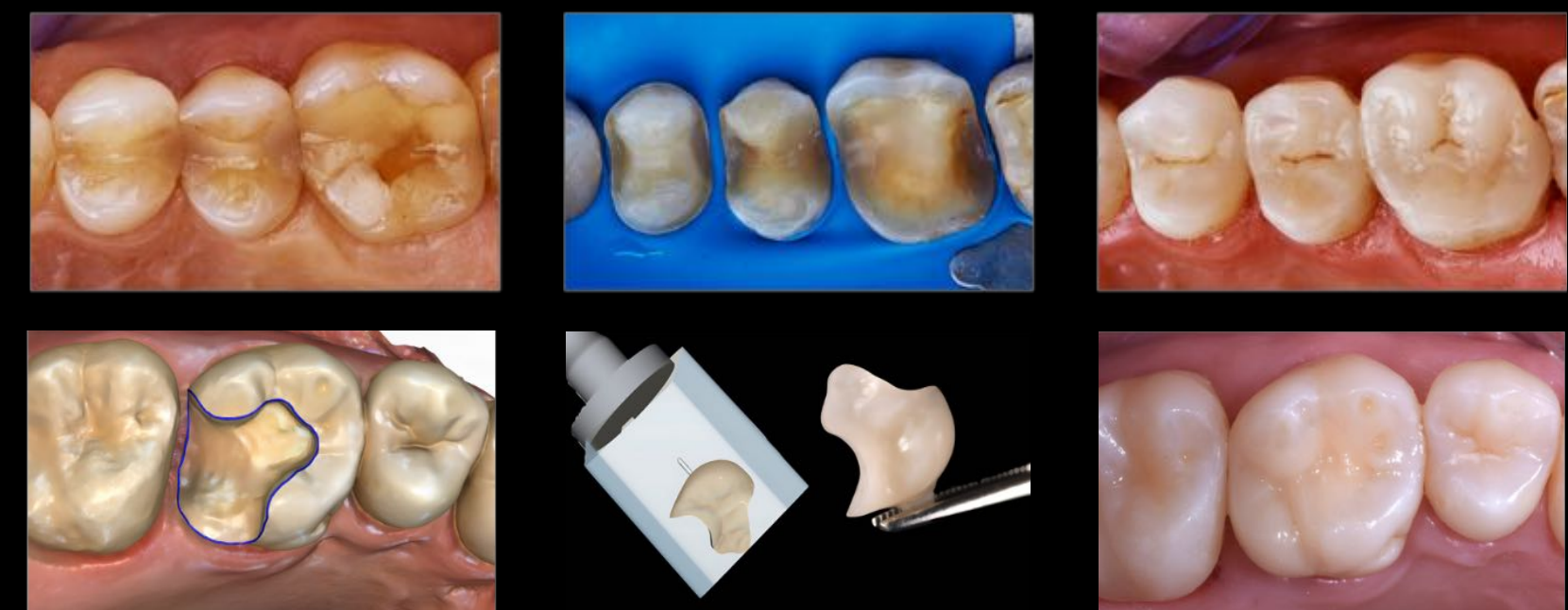


Dr Witowski

Worn dentition - Palatal Veneers



Dr Javier Tapia

GC recommends the use of **G-Multi PRIMER** for the pre-treatment of CERASMART. CERASMART blocks are CEREC® compatible and are listed in the latest CEREC® software update™

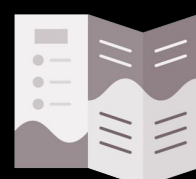

[View product brochure](#)



[View instructions for use](#)



OPTIGLAZE™ color

Are you looking for a simple solution to add characterisation to your everyday aesthetic creations?


Then OPTIGLAZE colour will impress you because it adds colour, surface gloss and high wear resistance to all your restorations. Ready to use, easy to handle and cutting on the polishing stage this product will save you valuable time.

You can choose from a wide variety of colours, both for internal and external characterisation, leading to perfect aesthetics.

The renowned nano-filler technology gives your restoration a high wear resistance, for long-lasting gloss with a high discolouration resistance.

Clinical Features & Benefits

- A wide variety of colours - adds natural gloss and brilliance
- Thin 25-50 um coating layer - long-lasting gloss
- High discolouration resistance & high wear resistance
- Excellent bonding durability

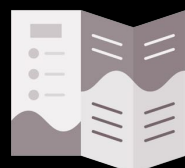

[View product brochure](#)



[View instructions for use](#)



Kit contains

- 15 x 2.6ml bottles Shades: A Plus, B Plus, C Plus, White, Ivory White, Yellow, Orange, Pink
- Orange, Pink, Red Brown, Olive, Lavender, Grey, Blue, Red
- 2 x 5ml bottles - Clear, Clear HV
- Disposable dispensing dish
- Flat brush & Round brush




Milled CERASMART™



Milled CERASMART™ with
OPTIGLAZE™ Color

Imagos: Dr Anthony Mak (Australia), Brad Grobler (Oral Dynamics NZ)

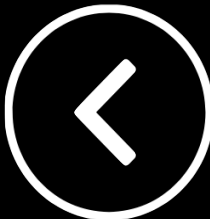

TEMPSMART®



Clinical Features & Benefits

High fracture toughness and flexural strength for clinical longevity.

Non-sticky - easy to handle

Low volumetric shrinkage - just 0.86% ensuring accurate fit

Requires minimal polishing with high gloss retention

Dual-cured for rapid setting with a three-unit bridge completed in just 3 minutes

Repairable with TEMPSMART® or a high strength flowable composite

Option to glazed with EQUIA® Forte Coat or G-COAT PLUS™

Option to characterised with OPTIGLAZE

48ml cartridge delivers approx. 64 temporary restorations

Available in A1, A2, A3, A3.5, B1 and BW shades

Clinical Indications

Single-unit crowns


Bridges, inlays & onlays

Veneers

Implants and long term temporisation

TEMPSMART® is an innovative, dual-cured, bis-acrylic composite temporary crown and bridge material using micro-filled resin (MFR) and nano-filler technologies.


It maintains a smooth surface after polishing because of this revolutionary combination. The high density polymer network for TEMPSMART® makes it a strong material ideal for any sized provisional. This product is ideal for temporary crowns, bridges, inlays, onlays and veneers, as well as long-term temporary restorations.

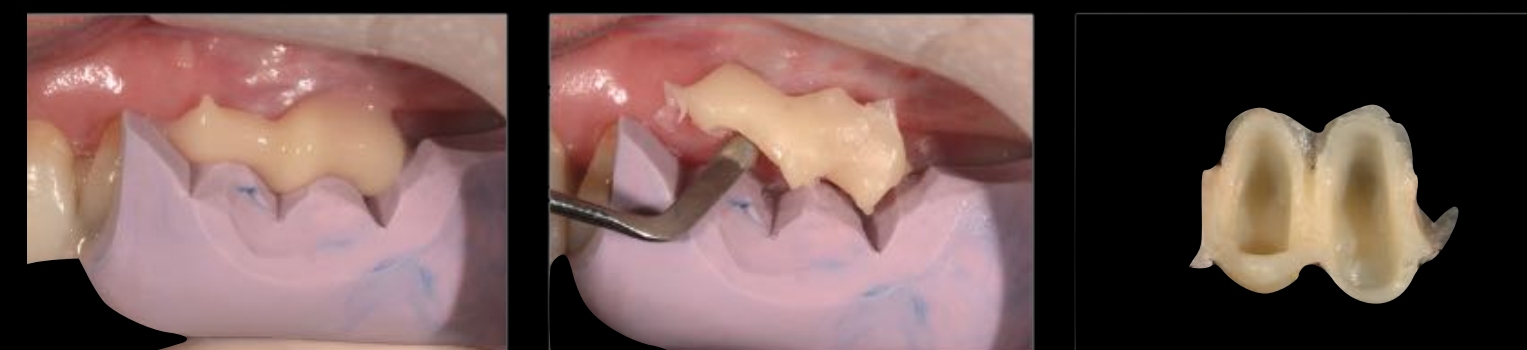

TEMPSMART® temporary restoration

Dr Anthony Mak


Dr Chee Chang

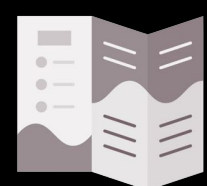

[View product brochure](#)



[View instructions for use](#)

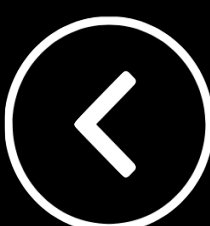

CAVITON

Hydraulic temporary restorative

Clinical Features & Benefits

Setting initiated by water and saliva

Excellent adaptability

Easy to remove


Clinical Indications

Temporary filling material for standard cavities

After endodontic treatments



[View instructions for use](#)

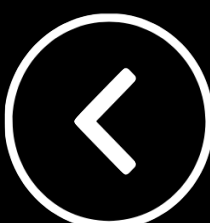

GC IMPRESSION & RELINE

Click a product image to view product information



COE Spacer Trays



GC AROMA FINE PLUS



GC EXACLEAR



EXA ADVANCED



EXAFLEX



GC Reline II Soft / Extra Soft



Soft-Liner



Exabite III



Fit Checker
Advanced White

ALWAYS READ THE LABEL AND FOLLOW THE DIRECTIONS FOR USE




COE Spacer Tray

Features & Benefits


- Disposable, strong, rigid plastic, excellent stability
- Handle will remain intact when removing the impression
- Unique, zigzag spacer, Anatomically designed
- Avoids pressure on lingual frenum and soft tissue
- Provides increased retention and maximum compression
- Perforate style/solid style
- Lessens the risk of cross contamination; no clean up or sterilisation
- Minimises contact of cusps to tray,
- Ensuring sufficient bulk of impression around the teeth

The right impression material and a good impression-taking technique go a long way to ensuring the best results; and the tray is the foundation of every accurate impression. As one of the world's leading manufacturers of impression trays and materials, GC offers a wide range of stock trays, as well as accessories and materials to make custom trays.

Full Arch Solid

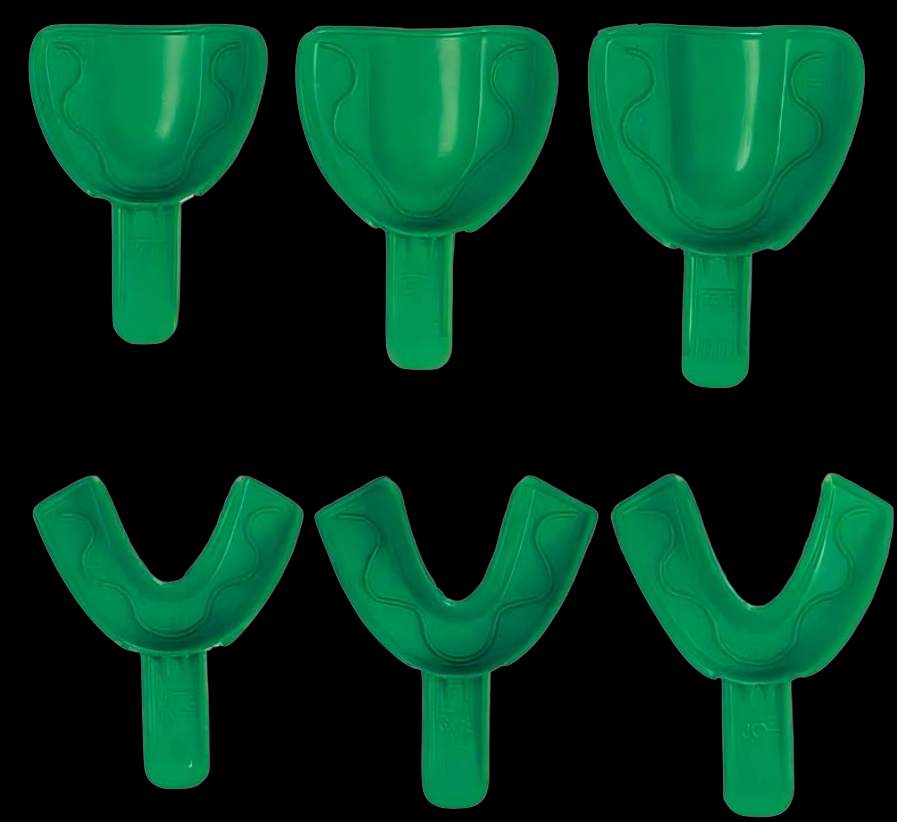

Full Arch Perforated

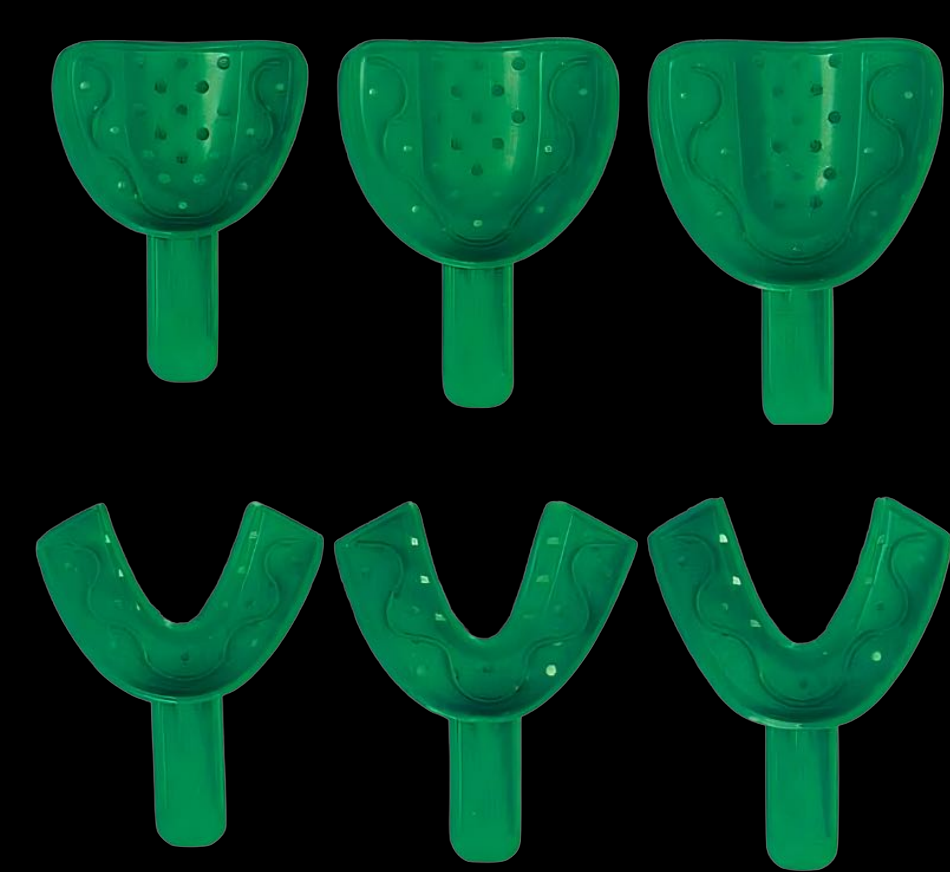


Quadrant



[View product brochure](#)



[View instructions for use](#)

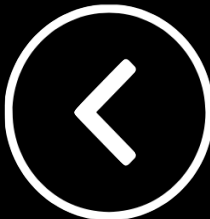

AROMA FINE PLUS

Based on GC's extensive knowledge of alginate materials, AROMA FINE PLUS has been designed to offer users final impressions of the highest quality combined with an easier and faster mixing.

Advanced formulation

Easy and rapid powder incorporation into water within just 3 seconds

Dust free formulation allowing fast and vigorous mixing for a smooth mixture without bubbles

Available in 2 setting types: Fast set (green 90 sec) and Normal set (pink 120 sec)

Ideal thixotropic property



Does not leak down the patient's throat during impression taking

Under pressure the alginate flows into the narrowest spaces for highest detail reproduction

Clinic and Laboratory optimised

Compatible with a wide range of Type IV stones for a better detail reproduction

Smooth model surface

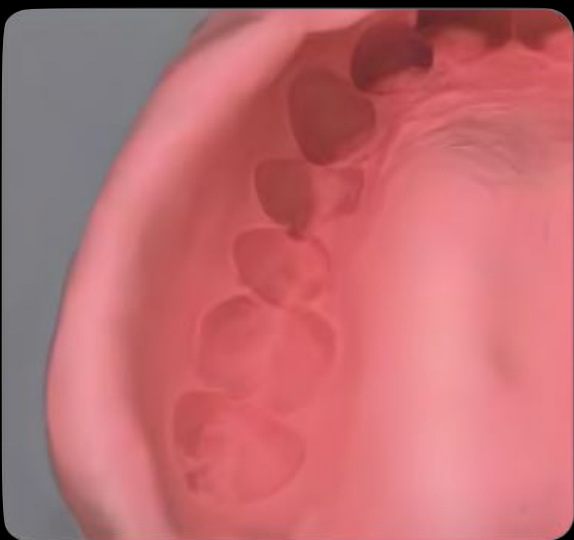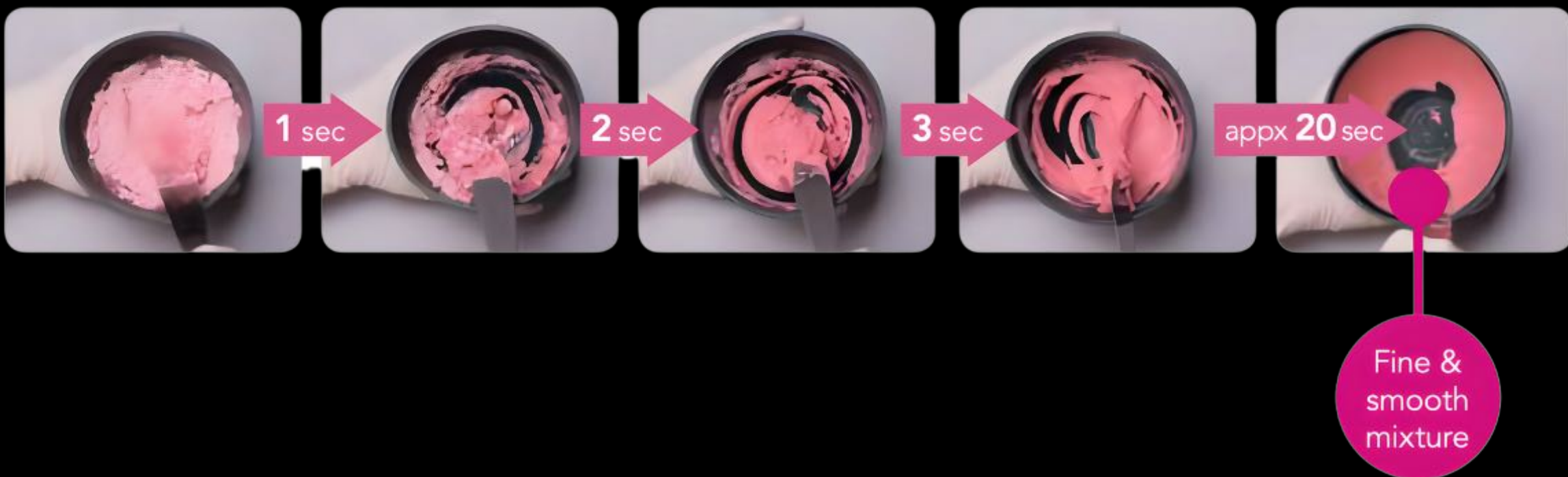


Easy-mixing of AROMA FINE PLUS results in a creamy, bubble-free consistency that captures details accurately.




AROMA FINE PLUS provides a balance of good drip resistance and fluidity to effectively provide coverage enabling detailed impression taking with minimal discomfort to patient.



AROMA FINE PLUS has good compatibility with hard stones producing smooth and precise detail reproduction on the gypsum model.

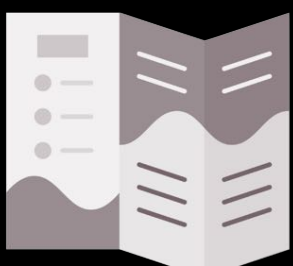

[View product brochure & Instructions for Use](#)

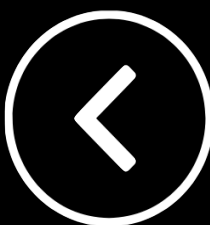

EXACLEAR




Clinical Indications

Injection Moulding Technique (IMT)
Clear stent fabrication

Features & Benefits

- High transparency - crystal clear for absolute control
- Very effective light-curing through the silicone
- Optimal consistency - balanced handling and flow under pressure
- Intra-oral setting time - 2 mins
- Extra-oral setting time- 7 minutes

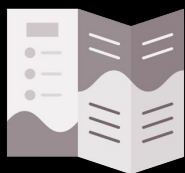

View product brochure



View instructions for use



View Injection
Moulding Guide



View Anterior Injection
Moulding Video




View Posterior Injection
Moulding Video

EXACLEAR™ is an innovative, clear silicone material that fulfils the needs of demanding aesthetic cases and facilitates consistent results. It's absolute transparency offers total visual control making it ideal for creating outstanding aesthetic results in complex clinical situations.

EXACLEAR™ offers the perfect balance of properties for ease of handling, staying in place when extruded but flowing nicely under pressure, resulting in perfect adaptation and reproduction of details. EXACLEAR™ is ideally suited for both chair side and laboratory use, the injection moulding technique and clear stent fabrication for creation of provisional and permanent restorations.



Dr. Anthony Tay

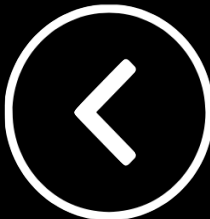

EXA Advanced™

Hydrophilic VPS Impression Material


Create accurate reproductions of the smallest details with EXA Advanced Hydrophilic Vinyl Polysiloxane Impression Material.

EXA Advanced is a significantly advanced addition reaction (platinum catalysed) VPS impression material.

It produces crisp, extraordinarily smooth and highly precise impressions for models, crowns, bridges and dentures. Available in Normal and Fast setting options.

Clinical Features & Benefits

- Rigidity and elasticity in balance for optimal results
- Resistant to tearing for an accurate reproduction of the finest preparation margins
- Hydrophilic to capture all details required for an accurate reproduction
- Elasticity, together with a pleasant smell and taste for the best possible comfort for patients and predictability for the clinician
- Recovers from the deformation that takes place upon removal from the patient's mouth
- Maintains its accuracy during storage and transportation, due to excellent dimensional stability

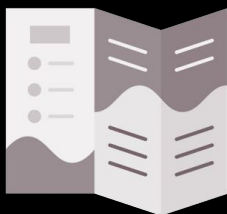

[View product brochure](#)




[View instructions for use](#)



Exa Advanced Normal & Fast Set	Working time at 23°C	Intraoral setting time
EXA Advanced Normal set (Injection and Regular)	2'00"	4'00"
EXA Advanced Fast set (Injection and Regular)	1'00"	1'30"

Double mix	● Heavy body	● Regular or ● Injection
Double arch – three function	● Heavy body	● Regular or ● Injection
Putty-wash Also 2-step putty wash	● Putty	● Regular ● Injection
Simultaneous Also 1-step putty wash or dimple technique	● Putty	● Regular or ● Injection

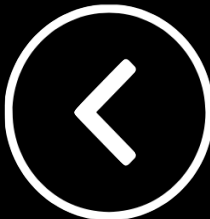

EXAFLEX®

Hydrophilic VPS Impression Material


Thixotropic impression material for all popular impression techniques.
The adapted flow, extended working time and the snap set make this an extremely easy to use A-silicone material. The reaction silicone also has a unique combination of innovative, clinical advantages for impression taking.




Clinical Features & Benefits

- Easy to mix, Optimised thixotropic properties
- Extended working time & Sharp setting curve
- Hydrophilic properties
- Extremely high precision
- Ideal affinity to modelling stone materials
- High recovery from deformation
- High dimensional stability
- High elasticity
- Very tear resistant
- Can be disinfected



2 Step Putty / Wash Technique

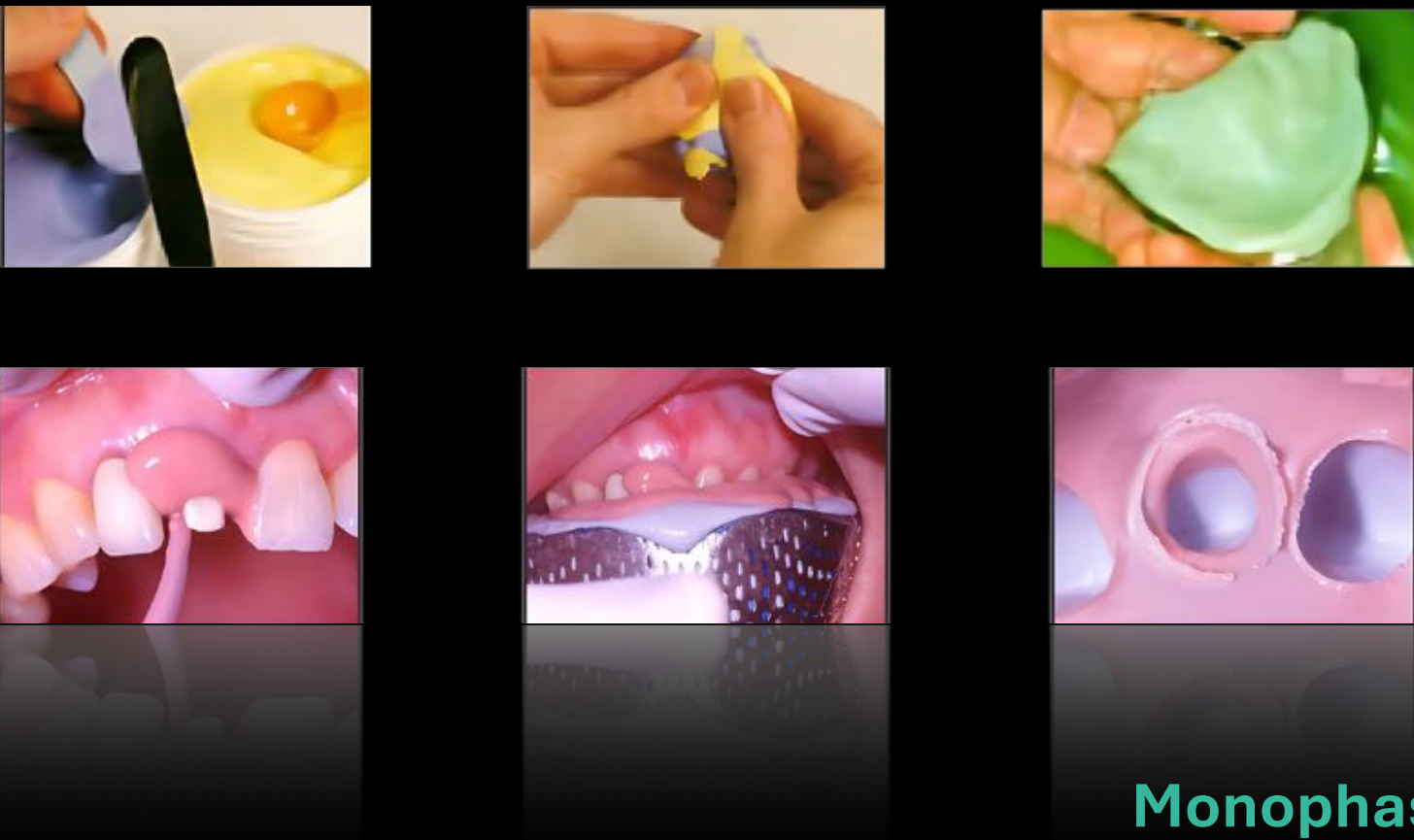

Double Impression Technique

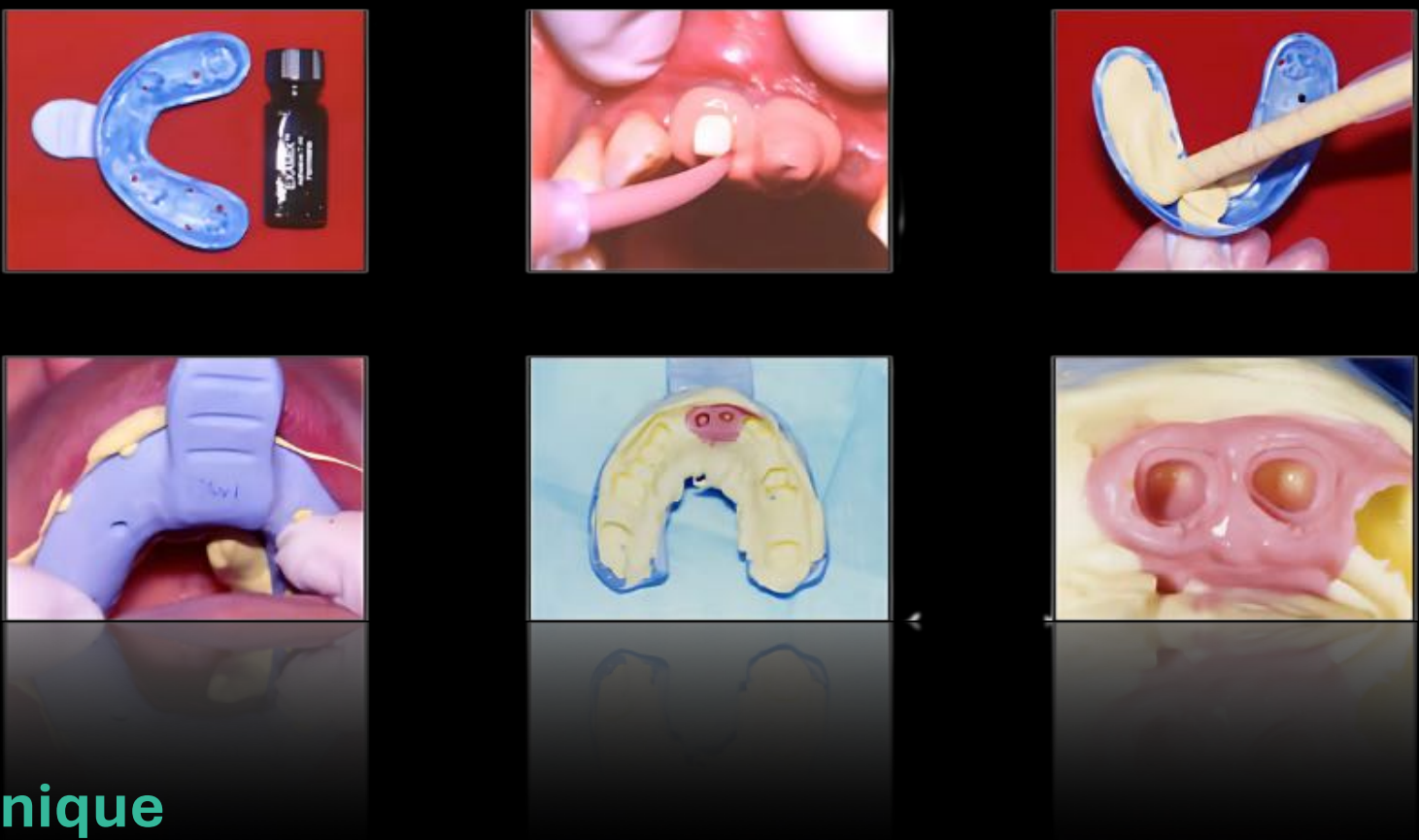

Monophase Technique

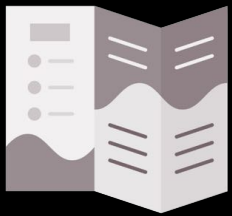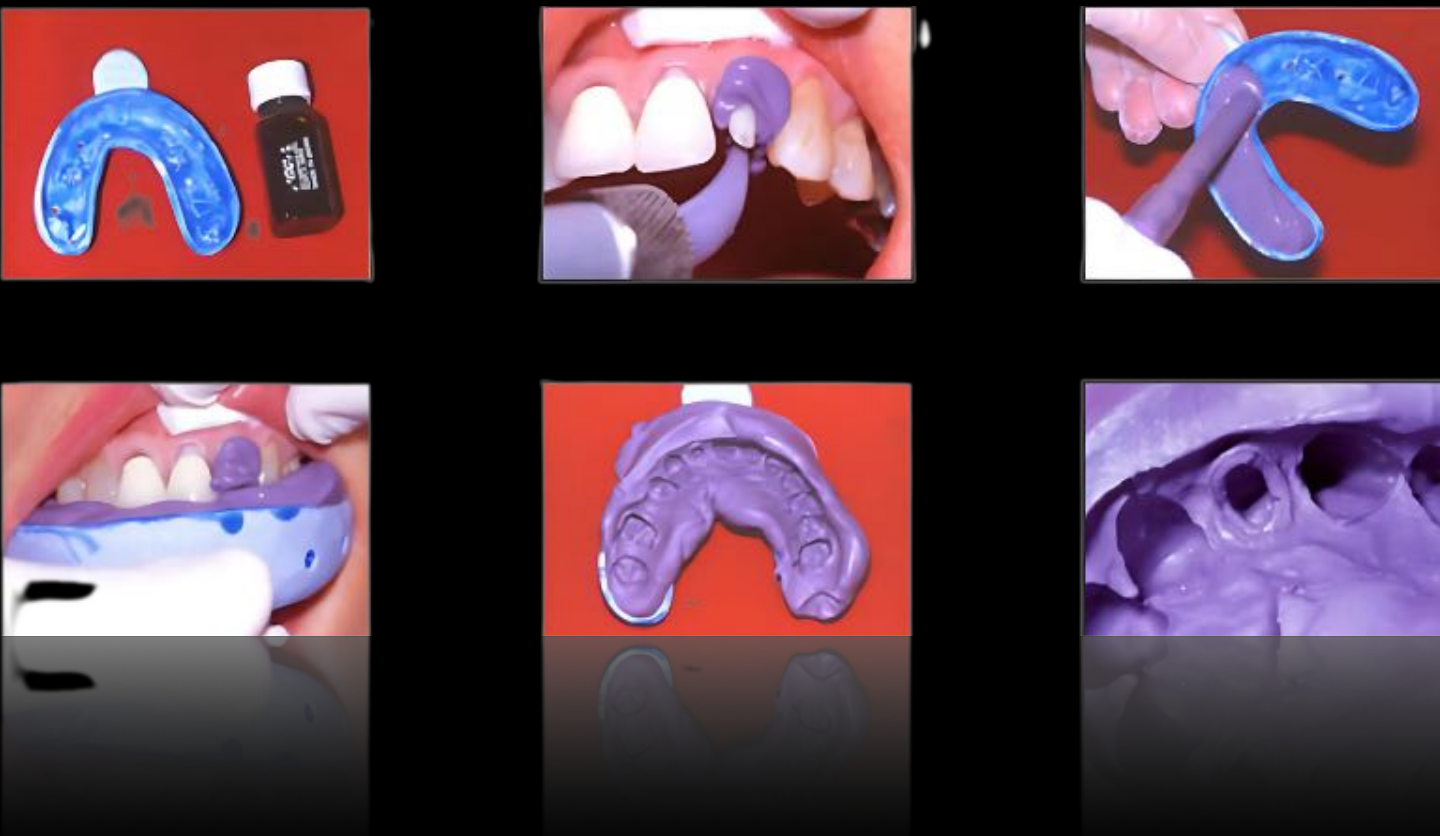


[View product brochure](#)



[View EXAFLEX
Instructions for use](#)



[View EXAFLEX PUTTY
Instructions for use](#)

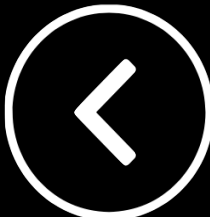

GC RELINE II Soft & Extra Soft

A-silicone reline material providing comfort and long-lasting elasticity.

GC Reline II Soft is indicated for sharp or bony alveolar ridges and undercuts. It is used to adapt dentures after extractions or implant surgeries.

Clinical Features & Benefits

Resilient: acts as a cushion between the denture and the oral tissue to buffer masticatory pressure and

solve anatomical problems

Can stay in place for weeks or months

Convenient cartridge delivery

Ideal viscosity, flowing under pressure to adapt to anatomical details

Long-lasting comfort & gentle relief: the material's surface remains elastic, smooth and clean over time

RELINE II Primer for Resin provides a long-lasting adhesion to the underlying denture

RELINE II Remover for Resin makes it possible to easily remove & replace the relining layer

RELINE II Finishing Material can be used to further seal the margins

Indications

Chairside and lab relining procedures

Relining of dentures - thin sensitive mucosa and resorbed ridges

Relining of dentures - sharp and bony ridges or bony undercuts

Relining of implant-supported dentures



Initial denture, to be relined following the placement of 4 implants.



Relieve the area to be relined and roughen the surface.



Apply GC RELINE II Primer on the area to be relined. Gently dry with air.



Easy cartridge delivery.



Apply GC RELINE II (Extra) Soft. Working time is 2'00".



Seat the denture in the mouth and muscle trim. Retain in situ for 5'00".



Remove the excesses using a scalpel or scissors. Trim with GC RELINE II POINT FOR TRIMMING.



Finish with GC RELINE II WHEEL FOR FINISHING.



Relined denture.

• Dr E. D'Incau


	GC RELINE II Extra Soft	GC RELINE II Soft
Material	A-Silicone	A-Silicone
Indications	• Tissue toning material • Sensitive mucosa • Pressure points • Normalising irritated tissues	• Implant cases • Pressure points • Bony undercuts • Pain reduction & increase of comfort • Long-lasting relining in case of missing supportive tissue
Time in mouth	Days/weeks	Weeks/months



[View product brochure](#)



[View instructions for use](#)

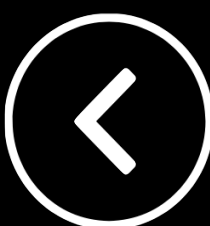

EXABITE III

A-silicone material with properties especially adjusted to the requirements of bite registration.



Clinical Features & Benefits

- EXABITE III is fast - it sets in just 30 seconds, with 15 seconds working time in the mouth
- Suitable for all bite registration techniques
- Open bite registration for standard crown and bridge indexes
- Closed bites for edentulous or partial patients
- Adapts well over wax bite blocks when taking denture registrations
- Suitable for CAD/CAM fabrications, as it is easy to powder and obtain a flawless optical image
- Creamy consistency, no resistance to biting
- No taste, no heat, no smell or odour

Easy to read & trim

EXABITE III makes it much easier to read the occlusal detail with its contrasting green colour. Light reflecting off the surface results in greater detail being seen.

Shore A hardness of 94 guarantees there will be no distortion or movement. This enables fast and easy trimming with no displacement.




Ideal consistency for closed bite technique



30 seconds setting time in the mouth

Dr Jim Ironside


Physical Properties	EXABITE III
Working time 23°C (sec)	35
Working time 37°C (sec)	15
Minimum time in mouth (sec)	30
Strain in Compression (%)	0.7
Recovery from deformation (%)	99.2
Linear dimensional change (%)	-0.3
Shore A Hardness	94

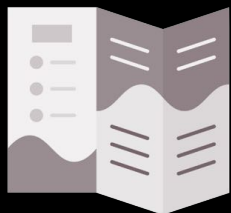

[View product brochure](#)



[View instructions for use](#)

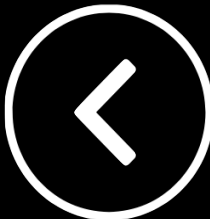

FIT CHECKER™ ADVANCED WHITE

Vinyl Polyether materials for checking pressure spots and the fit accuracy of prosthetic works.

Thanks to their very thin consistency, sharp setting time and hydrophilicity, both materials are very convenient to use and offer an accurate identification of misfits of prosthetic works of all kinds.

The ideal colour for acrylic and metal FIT CHECKER™ ADVANCED WHITE offers ideal consistency and high wettability for easy flow onto dentures and saliva-coated oral tissues.

Perfect for all kinds of prosthetic work, such as dentures and metal-based crowns and bridges.

Clinical Features & Benefits

Ideal consistency for easy flow onto dentures and oral tissues

Excellent wettability against saliva coated surfaces

Sufficient working time - up to 1min

Rapid setting time - 1 min in the mouth

A perfect colour gradation for clear and precise visualisation of poorly fitting areas

Easy adjustments - the set material does not get caught by the bur while relieving pressure spots.

Easy to peel off ensuring a fast, clean procedure.



Clinical Indications

Checking fit for indirect restorations - ceramics, metals, zirconia, acrylic

Evaluation of occlusal contacts

Fit of occlusal splints

Can be used for bite registration



Apply & Seat



Check & Adjust



Ideal for metal restorations




Seat



Check & Adjust




Bite Registration

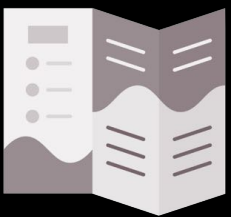

[View product brochure](#)



[View instructions for use](#)

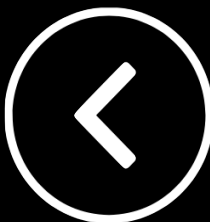

SOFT-LINER - TISSUE CONDITIONER

Soft relining material for short-term use, including functional impression taking using existing removable dentures.

- Tissue conditioning
- Temporary relining
- Functional impression taking

Contradictions


Avoid use of this product in patients with known allergies to methacrylate polymer, ethanol or phthalic ester-based plasticiser.


Clinical Features & Benefits

- Fast application
- No heat irritation
- Easy to trim
- Good retention
- Smooth surface minimises unwanted food debris retention
- Simple to remove from the denture at the end of the treatment



[View instructions for use](#)

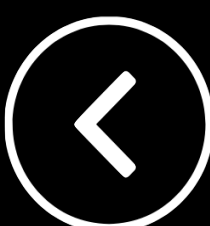

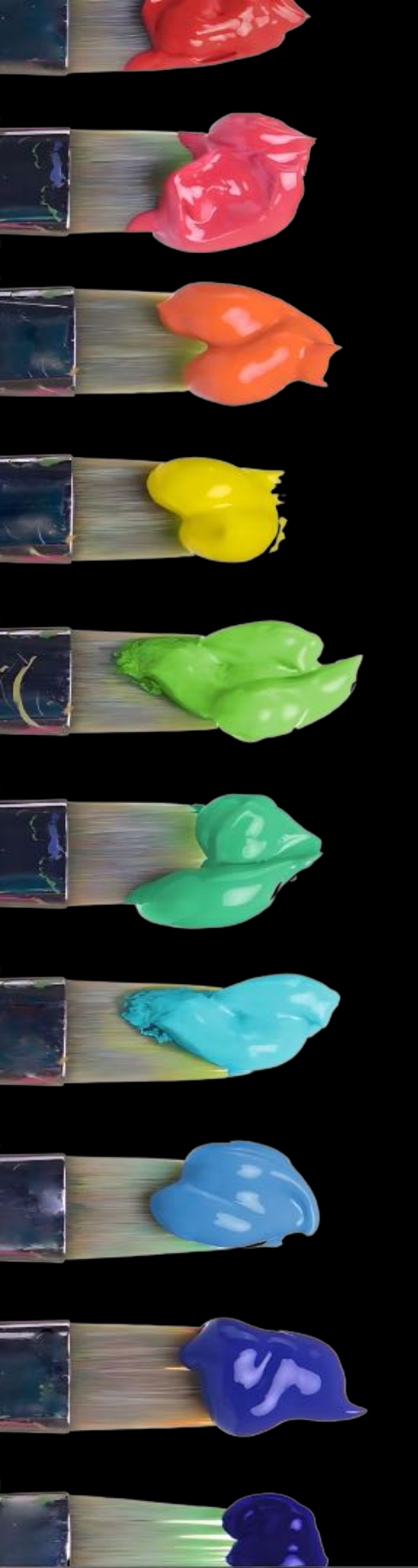


Subscribe to the GC CPD Monthly Newsletter

Webinars | Courses | Events



GC SMART GUIDE INSTRUCTIONS

1. Click to select a master product category

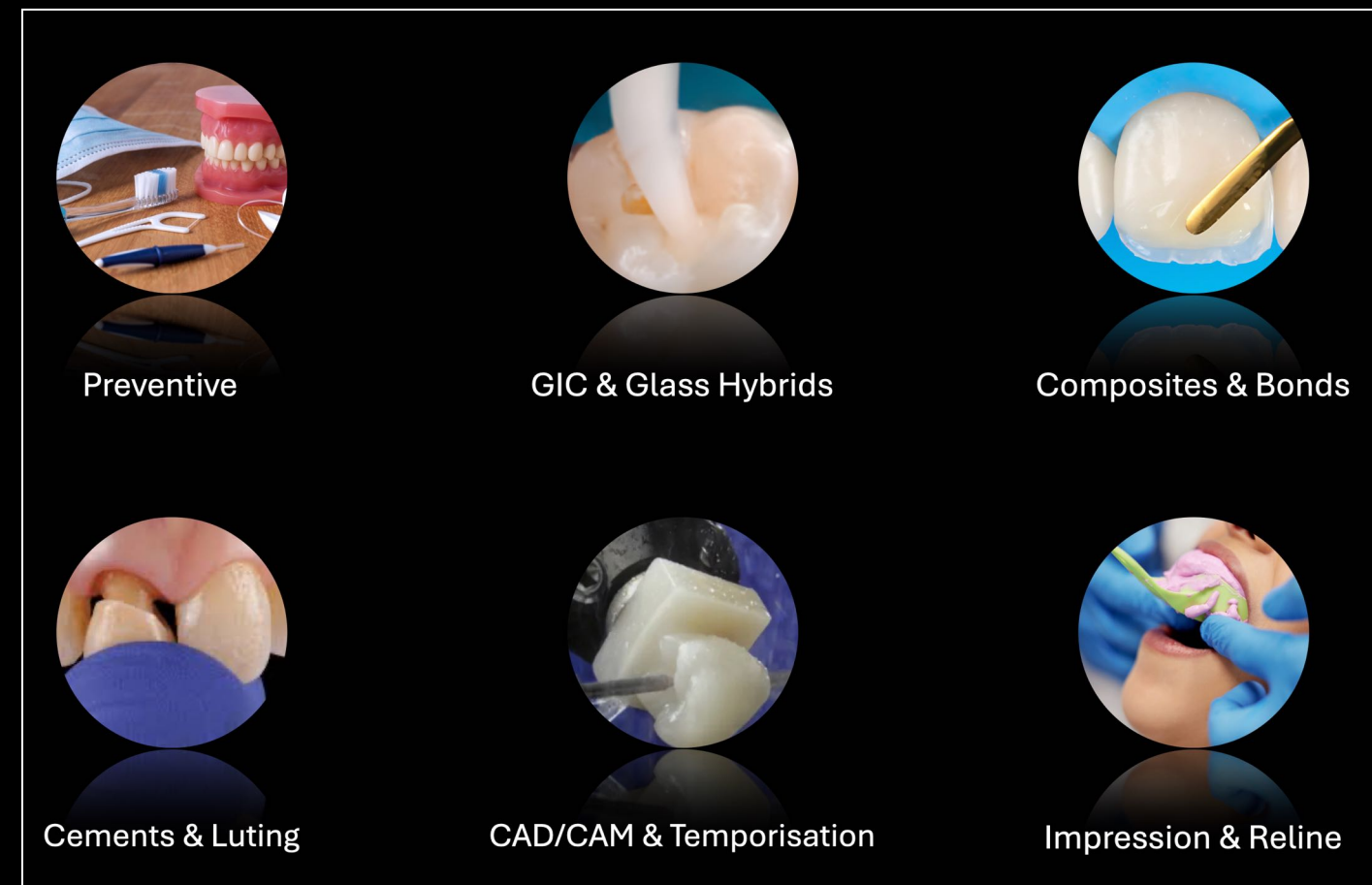

2. Click a product image to view product information



7. Click the product name to open the Henry Schein Online product order page



3. The selected product information will be displayed



4. Click the back button on any product page to go back to the product selection page for that category



5. Click the home button on any page to go back to the master product category selection page



6. Click the icons on the product information page to view further information, brochures and instructions for use



GC NZ SALES PRODUCT CONTACT

Click an email address to automatically open your e-mail client and send an email or click a phone number to call a product specialist

**Northland | North Shore | Wellington |
Christchurch | Hawkes Bay | Nelson
Marlborough**

Tamera Greenwood

tamera.Greenwood@gc.dental

021 689 016

Auckland | Waikato | BOP

Karyn Bedford

karyn.bedford@gc.dental

021 999 938

**Auckland | Christchurch | Otago | Central Otago
Southland | South Canterbury**

Colleen Coulter

colleen.Coulter@gc.dental

021 433200

GC Australasia Dental Pty Ltd (Head Office)

1753 Botany Road Banksmeadow, NSW 2019 Australia, info.australasia@gc.dental, www.gc.dental/australasia, +61 2 9301 8200

