KOOLINER™

Prior to use, carefully read the instructions for use. 

HARD CHAIRSIDE DENTURE RELINE

For use only by a dental professional in the recommended indications.

RECOMMENDED INDICATIONS


A temporary lining for acrylic dentures. For use in chairside procedures.

CONTRAINDICATIONS

Patients who have shown sensitivity to methacrylates. In case of allergy refer to a physician. Not intended for permanent lining.

DIRECTIONS FOR USE

- Preparation of the denture:**
Relieve and roughen the area of the denture to be relined. Clean and dry the denture thoroughly. Coat labial and buccal surfaces of the denture with COE LUBRICANT. Do not apply coating within 3mm (1/8 inch) of the peripheral border. If the denture has plastic teeth also protect them with COE LUBRICANT.
- Preparation of KOOLINER:**
Recommended powder / liquid ratio is 15ml powder to 6ml liquid. Pour liquid into the mixing cup and then add the powder slowly. Stir thoroughly for no more than 30 seconds and avoid the introduction of air bubbles.

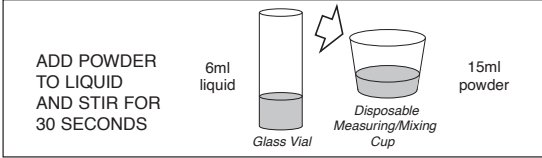- Application:**
After approximately 1-2 minutes, spread the mixture of KOOLINER over the area to be relined. Seat the denture in the manner of taking an impression and instruct the patient to close lightly into occlusion. After 3 minutes, instruct the patient to move lips and cheeks so that a muscle trimmed periphery is obtained. Remove the denture and rinse under cold water. Trim away excess material. Re-seat the denture and instruct the patient to close FIRMLY into occlusion, and to hold this position for 5 minutes. Remove the denture and rinse again in cold water.
- Finishing:**
Peak curing temperature: Approximately 40°C/104°F at 6-15 minutes when tested according to ADA/ANSI specification number 17. In thicker applications, peak temperature may exceed that stated above, possibly producing a hazardous condition in the mouth during curing. When curing is complete, trim away excess. For smoothing the edges, use a hot spatula.

STORAGE

Recommended for optimal performance; store at temperature of 4-25°C (39-77°F).

PACKAGES

345001 KOOLINER Professional Package

345002 KOOLINER Powder, 80 g

345091 KOOLINER Liquid, 55 mL

CAUTION

1. Patient should clean daily to remove food deposits and plaque. Recommend a commercially available denture cleaner and brush. Do not recommend toothpastes or hard bristle brushes as they may damage the denture liner or denture.
2. If patient notices damage to denture (e.g., chipping, delamination, distortion, etc.) or changes in tissue condition (e.g., inflammation, discomfort, allergic reaction, etc.) have patient discontinue use and return for evaluation and consultation.
3. Personal Protective Equipment (PPE) such as gloves, face masks and safety eyewear should always be worn.
4. Ensure good ventilation/exhaustion at the workplace. Keep ignition sources away.

CLEANING AND DISINFECTING RECOMMENDATION

MULTI-USE DELIVERY SYSTEMS: To avoid cross-contamination between patients the bottles and measuring devices require mid-level disinfection. Immediately after use inspect the bottles, measuring devices and label for deterioration. Discard if damaged. DO NOT IMMERSE. Thoroughly clean bottles and measuring devices to prevent drying and accumulation of contaminants. Disinfect with a mid-level registered healthcare-grade infection control product according to regional/national guidelines.

Some products referenced in the present IFU may be classified as hazardous according to GHS. Always familiarize yourself with the Safety Data Sheets available at: <https://www.gc.dental/europe> or for the Americas: <https://www.gc.dental/america>

They can also be obtained from your supplier.

Undesired effects- Reporting:
If you become aware of any kind of undesired effect, reaction or similar events experienced by use of this product, including those not listed in this instruction for use, please report them directly through the relevant vigilance system, by selecting the proper authority of your country accessible through the following link:
https://ec.europa.eu/growth/sectors/medical-devices/contacts_en as well as to our internal vigilance system: vigilance@gc.dental
In this way you will contribute to improve the safety of this product.

UK Responsible Person
GC UNITED KINGDOM Ltd.
Coopers Court Newport Pagnell, Buckinghamshire, MK16 8JS
United Kingdom

Last Revised: 02/2024

MANUFACTURED BY
GC AMERICA INC.
3737 West 127th Street, Alsip, IL 60803 U.S.A.
TEL: +1-708-591-0900
www.gc.dental/america

DISTRIBUTED BY
GC AUSTRALASIA DENTAL PTY. LTD.
1751 Botany Rise, Bardon, Queensland, NSW 2107, Australia
TEL: +61-2-9301-8200

DISTRIBUTED BY
GC AMERICA INC.
3737 West 127th Street, Alsip, IL 60803 U.S.A.
TEL: +1-708-591-0900
www.gc.dental/america

GC EUROPE N.V.
Flevostraat 399, Scharnestein, 3845 PA, The Netherlands
B-3001 Leuven, Belgium TEL: +32-16-74-10-00

GC SOUTH AMERICA
Rua Heliodora, 399, Sarambaia - São Paulo, SP, BRASIL
CEP: 02022-051 - TEL: +55-11-2925-0965
CNPJ: 08.279.968/0001-61
REGIST. TEC. Enca de Lina - CHOSP 100.866

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