Prior to use, carefully read the instructions for use.

EN

G-ænial™ Universal Injectable

UNIVERSAL LIGHT-CURED RADIOPAQUE INJECTABLE COMPOSITE

For use only by a dental professional in the indications for use.

INDICATIONS FOR USE

- Direct restorative for Class I, II, III, IV and V cavities
- Fissure sealant
- Sealing hypersensitive areas
- Repair of (indirect aesthetic restorations, temporary crown & bridge, defect margins when margins are in enamel
- Blocking out undercuts
- Liner or base
- Core Build-up

CONTRAINDICATIONS

- Direct pulp capping.
- Avoid use of this product in patients with known allergies to methacrylate monomer or methacrylate polymer.

PRODUCT DESCRIPTION

G-ænial Universal Injectable is a light-cure, radio-opaque restorative material to be used intra-orally and classified as a Type 1 and Class 2 (Group 1) per ISO standard 4049.

This material has a radiopacity equivalent to 2.5-3.0 mm of aluminum (dentin = 1 mm, enamel = 2 mm).

The particle size of inorganic fillers range is 0.01-0.5 µm. The total amount of inorganic filler is approximately 46 vol%.

DIRECTIONS FOR USE




- Shade Selection
 Select shade from 16 shades of XBW, BW, A1, A2, A3, A3.5, A4, B1, B2, CV, CVD, AO1, AO2, AO3, JE, AE.
 A, B, C, D shades are based on Vita® Shade.

**Vita® is a registered trademark of Vita Zahnfabrik, Bad Säckingen, Germany.*

- Cavity Preparation
 Prepare cavity using standard techniques. Dry by gently blowing with oil free air. Rubber dam is recommended to isolate the prepared tooth from contamination with saliva, blood or sulcus fluid.
 Note:
 For pulp capping, use calcium hydroxide.

- Bonding Treatment
 For bonding G-ænial Universal Injectable to enamel and/or dentin, use a light-cured bonding system such as G-Premio BOND, G-ænial Bond or G-BOND. Follow manufacturer's instructions.
- Placement of G-ænial Universal Injectable

- Dispensing from a syringe
 1. Hold the syringe upright and remove the wing cap by turning counterclockwise (Fig. 1).
 2. Promptly and securely attach the dispensing tip to the syringe by turning clockwise.
 3. Place the light protective cap until use.
 4. Remove cap from the dispensing tip.
 Note:
 Take care not to attach the dispensing tip too tightly. This may damage its screw. In order to ensure a light connection make sure that the threads are free of remnant material.
 Prior to dispensing the material into the cavity, make sure to check the following:
 a. Dispensing tip is securely attached to the syringe.
 b. Any possible air from the dispensing tip is removed by gently pushing forward the syringe plunger with the tip pointing upwards until material reaches the opening of the tip (Fig. 2).
 5. Place the dispensing tip as close as possible to the prepared cavity, and slowly push the plunger to extrude material.
 Alternatively, dispense material onto a mixing pad and transfer to the cavity using a suitable instrument.
 Note:
 a. If plugged, remove the dispensing tip and extrude a small amount of material directly from the syringe.
 b. Minimise exposure to ambient light. Ambient light can shorten the manipulation time.
 c. After use, immediately remove and dispose of the dispensing tip and tightly close the syringe with the wing cap.
- Dispensing from a Unitip
 Insert the G-ænial Universal Injectable Unitip into the Unitip APPLIER II or equivalent. Remove the cap and extrude material directly into the prepared cavity. Use steady pressure (Fig. 3). Maintain pressure on the applier handle while removing the Unitip APPLIER II and Unitip from the mouth. This prevents the Unitip from coming loose from the applier.
 Note:
 a. The material may be difficult to extrude immediately after removing from cold storage. Prior to use, let material stand for a few minutes at room temperature.
 b. For first use, apply pressure slowly and gradually to control extrusion of material.

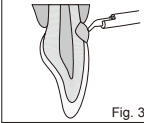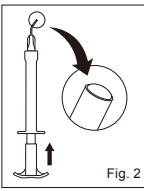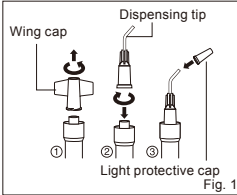- Finishing and Polishing
 Finish and polish using standard techniques.

Clinical Hints

In order to inject effectively, use the surface tension of the material to ensure uniformity across the entire surface of the restoration during build up. Once the required amount has been injected, release the pressure on the plunger and withdraw the syringe in a direction perpendicular to the surface. This will allow the material to separate from the dispensing tip and provides a smooth surface over the restoration.

- Light-curing
 Light-cure G-ænial Universal Injectable using a light curing unit. Keep light guide as close as possible to the surface. Refer to the following chart for Irradiation Time and Effective Depth of Cure.

Irradiation Time and Effective Depth of Cure (ISO 4049)	10 sec. (High power LED) (>1200 mW/cm ²)
	20 sec. (Halogen/LED) (700-1200 mW/cm ²)
A1, A2, A3, B1, B2, JE, AE	2.5 mm
XBW, BW, A3.5, A4, AO1, AO2, CV, CVD	2.0 mm
AO3	1.5 mm

The effective wavelength range of each dental curing unit must cover 450-480 nm.

- Note:
- The material should be placed and light-cured in layers. For maximum layer thickness, please consult above table.
 - Lower light intensity may cause insufficient curing or discoloration of the material.

- Finishing and Polishing
 Finish and polish using standard techniques.

SHADES

XBW, BW, A1, A2, A3, A3.5, A4, B1, B2, CV, CVD, AO1, AO2, AO3, JE, AE

STORAGE

Recommended for optimal performance, store in a cool and dark place (4-25°C / 39-77°F).

PACKAGES

- Syringes
 1. Syringe 1.7 g (1.0 mL) x 1, Dispensing Tip Needle x 5, Dispensing Tip Long Needle x 5, Light-protective cap x 1
 2. Syringe 1.7 g (1.0 mL) x 2, Dispensing Tip Needle x 10, Dispensing Tip Long Needle x 10, Light-protective cap x 2
 3. Dispensing Tip Refill, Dispensing Tip Long Needle x 30, Light-protective cap x 2

II. Unitips

- Unitip 0.27 g (0.16 mL) x 15

CAUTION

- In case of contact with oral tissue or skin, remove immediately with cotton or a sponge soaked in alcohol. Flush with water. To avoid contact, a rubber dam and/or COCOA BUTTER can be used to isolate the operation field from the skin or oral tissue.
- In case of contact with eyes, flush immediately with water and seek medical attention.
- Take care to avoid ingestion of the material.
- Prior to starting the filling process, partially extrude paste outside of the patient's mouth and away from the patient to ensure paste is flowing properly.
- Do not mix with other similar products.
- Dispensing tips are for single use only. To prevent cross contamination between patients, do not reuse the dispensing tip. The dispensing tip cannot be sterilised in an autoclave or chemically.
- For infection control reasons, Unitips are for single use only.
- Wear plastic or rubber gloves during operation to avoid direct contact with air inhibited resin layers in order to prevent possible sensitivity.
- Wear protective eye glasses during light curing.
- When polishing the polymerised material, use a dust collector and wear dust mask to avoid inhalation of cutting dust.
- Avoid getting material on clothing.
- In case of contact with unintended areas of tooth or prosthetic appliances, remove with instrument, sponge or cotton pellet before light-curing.
- Do not use G-ænial Universal Injectable in combination with eugenol containing materials as eugenol may hinder G-ænial Universal Injectable from setting.
- Do not wipe the dispensing tip with excessive alcohol as this may cause contamination of the paste.
- Dispose of all wastes according to local regulations.
- In rare cases the product may cause sensitivity in some people. If any such reactions are experienced, discontinue the use of the product and refer to a physician.
- Personal protective equipment (PPE) such as gloves, face masks and safety eyewear should always be worn.

Some products referenced in the present IFU may be classified as hazardous according to GHS. Always familiarise yourself with the safety data sheets available at:
 https://www.gc.dental/america
 They can also be obtained from your supplier.

CLEANING AND DISINFECTING:
 MULTI-USE DELIVERY SYSTEMS: to avoid cross-contamination between patients this device requires mid-level disinfection. Immediately after use inspect device and label for deterioration. Discard device if damaged.
 DO NOT IMMERGE: Thoroughly clean device to prevent drying and accumulation of contaminants. Disinfect with a mid-level registered healthcare-grade infection control product according to regional / national guidelines.

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