GC Fuji LINING LC PASTE PAK

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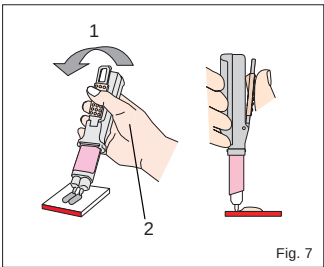
Fig. 7

GC Fuji LINING LC PASTE PAK RADIOPAQUE LIGHT CURED GLASS IONOMER LINING CEMENT

For use only by a dental professional in the recommended indications.

GC Fuji LINING LC PASTE PAK is the first radiopaque, light cured resin-modified glass ionomer lining cement available in paste-paste presentation. GC Fuji LINING LC PASTE PAK is available in a cartridge pack and can be dispensed easily and correctly using a special dispenser called Paste Pak Dispenser.

RECOMMENDED INDICATIONS

As a base or liner in prepared cavities.

CONTRAINDICATIONS

1. Pulp capping.
2. In rare cases the product may cause sensitivity to some persons. If such reactions are experienced, discontinue the use of the product and refer to a physician.

IDENTIFICATION OF PARTS

- (Paste Pak Cartridge)
1. Cover
 2. Inner plunger
 3. Tip (Paste Pak Dispenser)
 4. Piston
 5. Finger grip Long
 6. Lever
 7. Slide block
 8. Release lever (Fig. 1)

Paste / Paste Ratio (g/g)	1,0 / 1,1
Mixing Time (sec.)	10"
Working Time (min., sec.)	2'15"
Light Curing Time (sec.) (Halogen)	20"
Depth of Cure (mm)	1,3

Test conditions : Temperature (23 +/-1°C)
Relative humidity (50 +/-10%)
ISO 9917-2 : 1998 (E) (Light-activated cements)(Type I)

DIRECTIONS FOR USE

1. **PASTE PAK CARTRIDGE LOADING**
 - a) Make sure that the piston is completely retracted into the dispenser. (Fig. 2)
Note :
To retract the piston, hold both the slide block and release lever (on the underside of the dispenser) with fingers and pull them back together.
b) Load the cartridge into the Paste Pak Dispenser. (Fig. 3)
1) Ensure the ▲ mark on the cartridge lines up with the ▲ mark on the Paste Pak Dispenser as illustrated.
2) Then turn the cartridge fully to the right position.
Note :
1) If the cartridge is stored in a refrigerator, before use leave it at room temperature (23°C / 73.4°F) for 30 minutes.
2) To remove the cartridge, retract the piston fully and turn the cartridge back to the marked positions.
c) After cartridge loading, slide the release lever forward until it stops. (Fig. 4)
1) Push the front of the release lever.
2) Do not push the upper slide.

2. **DISPENSING**
 - a) Remove the cartridge cover. (Fig. 5)
b) Depress the lever to dispense the required amounts of pastes onto the mixing pad.
Note :
1) Before dispensing the first material from a new cartridge or when dispensing after a long interval, bleed the pastes in order to prevent the incorporation of air bubbles into the material.
2) The finger grip can be moved forwards or backwards to change the quantity of material to be dispensed.
Position of Finger Grip Long for lining. (Fig. 6)
c) While depressing the lever, level and cut off the extruded material from the tips of the cartridge by moving the cartridge and dispenser to an upright position on the mixing pad. (Fig. 7)
1. Level and cut off material.
2. Depressing the lever.
- When the lever is released after dispensing, residual pastes at the tips are retracted into the cartridge. Replace the cartridge cover.
Note :
1) If any surplus paste is retained on the cartridge tips after dispensing, remove with gauze, etc.
2) The amount of residual material in the cartridge can be confirmed by the position of the slide block on the underside of the dispenser. As the cartridge has a translucent tip on one of the barrels, the residual amount can also be confirmed by the position of the inner plunger seen through the tip.
3) Do not drop the dispenser with the cartridge. Otherwise the cartridge may be damaged.

3. **MIXING**

After dispensing, incorporate and spread the material out in a thin layer on the mixing pad using a plastic spatula. Mix thoroughly with lapping strokes, for 10 seconds. Be careful not to incorporate air bubbles.
Note :
1) The working time is 2 minutes 15 seconds from the start of mixing at 23°C (73.4°F).
Higher temperatures will shorten working time.
2) Do not mix pastes with other materials.

4. **RESTORATIVE TECHNIQUE**
 - a) Prepare tooth using standard techniques. Extensive mechanical retention is not necessary.
 - b) Wash and dry but DO NOT DESICCATE. For pulp capping use calcium hydroxide.
 - c) Mix the cement (as above).
 - d) Remove surface moisture, but DO NOT DESICCATE.
 - e) Transfer cement to the preparation using a syringe or suitable placement instrument, covering dentine up to the dentino-enamel junction.
 - f) Light cure with a halogen, plasma arc or LED light curing device*. Place light source as closely as possible to the

* Irradiation time (Halogen)	GC Co-bee, Coe Lunar TA	20 sec.
(Plasma arc)	GC Flipo	10 sec.
(LED)	GC e-Light (Mode : Fast cure)	18 sec.

- cement surface.
Note :
If application exceeds 2 mm in depth, use a layering technique.
g) Proceed with standard composite or amalgam restorative techniques.

- Cleaning / Sterilization of Paste Pak Dispenser**
1. Clean the dispenser and the finger grip with gauze soaked in alcohol.
 2. The dispenser may be steam autoclaved at 121-135°C (250-275°F) for 10-30 minutes or disinfected using alcohol. (Please remove the finger grip and the cartridge before autoclaving the dispenser.) The use of various chemical disinfectants / sterilants may cause damage to the dispenser and therefore is not recommended.
Note :
After autoclaving, the release lever of the dispenser may not easily be moved forwards and backwards by fingers. In such a case, slide the release lever forward until it stops. Apply a very small quantity of petroleum jelly onto the entire surface of the piston. Do not apply excess petroleum jelly. Then move the release lever back and forth a few times to ensure smooth movement.

SHADE

Universal

STORAGE

Store in a cool and dark place (4 - 25°C) (39.2 - 77.0°F)
(Shelf life : 2 years from date of manufacture)

PACKAGE

Paste Pak Cartridge (7.0g / 4.7mL) (1) with mixing pad (No.22)
Option : Paste Pak Dispenser (1)

CAUTION

1. In case of contact with oral tissues or skin, remove immediately with cotton, gauze or a sponge soaked in alcohol. Flush with water.
2. In case of contact with eyes, flush immediately with water and seek medical attention.
3. In rare cases metal may cause sensitivity to some persons. If such reactions are experienced, discontinue the use of Paste Pak Dispenser (nickel-plated aluminum) and refer to a physician.
4. Store Paste Pak Dispenser away from chemicals such as iodine tincture containing halogen element. Otherwise the dispenser may discolor.

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