

initial™ Zirconia Disk Multilayer Elite

The top choice for the highest standards.

initial™ Zirconia Disk Multilayer Elite

The top choice for the highest standards.

Benefits at a glance

- Superior aesthetics with a **combined shade and translucency gradient**.
- **High strength** throughout all layers, making it suitable for a wide range of indications.
- Limits the required finishing with a **framework in the right colour tone**.
- **Optimal synergy** with Initial Zirconia Coloring Liquid and the IQ ONE SQIN concept for high-end individualisation.
- All 16 V-shades and 3 bleach types are available in 14, 16, 20 & 25 mm thick disks with a diameter of 98.5 mm.



Perfect synergies



Initial Zirconia Coloring Liquid



Initial IQ ONE SQIN



Initial Zr-FS



GC EUROPE N.V.

Head Office
Researchpark,
Haasrode-Leuven 1240
Interleuvenlaan 33, B-3001 Leuven
Tel. +32 16 74 10 00
Fax. +32 16 40 48 32
info.gce@gc.dental
<https://www.gc.dental/europe>

GC UNITED KINGDOM Ltd.

Coopers Court Newport Pagnell
Buckinghamshire
MK16 8JS
United Kingdom
Tel. +44 1908 218 999
Fax. +44 1908 218 900
info.uk@gc.dental
<https://www.gc.dental/europe/en-GB>

z L L F E N 7 9 3 0 9 / 2 4



The top choice for the highest standards

Engineered to deliver the perfect balance of high strength and exceptional aesthetics, Initial Zirconia Disk Multilayer Elite is the ideal choice for both anterior and posterior restorations.




SUPERIOR AESTHETICS
True-to-nature shade and translucency gradient for anterior and posterior restorations.

HIGH STRENGTH
Classified as a Class V according to DIN EN ISO 6872 with a mean value of > 1.100 MPa.

OMNI-FUNCTIONAL
Versatile and suitable for a broad range of clinical indications. Fast sintering or regular sintering.


Natural gradient

Stunningly natural: a seamless gradient from cervical to incisal, in shade as well as translucency gradients, mimicking the natural tooth structure in the best possible way.

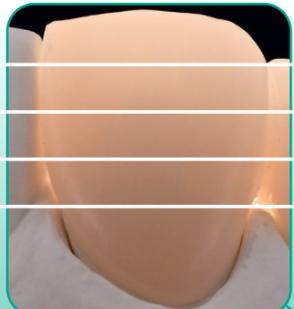Translucency >49% Strength 800 MPa		5Y INCISAL 20%
		4-5Y INTERLAYER 1 15%
		4-5Y INTERLAYER 2 15%
		4-5Y INTERLAYER 3 15%
Translucency >45% Strength 1.150 MPa		4Y CERVICAL 35%

Physical properties

- Zirconium dioxide Type II Class V DIN EN ISO 6872
- Flexural strength: typical mean value ≥ 1.100 MPa
- High translucency gradient from >45% to >49%
- Fast sintering for crowns and bridges up to 3 units

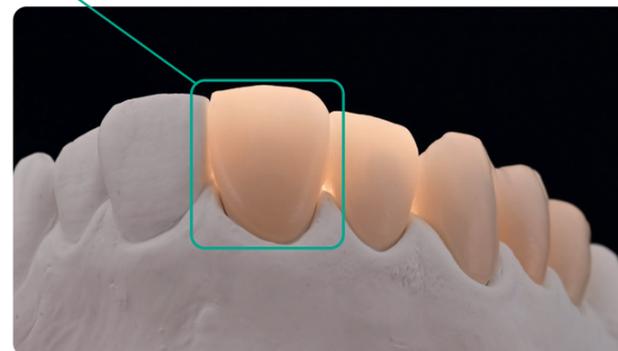


Made for qualitative productivity

Suitable for high productivity without compromising on quality.

SCAN & DESIGN with Aadvia Lab Scan 3
Gesture-controlled scanner with high precision and smart scan features



2D INNER STAIN with Initial Zirconia Coloring Liquid
Refine with effect colours for the ultimate sophistication.

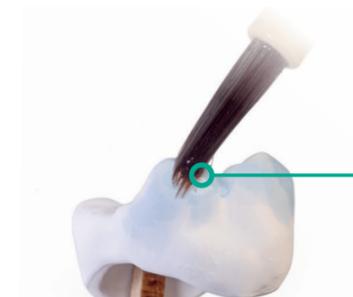

MICROLAYERING with Initial IQ ONE SQIN
3D PAINT with Initial IQ Lustre Pastes ONE
Colour and fluorescence - either for purely monolithic work or with micro-cutback in combination with IQ SQIN.



FORM & TEXTURE with Initial IQ SQIN
Vitality, natural texture and gloss within a microlayer.

CLASSIC VENEERING with Initial Zr-FS
Unsurpassed optical effects and high stability after multiple firings for frameworks with classic cutback.



CEMENT with FujiCEM Evolve
The preferred choice for zirconia, easy to use and ideal for juxtagingival or subgingival margins.

