

Essentia[®] Universal shade

Universal resin composite

One shade fits all



Since 1921
100 years of Quality in Dental




Essentia® Universal shade from GC

One shade fits all


With **Essentia Universal shade**, GC throws the door to simplification wide open. Essentia Universal exhibits an excellent chameleon effect and **seamlessly blends with the surrounding tooth structure, regardless of the shade of the cavity.**

Available in **three viscosities - all strong enough to be used in all cavity classes** - Essentia Universal is the ultimate solution for posterior restorations.

Thanks to a proprietary filler treatment and to the use of ultra-fine particles, **the physical properties of Essentia HiFlo & Essentia LoFlo are even higher** than those of Essentia Universal (paste type), confirming that all products can be safely used **for all indications.**

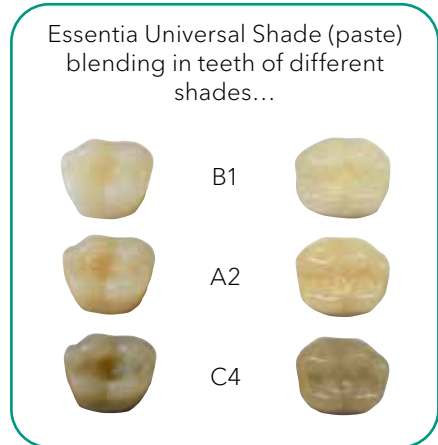

Flexural strength & 3-body wear of the three viscosities of Essentia Universal shade (Source: GCC R&D data)

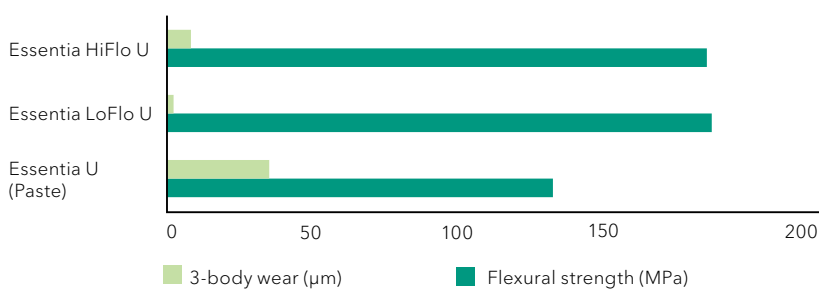

Preparation of the small Class I cavities



Restoration using **Essentia HiFlo Universal**



Post-operative view showing great blending with natural teeth




Air abrasion of the non-carious lesions



Restoration with **Essentia LoFlo Universal**




Post-operative view after polishing

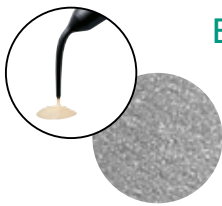

Essentia HiFlo is a fluid composite with a **high wettability**, perfect for easy adaptation to preparations and flow into narrow cavities. It offers great handling properties together with superior strength and a high radiopacity of over 300% facilitating the diagnosis of secondary caries.



Essentia LoFlo is an injectable composite with a **low flowability**, optimal when a thixotropic material is needed - for instance for cervical restorations. It promises minimal wear, ensuring long-lasting restorations even on the occlusal surface.

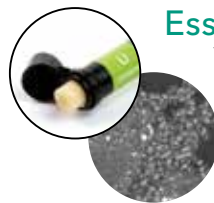

Essentia is a paste composite with a simplified 7-shade concept. Its Universal shade displays a great blending effect together with a **packable consistency**: perfect for mono-shade posterior restorations.



A direct approach was selected for this challenging post-endodontic case



Restoration with **everX Posterior** as a dentine core and **Essentia Universal** as the outer shell



Six-month recall showing good shade integration into the surrounding tissues

Essentia® Universal shade from GC

One shade fits all



Essentia HiFlo

Universal light-cured composite
with high flowability



Essentia LoFlo

Universal light-cured composite
with low flowability



Essentia

Universal light-cured composite
with a packable consistency

- 10003383 Essentia HiFlo Syringe, U (1 x 2mL)
- 10003385 Essentia LoFlo Syringe, U (1 x 2mL)
- 10003524 Essentia Syringe, U (1 x 2mL)
- 10003549 Essentia Unitips, U (15 x 0.16mL)

Also discover G-Premio BOND GC's universal bonding featuring zero compromises

- 10001457 G-Premio BOND Kit 5mL
- 10001461 G-Premio BOND Bottle Refill 5mL
- 10001458 G-Premio BOND 3-Bottle Pack
- 10001242 G-Premio BOND Unit Doses, 50 pieces (0.1mL)



GC EUROPE N.V.
Head Office
Researchpark,
Haasrode-Leuven 1240
Interleuvenlaan 33, B-3001 Leuven
Tel. +32 16 74 10 00
Fax. +32 16 40 48 32
info.gce@gc.dental
<https://www.gc.dental/europe>

GC UNITED KINGDOM Ltd.
Coopers Court Newport Pagnell
Buckinghamshire
MK16 8JS
United Kingdom
Tel. +44 1908 218 999
Fax. +44 1908 218 900
info.uk@gc.dental
<https://www.gc.dental/europe/en-GB>