



# Metal Primer Z from GC

The hero for  
all your metal based  
frameworks  
One Step Resin to Metal  
Bonding Agent


**GC**

## Looking for a bonding agent to create strong connections between metal frameworks and any type of resin based material?

With Metal Primer Z, GC brings you the successor of the proven Metal Primer II, the choice of many dental professionals since many years. With its 2 functional monomers, Metal Primer Z will create even stronger connections.

The bond to precious alloys is improved by the **MDTP** monomer.

The second monomer, **MDP**, enables a strong bond to non-precious alloys and even to zirconia for certain indications.

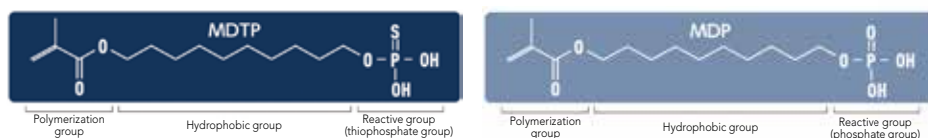

Metal Primer Z stands out, because of these features:

- A quick and easy procedure
  - No need for any additional equipment
- No heat treatment required
  - Avoiding unnecessary T.E.C. variations
- Creates a permanent bond, even in a wet environment
  - Can be applied directly in the mouth
- Universal in use
  - Is used with all dental alloys and acrylics



009295 GC Metal Primer Z, 5ml

Whenever you have to work with metal based frameworks and resin based materials, choose this all-round hero of all time.

- Application of (opaque) resin based materials to metal & zirconia for crowns and bridges
- Bonding or repair of dentures with metal bases, clasps or attachments



Metal Primer Z step by step



Sandblasting

Application

Resin/Opaque

**GC EUROPE N.V.**  
 Head Office  
 Researchpark  
 Haasrode-Leuven 1240  
 Interleuvenlaan 33  
 B-3001 Leuven  
 Tel. +32.16.74.10.00  
 Fax. +32.16.40.48.32  
 info@gceurope.com  
 http://www.gceurope.com

**GC UNITED KINGDOM Ltd.**  
 12-15, Coopers Court  
 Newport Pagnell  
 UK-Bucks. MK16 8JS  
 Tel. +44.1908.218.999  
 Fax. +44.1908.218.900  
 info@uk.gceurope.com  
 http://uk.gceurope.com

**GC**