

MI Varnish

treats hypersensitive teeth
and provides enhanced long term protection

Support, replenish
and protect with
the stroke of a brush



Support, replenish and protect with the stroke of a brush

GC introduces MI Varnish, the latest addition to your minimum intervention dentistry arsenal.

Offering superior coverage of tooth surfaces, MI Varnish **does substantially more** than your average fluoride varnish.


By delivering a powerful dose of fluoride, MI Varnish delivers an **added booster effect** of calcium and phosphate ions through its patented **Recaldent™** ingredient.

How MI Varnish works

Thanks to the world renowned Recaldent™ MI Varnish gives patients a head start in caries protection.


Without MI Varnish

In the acidic environment of the mouth calcium, phosphate and hydroxyl ions dissolve from the tooth. This will cause the tooth to become weaker and more vulnerable.

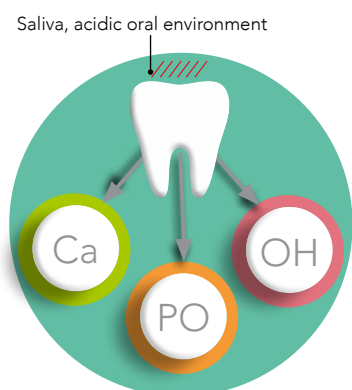


With MI Varnish

Fluoride, calcium and phosphate combine and form a strong protective coating. The neutral pH of 6,6 enhances enamel acid resistance and inhibits demineralisation.

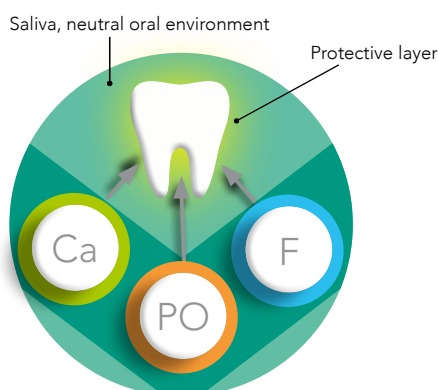



Coating with depot effect

In an acidic environment, calcium and fluoride dissolve from the covering layer. As a result, the teeth are no longer under direct attack by acids.

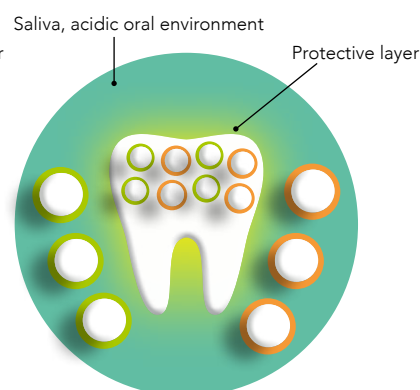

Fluoride varnishes can prevent hypersensitivity in various clinical situations:



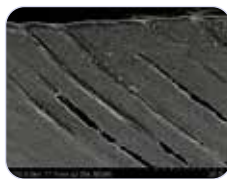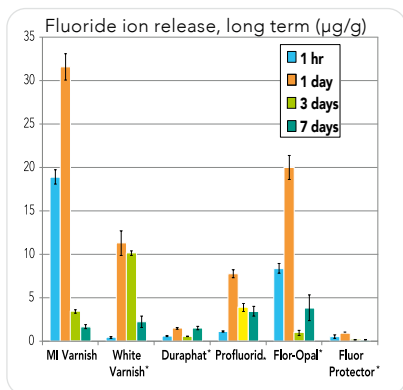

- MI Varnish penetrates and blocks dentinal tubules to give relief from hypersensitive teeth
- The high initial fluoride release minimises sensitivity and strengthens enamel

* is not a registered trademark of GC Europe



- ✓ The simple application requires no preparation or prophylaxis
- ✓ Comfortable and quick procedure completed within minutes

- ✓ Adheres quickly to the applicator brushes and flows easily in hard to reach areas
- ✓ Covers all tooth surfaces evenly and effectively



- ✓ MI Varnish sets when in contact with water or saliva, without clumping or coagulating
- ✓ Pleasant creamy texture, comfortable film thickness, even with multiple layers

- ✓ Available in 2 great tasting flavours, mint & strawberry
- ✓ Pleasant taste for patients



- ✓ Neutral shade with natural translucency
- ✓ No colour change when applied to the tooth

MI Varnish, achieving mutual benefits

MI Varnish forms part of a total MI treatment strategy



GC offers dentists the necessary tools to manage a minimum intervention continuum of care with a full suite of tailored products, ranging from diagnosis to prevention and treatment. **While MI Varnish provides high level protection and reduced sensitivity at the chairside, MI Paste can be used as follow up treatment for a patient's home care regime, extending the remineralisation effect.**



In office

Professional application of MI Varnish at intervals, determined by indication and patient needs.



At home

Daily application of Tooth Mousse or MI Paste Plus provides continued remineralisation and caries protection

RECALDENT™ and RECALDENT™ logo are trademarks of Mondelēz International group used under license. RECALDENT™ is derived from milk casein. Do not use on persons with a milk allergy. In case of allergic reaction: stop use, rinse mouth with water, and seek medical advice.

GC EUROPE N.V.

Head Office
Researchpark,
Haasrode-Leuven 1240
Interleuvenlaan 33, B-3001 Leuven
Tel. +32 16 74 10 00
Fax. +32 16 40 48 32
info.gce@gc.dental
<https://www.gc.dental/europe>

MI Treatment Ordering information

MI Varnish Intro pack

10003389 MI Varnish Intro pack,
10 unit doses MI Varnish (5 Strawberry + 5 Mint),
0.50 ml per unit dose, and 20 brushes

MI Varnish Refill pack, 35 unit doses, 0.50 ml per unit dose, and 50 brushes

10003390 Strawberry
10003391 Mint

MI Varnish Clinic pack, 100 unit doses

MI Varnish, 0.50 ml per unit dose and 100 brushes

10003392 Strawberry
10003393 Mint

GC MI Paste Plus

10003270 GC MI Paste Plus Promotional 5-Pack Box with
5 Tubes, 40 g (35 ml) per Tube, 1 of each flavour
(Melon, Strawberry, Tutti-Frutti, Mint, Vanilla)
10001749 GC MI Paste Plus, 10-Pack Assorted Box with
10 Tubes, 40 g (35 ml) per Tube, 2 of each flavour
(Strawberry, Vanilla, Tutti-Frutti, Mint, Melon)
40 g (35 ml) per Unit

GC MI Paste Plus, 10-Pack Box with 10 Tubes of a single flavour, 40 g (35 ml) per Tube

10001750 Mint
10001751 Strawberry
10001752 Melon
10001753 Vanilla
10001754 Tutti Frutti

Tooth Mousse

10003205 GC Tooth Mousse, Promotional 5-Pack Box with
5 Tubes, 40 g (35 ml) per Tube, 1 of each flavour
(Melon, Strawberry, Tutti-Frutti, Mint, Vanilla)
10001743 GC Tooth Mousse, 10-Pack Assorted Box of
10 Tubes, 40 g (35 ml) per Tube, 2 of each flavour
(Melon, Strawberry, Tutti-Frutti, Mint, Vanilla)

GC Tooth Mousse, 10-Pack Box with 10 Tubes of a single flavour, 40 g (35 ml) per Tube

10001745 Strawberry
10001746 Melon
10001747 Mint
10001748 Tutti Frutti
10001744 Vanilla

GC UNITED KINGDOM Ltd.

Coopers Court Newport Pagnell
Buckinghamshire
MK16 8JS
United Kingdom
Tel. +44 1908 218 999
Fax. +44 1908 218 900
info.uk@gc.dental
<https://www.gc.dental/europe/en-GB>