

Chic is the answer for doctors who want an impressive ceramic bracket option for their patients without compromising the performance or strength of the bracket. Chic is manufactured by Ceramic Injection Molding (CIM) creating uniformity, strength, and a true mechanical base. Its translucency compliments the patient's tooth shade while offering optimum stain resistance.




## Design Features

- Translucent ceramic blends with enamel
- Stain resistant
- 100% mechanical base
- Strong enough to rebond
- Easy debonding with debonding indent
- Engineered for strength and aesthetics

**Hi** Aesthetic

# Chic.

## Ceramic Standard Brackets

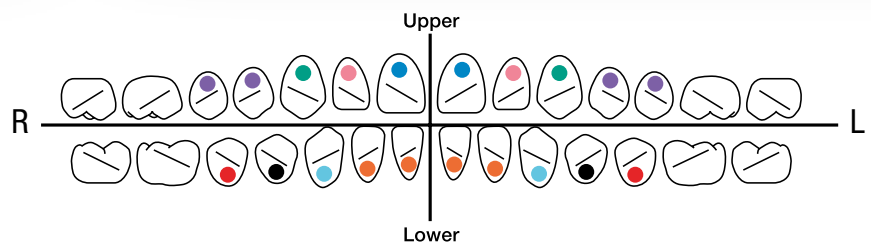

**Experts in Orthodontics**

GC Orthodontics Europe GmbH

[www.gcorthodontics.eu](http://www.gcorthodontics.eu)

# Chic.

**Rhomboid bracket**  
for easy bracket placement

**Ceramic Injection Molding**  
for uniformity, strength & precision

**Triple chamfered slot**  
for easy wire engagement and to eliminate arch wire binding & notching



**Torque in base**

**Long axis line**  
for easy bracket placement; easily removed with first brushing


**Translucent ceramic**  
for excellent aesthetics

**Milled slot**



**Excellent under tie-wing area**  
for easy ligation

**Molded mechanical base**  
for greater bracket retention and predictable debonding

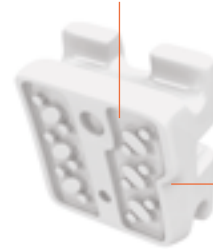

**Debonding indent**  
at the edge of the bracket for easier debonding

## Bracket combinations

You can combine Chic Ceramic brackets with LEGEND™ Rhodium brackets for a complete aesthetic result.



Debonding Chic



**Chic is manufactured by TOMY Inc.** in Yumoto, Japan. Orthodontists around the world have used and adored TOMY brackets since their first bracket system launched in December of 1961. As of today, TOMY has grown to be one of the most respected bracket manufacturers in the world.

Chic is available in a variety of prescriptions.

## Bracket & Wire Combination

Pair Chic with our Initalloy, BIO-ACTIVE, and BIO-EDGE thermal NiTi arch wires for increased patient comfort and treatment efficiency. Rhodium-coated arch wires also available for a completely discrete orthodontic bracket system.



99-14ET-1300 REV01-09.2019



**GC Orthodontics Europe GmbH**  
Harkortstraße 2 · 58339 Breckerfeld · Germany  
Phone: +49 2338 801-888 · info.gc.germany@gc.dental