


LEGEND™ mini Rhodium features the same characteristics as the standard bracket but with enhanced aesthetics. LEGEND™ mini is manufactured by metal injection (MIM) and its slot undergoes a CNC (Computer Numerical Controlled) process guaranteeing an exact slot size of .018" x .018" or .022" so that the bracket information is precisely expressed. Its anatomical mesh base ensures maximum bracket retention on the tooth.

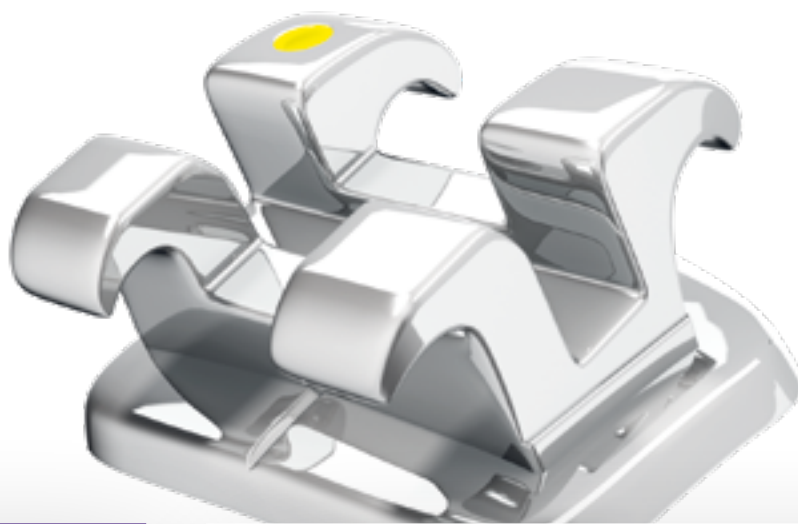

## Design Features

- Made from aerospace strength steel
- Tumbled for a highly polished and better looking bracket
- Rhodium coated for better aesthetic
- Laser-fused hook options designed for malleability and strength
- Anatomically correct bonding pad for each tooth
- In/Outs compatible between medium & mini

**Rh** Rhodium


# LEGEND mini Rhodium

Standard Metal Bracket with Rhodium



With Rhodium

Without Rhodium

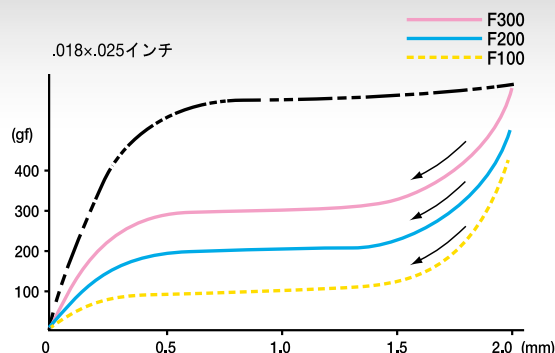

Experts in Orthodontics

GC Orthodontics Europe GmbH


[www.gorthodontics.eu](http://www.gorthodontics.eu)

# LEGEND mini Rhodium

**Rhomboid bracket**  
for easy bracket placement

**Long axis line**  
for accurate positioning

**Aerospace steel**

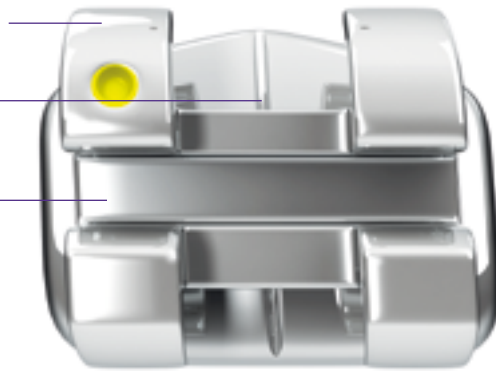

**Triple chamfered slot**  
for easy wire engagement  
and to eliminate arch wire  
binding & notching

**CNC milled slot**

**Smooth swept tie wings**  
engineered to minimize  
occlusal interference

**Large undercuts**  
for easy ligation

**Torque in Base**



**Micro-etched, 80 gauge  
mesh pad**  
for maximum bond  
retention

**Compound contour pad**  
for an anatomically  
correct fit for each tooth

## Bracket combinations

You can combine LEGEND™ Rhodium brackets with Chic Ceramic brackets for a complete aesthetic result.

**LEGEND™ is manufactured by TOMY Inc.** in Yumoto, Japan. Orthodontists around the world have used and adored TOMY brackets since their first bracket system launched in December of 1961. As of today, TOMY has grown to be one of the most respected bracket manufacturers in the world. LEGEND is available in standard and mini sizes, in metal or with a rhodium coating, in a variety of prescriptions.



## Bracket & Wire Combination

Pair LEGEND mini Rh with our Initalloy, BIO-ACTIVE and BIO-EDGE thermal NiTi arch wires for increased patient comfort and treatment efficiency. Rhodium-coated arch wires also available for a completely discrete orthodontic bracket system.



99-38ET-1300 REV01-09-2019



**GC Orthodontics Europe GmbH**  
Harkortstraße 2 · 58339 Breckerfeld · Germany  
Phone: +49 2338 801-888 · info.gc.germany@gc.dental