

# Fit and comfort for an active life



RELIN<sup>™</sup> II Soft &  
RELIN<sup>™</sup> II Extra Soft  
from GC

Soft silicone relining materials  
in cartridge delivery

**GC**

# Fit, function and comfort to keep on smiling

Irritation, extractions or implant placement will lead to changes in soft and hard tissues and **it will be needed to adapt the dentures to those new local conditions.**

Because elderly patients are getting healthier and stay active longer, they want to be able to continue smiling, speaking and enjoying eating at all times; even when wearing dentures.

**Reline II Soft and Extra Soft from GC** have been designed to create a soft layer between the denture and the mucosa which will **maintain comfort and function** for weeks and up to several months.



*Pictures courtesy of Dr D. Garcia-Baeza, Spain*

## An easy-to-apply silicone liner in cartridge delivery

The material's **ideal viscosity** makes it **flow under pressure, adapting nicely** to all anatomical details.

Its **improved trimming ability** makes the ultimate finishing of margins even easier.



Improved trimming ability  
for a better finition of margins

## Long lasting comfort and gentle relief for your patients\*

With RELINE II, temperature will not increase during setting, resulting in more comfort for your patients.

The material **remains smooth and odorless** and maintains its **dimensional stability and elasticity** for weeks or months. The result? A **durable\* & convenient solution** for denture patients.



Protection of margins with Reline II  
finishing material

## Strong adhesion, yet easy to maintain and adapt\*

The enhanced composition of **GC RELINE II Primer for resin** provides **strong adhesion to the denture\*** for a longer time. However, this strong adhesion does not make it any less flexible.

You can **add material within 3 months** after the first application, to adapt to changes of the mucosa and bone. The new **GC RELINE II Remover for resin** makes it possible to remove the relining layer and replace it very quickly.



**GC RELINE II  
Primer  
for resin**




**GC RELINE II  
Remover  
for resin**

# in all circumstances

From resilience to elasticity:  
choose the product that best fits the clinical situation

Featuring a Hardness of 32 (Shore A), **GC RELINE™ II Extra Soft** is a resilient material. It's particularly helpful to **normalise** irritated tissues after extraction or implant surgery. It can stay in place for **days or weeks\***.

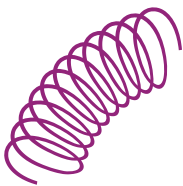

With a Shore A Hardness of 48, **GC RELINE™ II Soft** will maintain its elasticity in time and, much **like a cushion**, **spreads the load pressure** in case of anatomical problems. It is also ideal to reline **implant-supported dentures** after the healing phase. It can stay in place for **weeks or months\***.



Straightforward application & handling, throughout the entire procedure



Initial denture, to be relined following the placement of 4 implants.



Relieve the area to be relined and roughen the surface.



Apply GC RELINE II Primer on the area to be relined. Gently dry with air.



Easy cartridge delivery.



Apply GC RELINE II (Extra) Soft. Working time is 2'00".



Seat the denture in the mouth and muscle trim. Retain in situ for 5'00".



Remove the excesses using a scalpel or scissors. Trim with GC RELINE II POINT FOR TRIMMING.



Finish with GC RELINE II WHEEL FOR FINISHING.



Reline denture.

# RELINE™ II Soft & RELINE™ II Extra Soft

## Fit and comfort for an active life

Discover the all-round family of reline materials from GC



	GC TISSUE CONDITIONER	GC RELINE II Extra Soft	GC RELINE II Soft
Material	Acrylic Resin	A-Silicone	A-Silicone
Indications	<ul style="list-style-type: none"> <li>• Pressure points</li> <li>• During implant procedures</li> <li>• Soft relining</li> <li>• Tissue recovery after trauma/surgery</li> <li>• Functional impressions</li> </ul>	<ul style="list-style-type: none"> <li>• Tissue toning material</li> <li>• Sensitive mucosa</li> <li>• Pressure points</li> <li>• Normalising irritated tissues</li> </ul>	<ul style="list-style-type: none"> <li>• Implant cases</li> <li>• Pressure points</li> <li>• Bony undercuts</li> <li>• Pain reduction &amp; increase of comfort</li> <li>• Long-lasting relining in case of missing supportive tissue</li> </ul>
Time in mouth	Days/weeks/month	Days/weeks	Weeks/months

### Available packages of GC RELINE II Soft / Extra Soft

#### Intro Kit

##### GC RELINE II Soft / Extra Soft Intro Kit:

- 1x GC RELINE II Soft / Extra Soft Cartridge
- 1x GC RELINE II Primer for Resin
- 1x GC RELINE II Point for Trimming
- 1x GC RELINE II Wheel for Finishing
- 1x GC RELINE II Remover for Resin
- 6x Mixing Tips
- 3x Disposable Brushes
- 1x Plastic Spatula

#### Refills

- GC RELINE II Soft Cartridge, 48 ml
- GC RELINE II Extra Soft Cartridge, 48 ml
- GC RELINE II Primer for Resin, 10 ml
- GC RELINE II Remover For Resin, 5.5 ml
- GC RELINE II Finishing Material, 9 ml
- GC RELINE II Point for Trimming x 3
- GC RELINE II Finishing Wheel x 12
- GC RELINE II Disposable Brushes x 6

*This product is not available for purchase by the general public. Always read the label and follow the directions for use.*

\* Daisuke et al (2015). Peel Bond Strength of Silicone Denture Relining Materials. J Dent Res 94 (Spec Iss A):3487.  
<https://iadr.abstractarchives.com/abstract/15iags-2108504/peel-bond-strength-of-silicone-denture-relining-materials>

**GC Asia Pte Ltd**  
 5 Tampines Central 1,  
 Tampines Plaza 2, #06 01,  
 Singapore 529541  
 Tel: +65 67979863  
[gcasia@gc.dental](mailto:gcasia@gc.dental)  
<https://sea.gcasiadental.com>

